



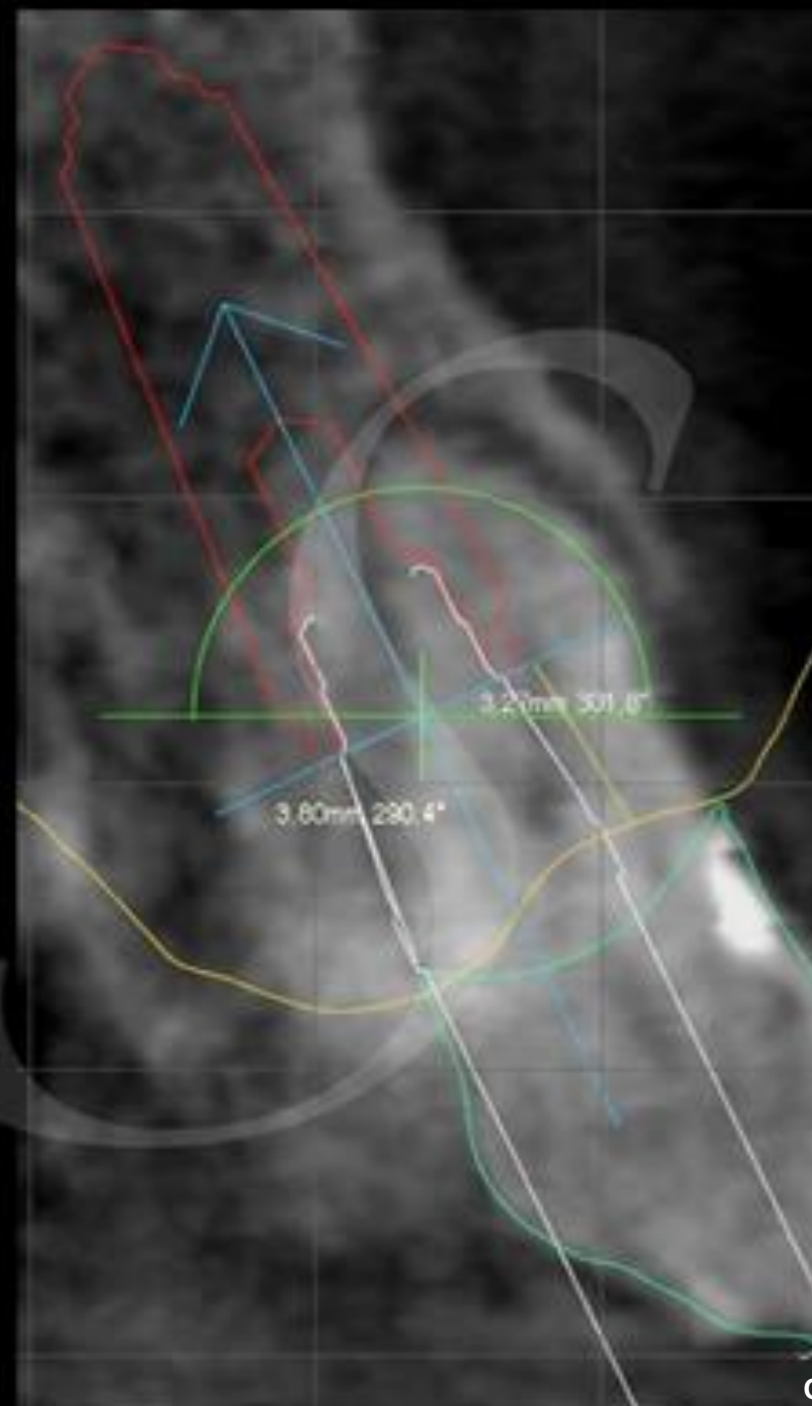


#11 Implant placement Case, Surgery using R2GATE Guide and ONE-DAY Prosthesis (Zirconia customized abutment and PMMA Crown).

Dr. Laurent Sers

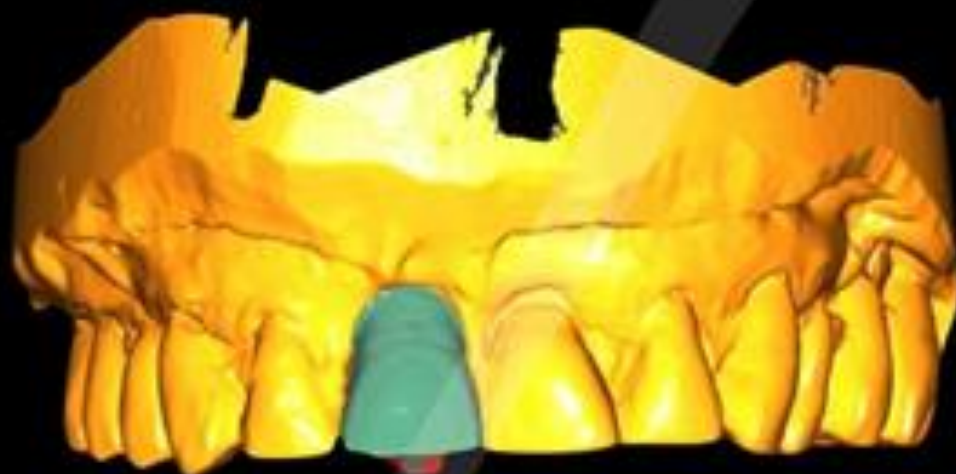

Case by Dr. Laurent Sers

Clinical Case: Diagnosis the implant by Top-Down Treatment with R2GATE.

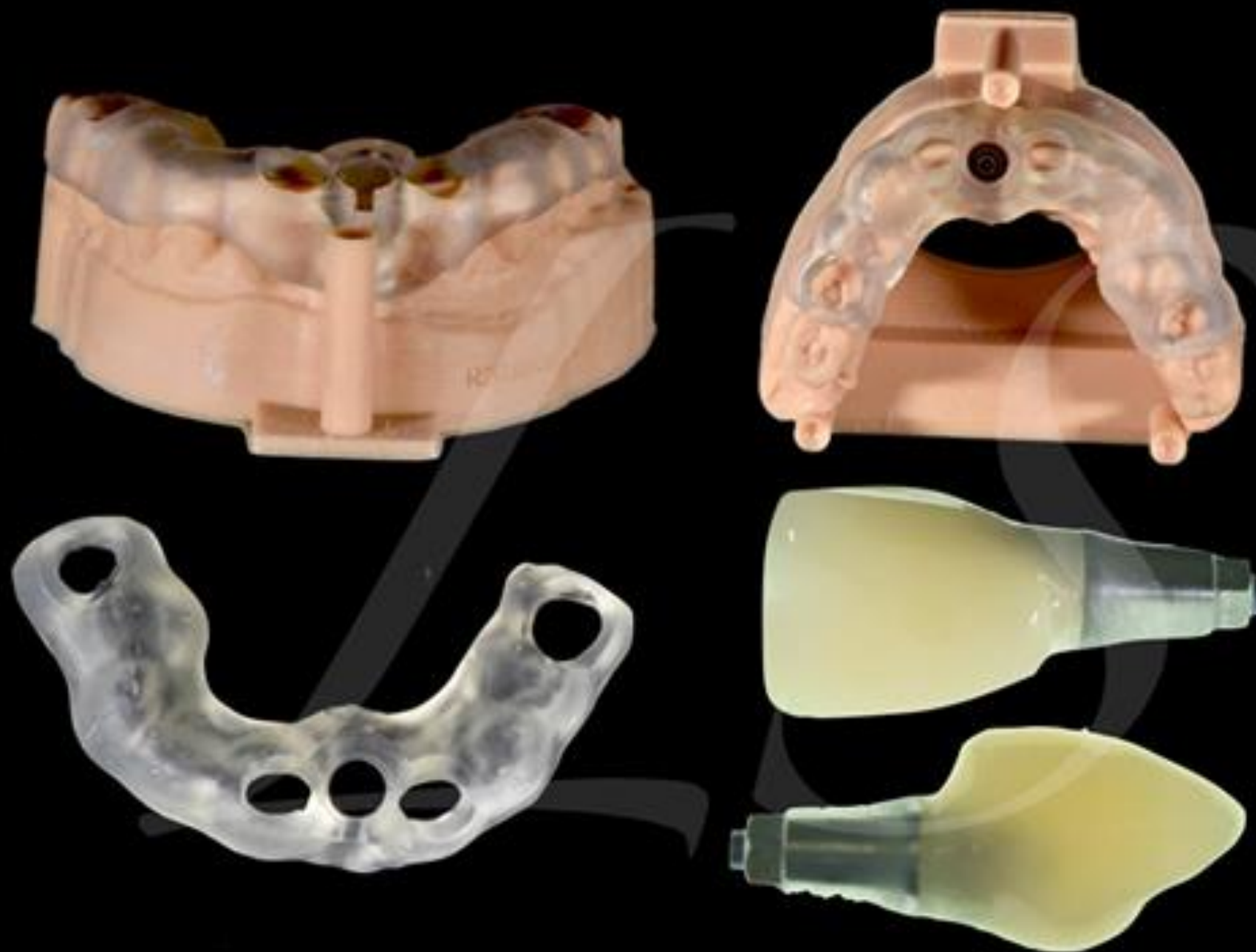

Case by Dr. Laurent Sersr

Clinical Case: Surgical guide , digital models and Prefabricated printed provisional resin crowns.




Case by Dr. Laurent Sers

Clinical Case: # 11 tooth extraction progress




Case by Dr. Laurent Sers

Clinical Case: Put the R2GATE Guide on the teeth

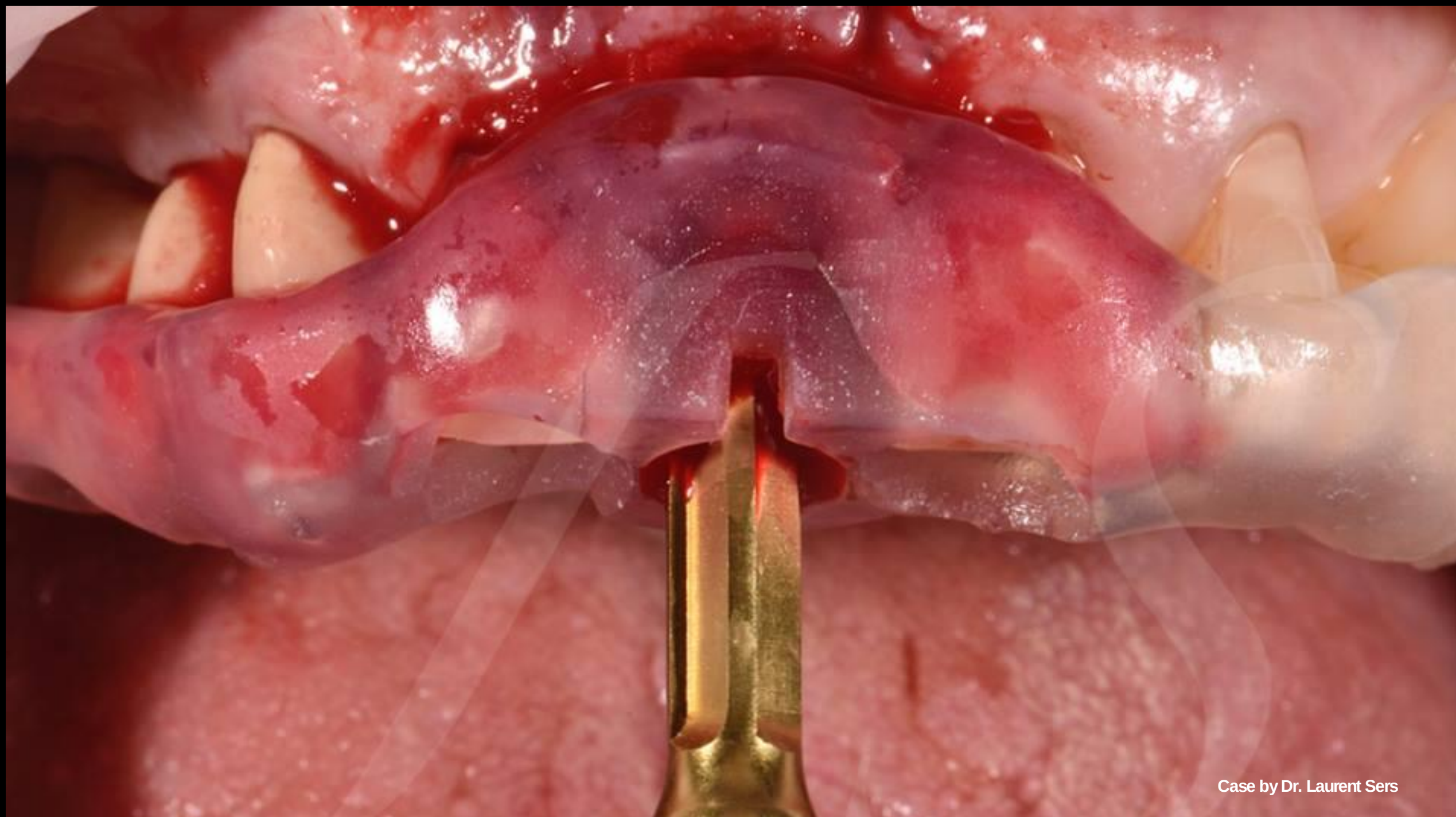

Case by Dr. Laurent Sers

Clinical Case: # 11, Final drilling 3.8x13mm

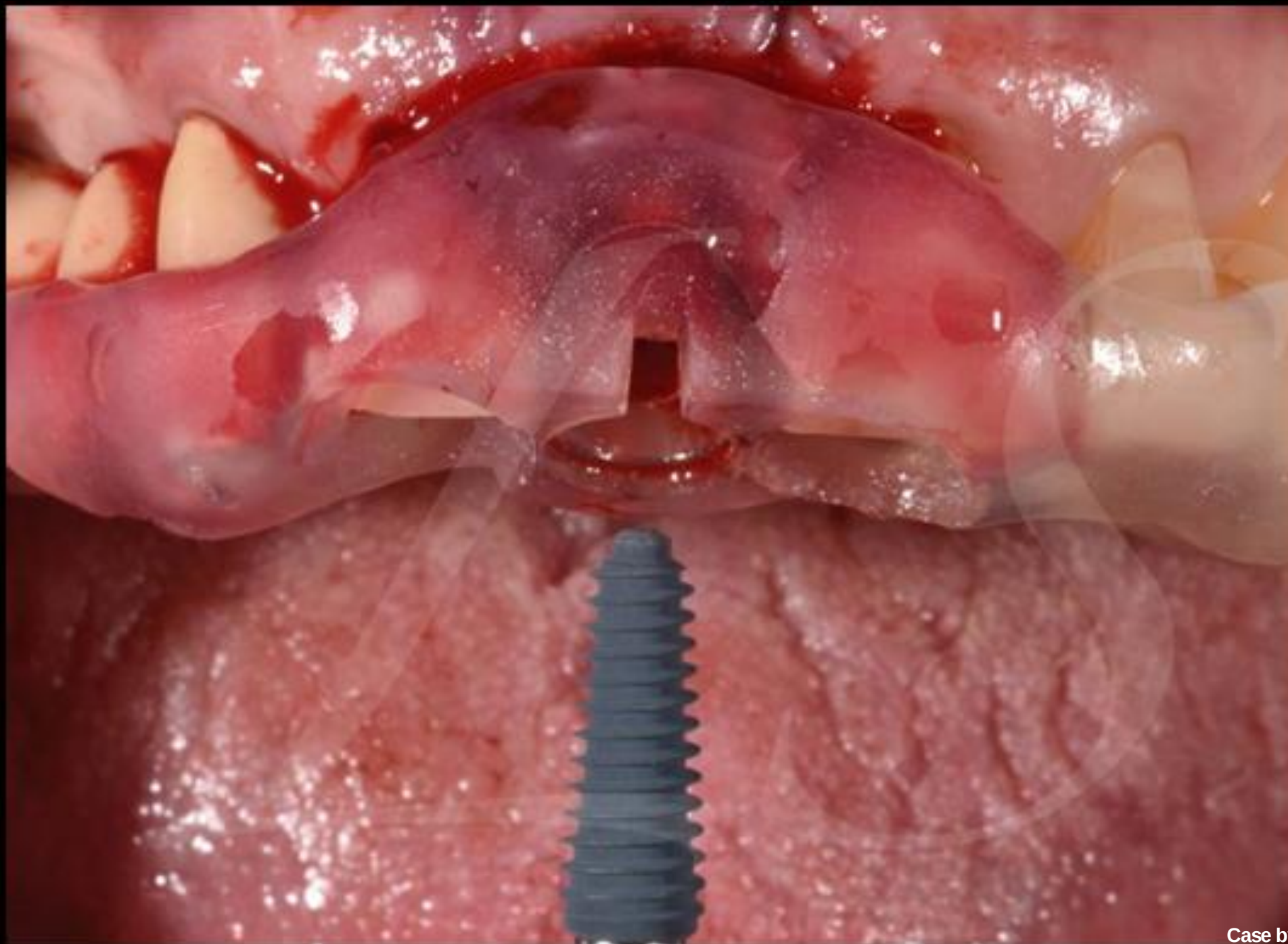


Case by Dr. Laurent Sers

Clinical Case: #11, AnyRidge 4.0x13mm Implant Placement



Case by Dr. Laurent Sers

Clinical Case: ONE-DAY Prosthesis (PMMA SCR P TYPE Crown)

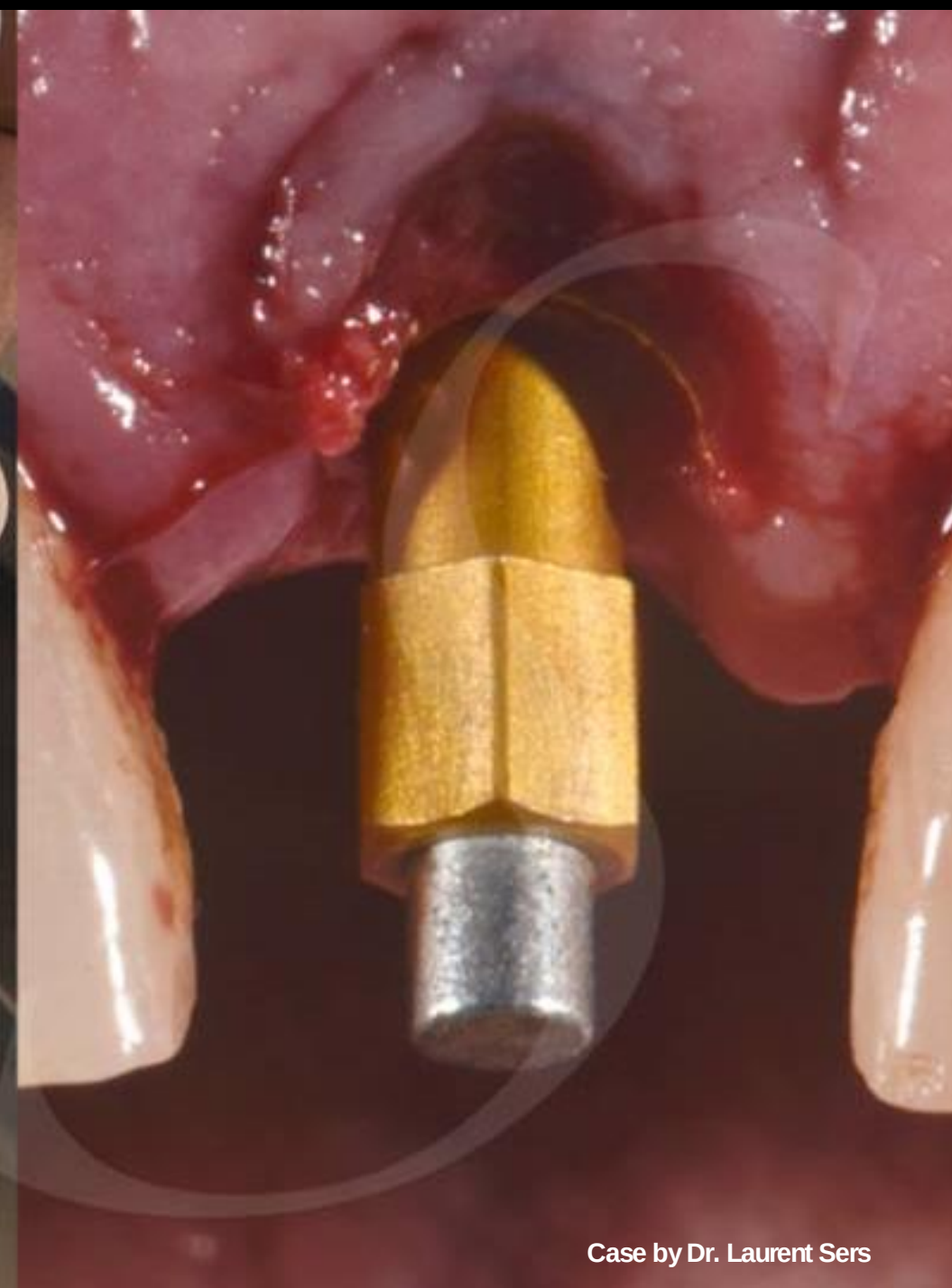

Case by Dr. Laurent Sers

Clinical Case: Measuring ISQ adequate for immediate loading (75 recommended)



Case by Dr. Laurent Sers

Clinical Case: Set the PMMA Crown and fix the gingiva to crown to prevent gingival retraction



Case by Dr. Laurent Sers

Clinical Case: Healing 2 weeks later

