





#11,21 Implant placement Case, Surgery using R2GATE Guide and ONE-DAY Prosthesis (Zirconia customized abutment and PMMA Crown).

Dr. Kwang Bum Park

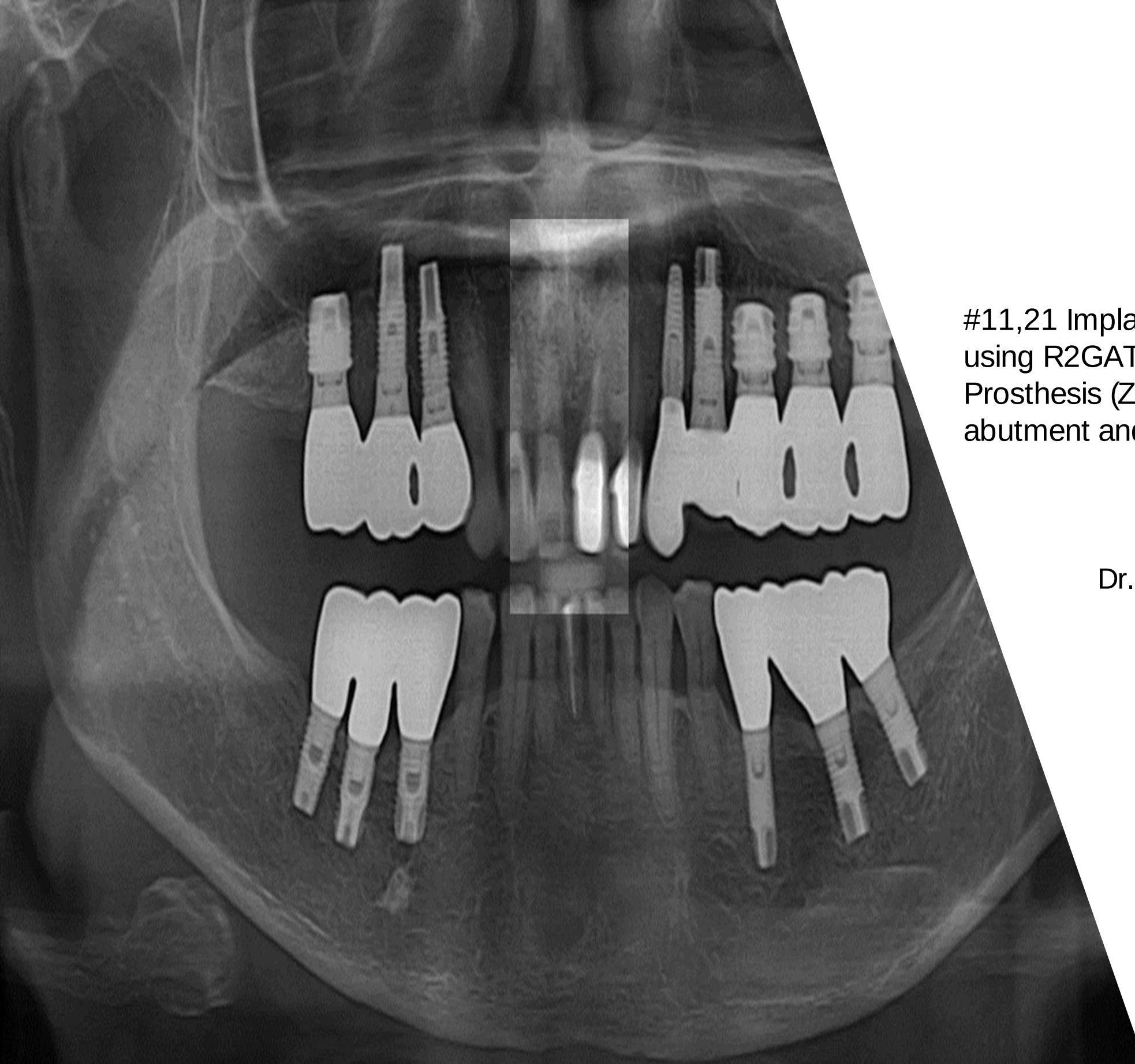

Case by Dr. Kwang Bum Park

Clinical Case: Photo to Before an extraction teeth

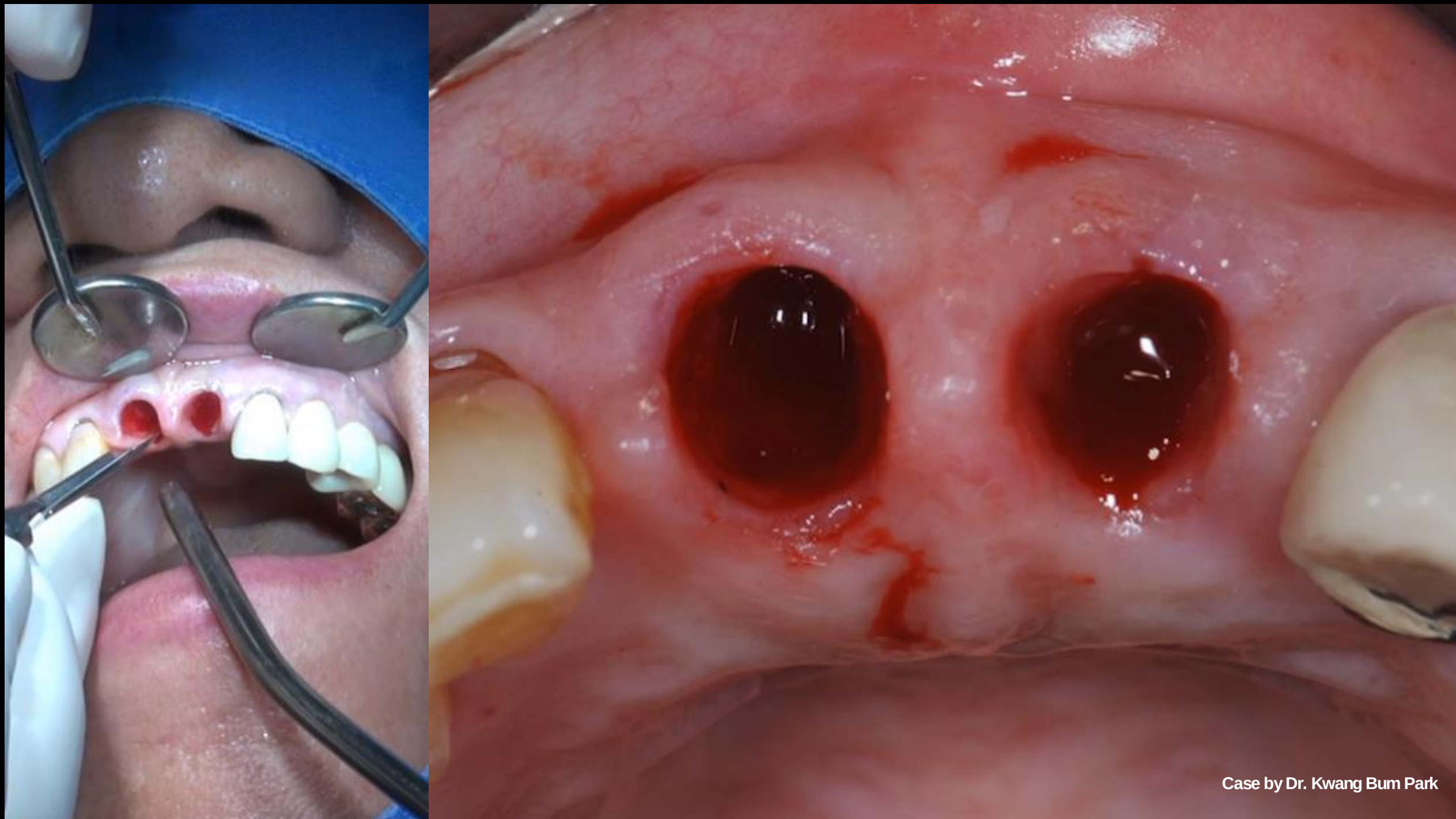

Case by Dr. Kwang Bum Park

Clinical Case: Photo to after an extraction teeth





Case by Dr. Kwang Bum Park

Clinical Case: Put the R2GATE Guide on the teeth




Clinical Case: Drilling according to drilling sequence using bone quality through digital eye: Drilling sequence(Initial drill, Second drill, 2.0x7mm, 2.0x8.5mm, 2.0x10mm, 2.0x13mm, 2.5x13mm, 2.8x13mm, 3.3x8.5mm)

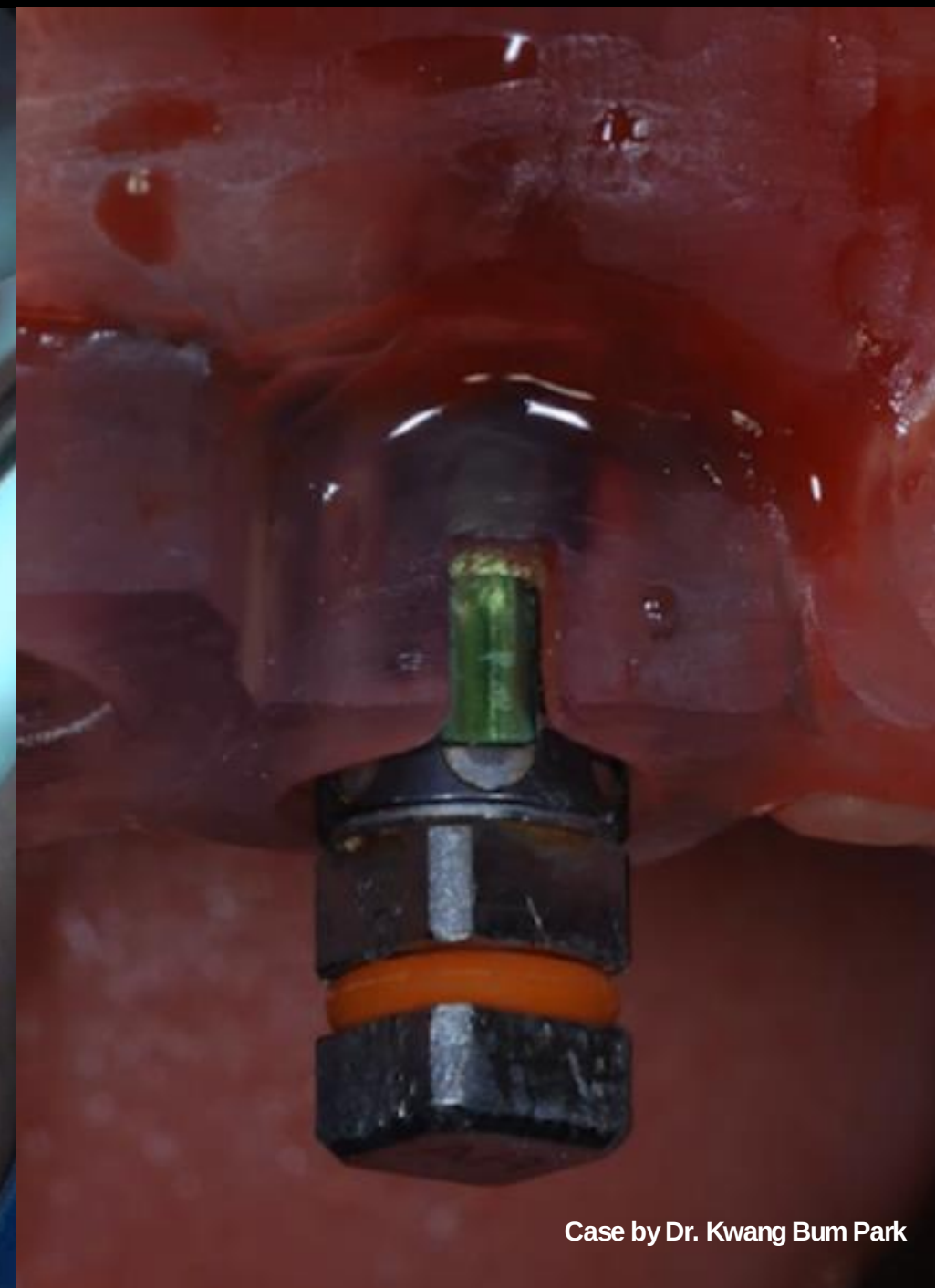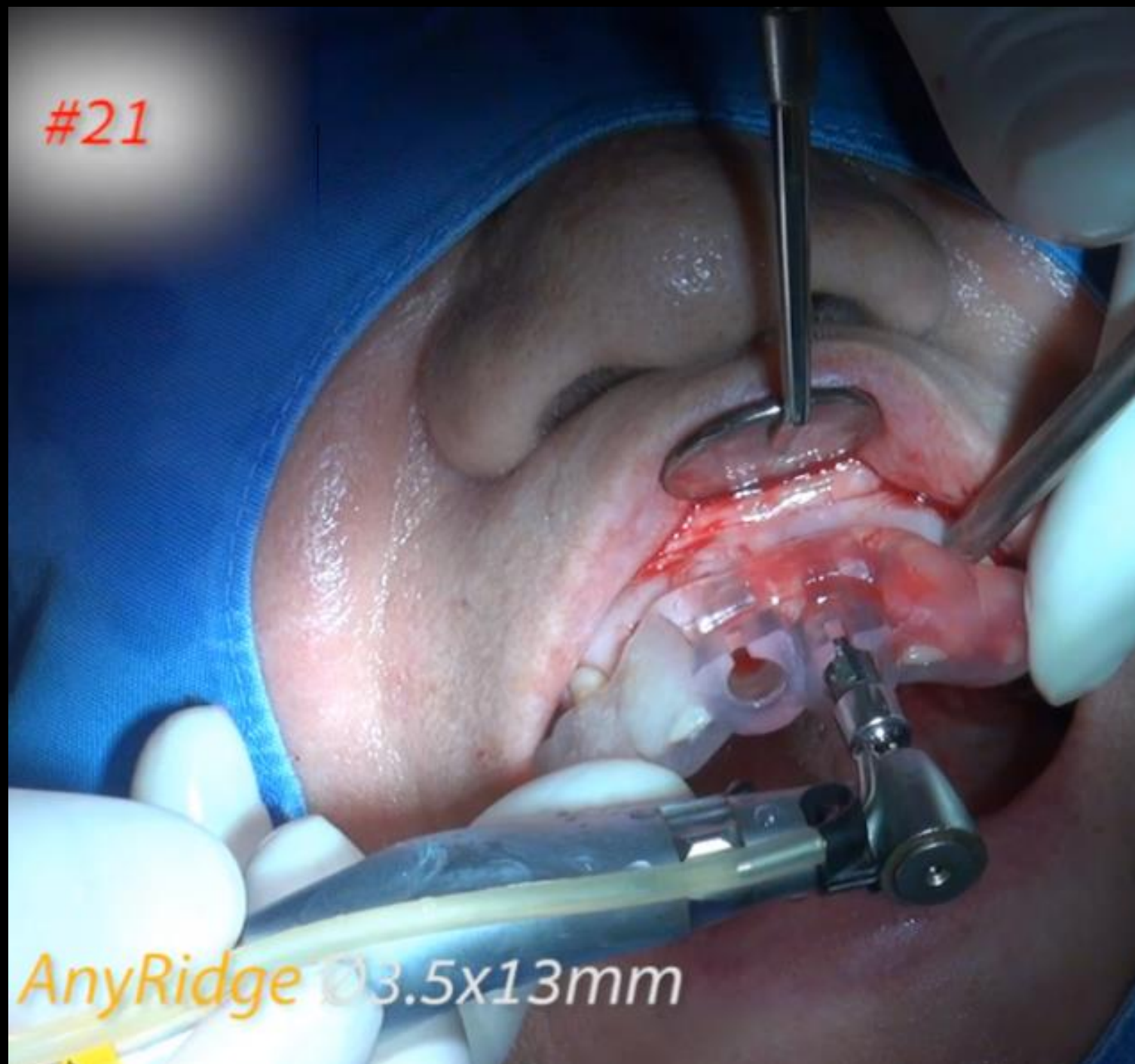

Case by Dr. Kwang Bum Park


Clinical Case: # 21, AnyRidge 3.5x13mm Implant Placement

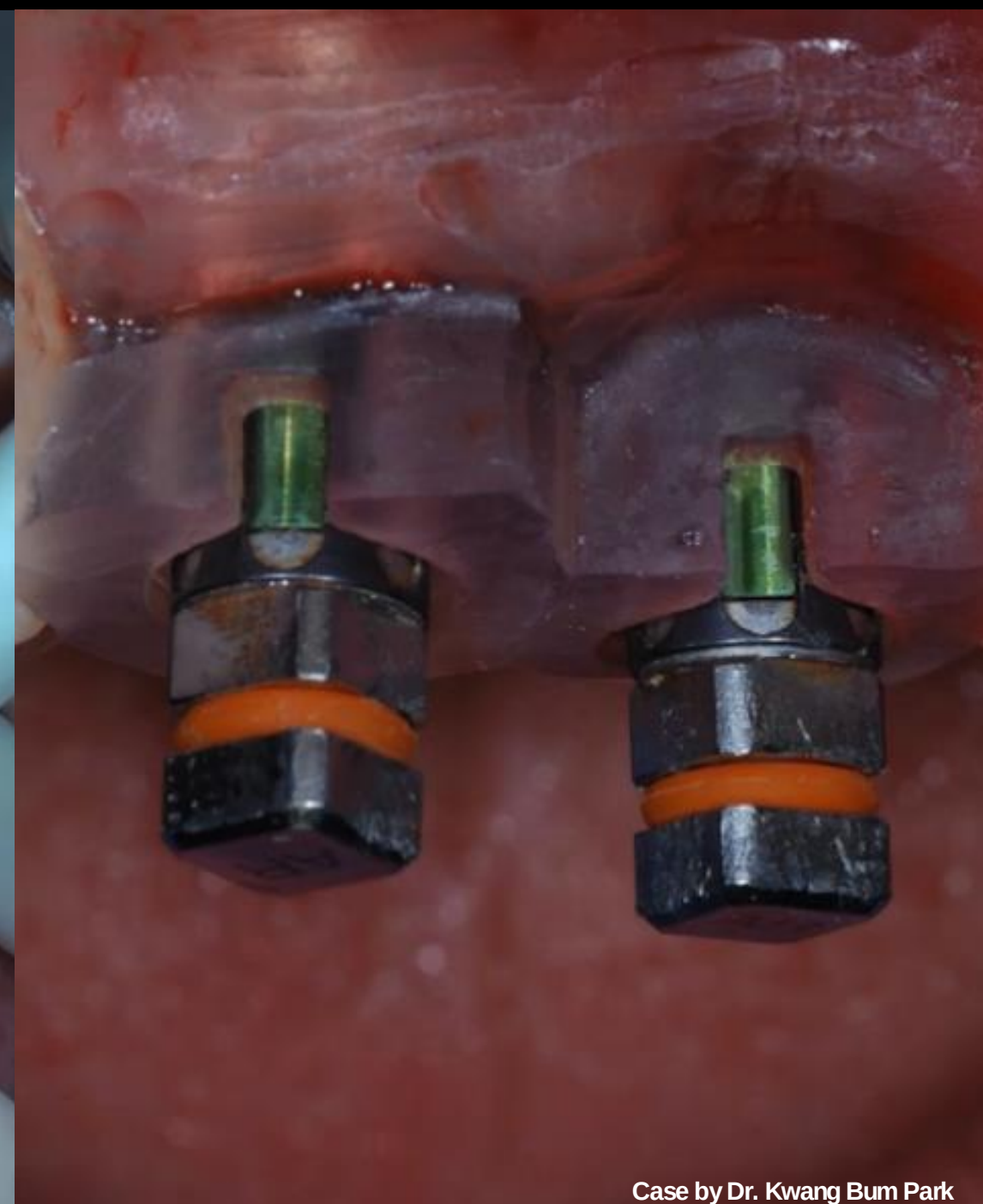

Clinical Case: # 11, AnyRidge 3.5x13mm Implant Placement



Clinical Case: # 21, Check the ISQ (B: 63, L: 70)

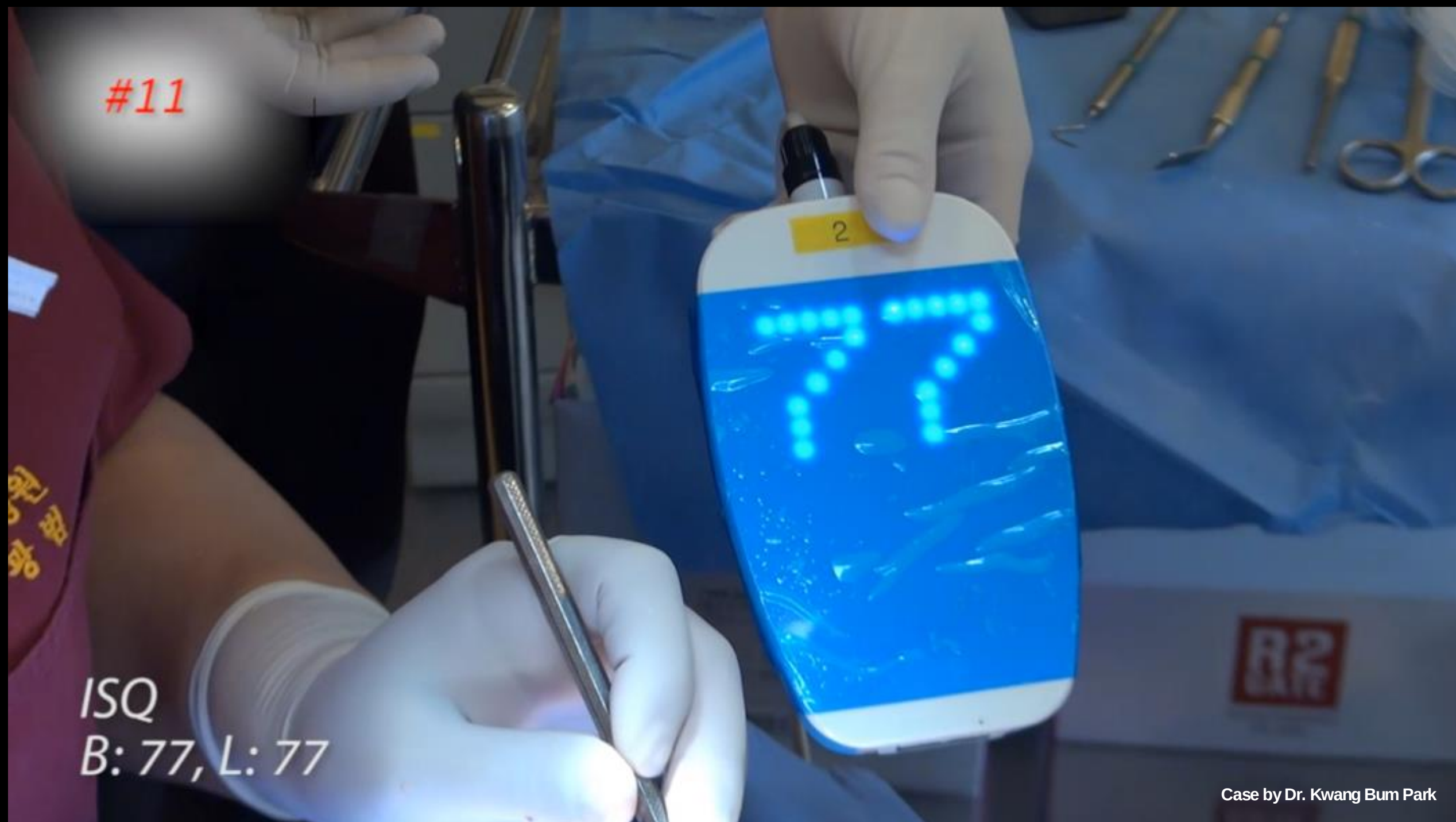

Clinical Case: # 11, Check the ISQ (B:77, L:77)



Case by Dr. Kwang Bum Park

Clinical Case: GBR using Mega-Oss




Case by Dr. Kwang Bum Park

Clinical Case: Set the Zirconia customized abutment and PMMA Crown



Case by Dr. Kwang Bum Park

Clinical Case: Set the PMMA Crown

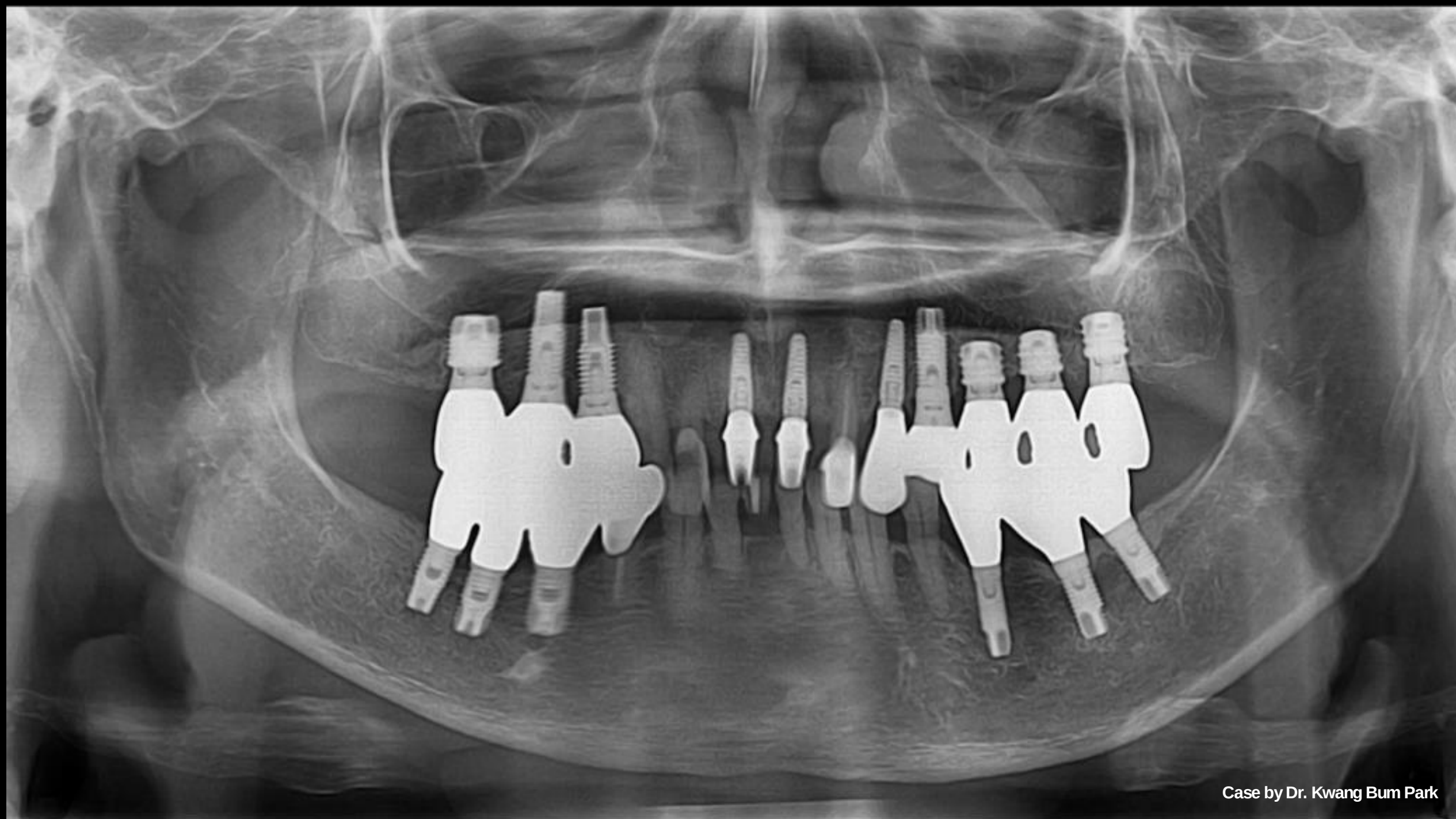

Case by Dr. Kwang Bum Park

Clinical Case: Panoramic picture after surgery

