

AnyRidge
2nd edition basic guidance

Dr. Miguel Stanley

6 MONTHS FOLLOW-UP

Case by Dr. Miguel Stanley

Fig. 11.



Fig. 12.

Case by Dr. Miguel Stanley

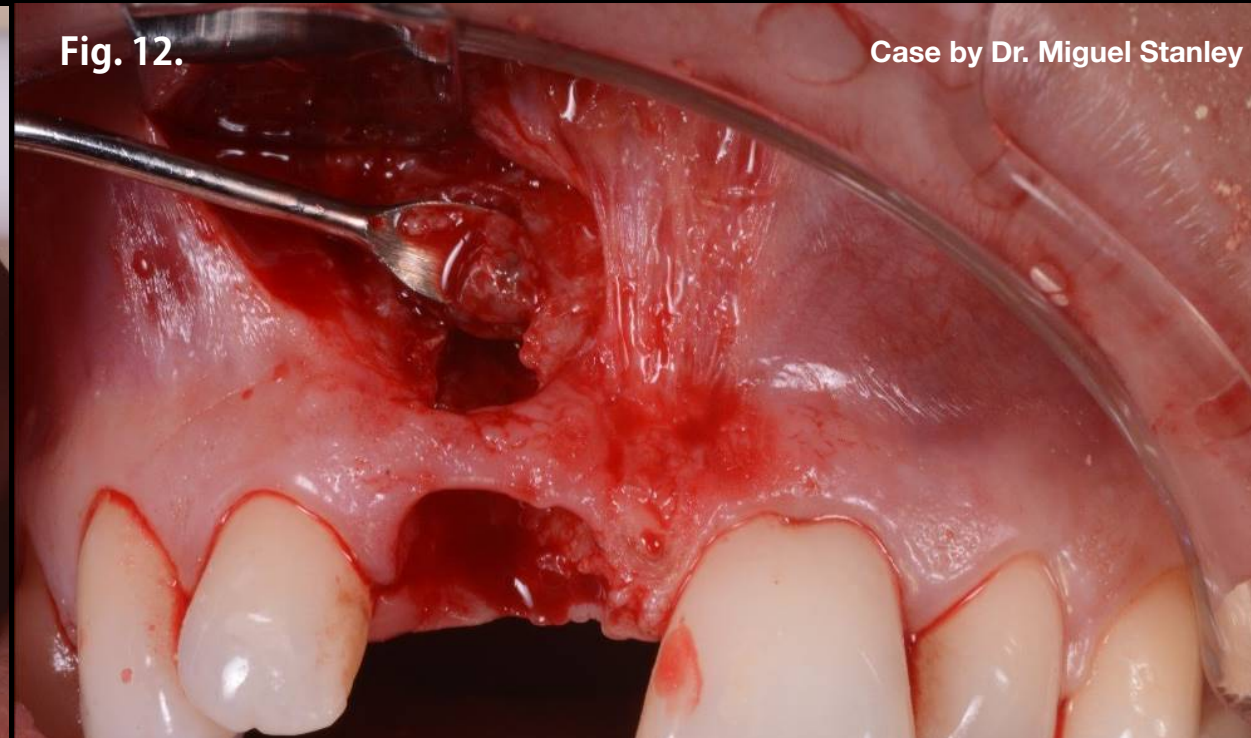


Fig. 13.



Fig. 11. Loss of buccal plate and vertical root fracture.
Surgery performed with aesthetic buccal flap.
Fig. 12. Cyst removal using curettes
Fig. 13. Cyst removed

Intra-oral clinical photos of the surgical procedure 17/11/2015 (Fig. 11-13)

Fig. 14.

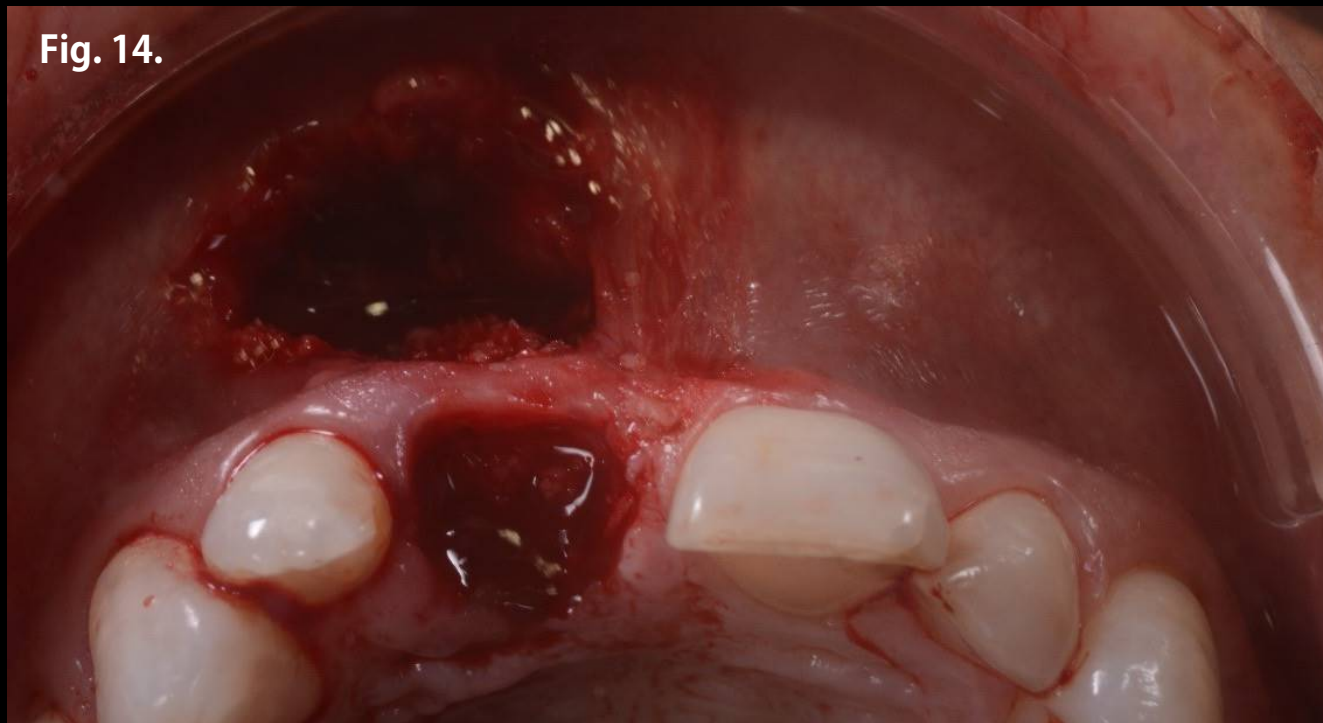


Fig. 15.

Case by Dr. Miguel Stanley



Fig. 16.



Fig. 14. Occlusal view during surgery

Fig. 15. Bone graft (Gen-os, Osteobiol) mixed with i-PRF

Fig. 16. Application of bone graft

Application of bone graft (Gen-os, Osteobiol) for ridge augmentation mixed with i-PRF (Fig. 14 and 16)



Fig. 17.



Fig. 18.



Fig. 19.

Case by Dr. Miguel Stanley

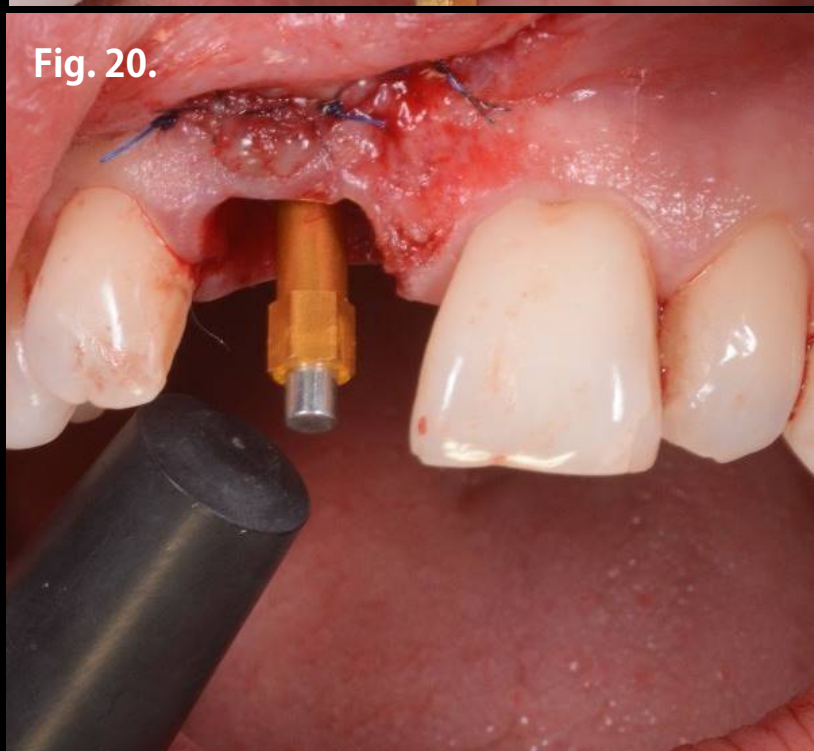


Fig. 20.

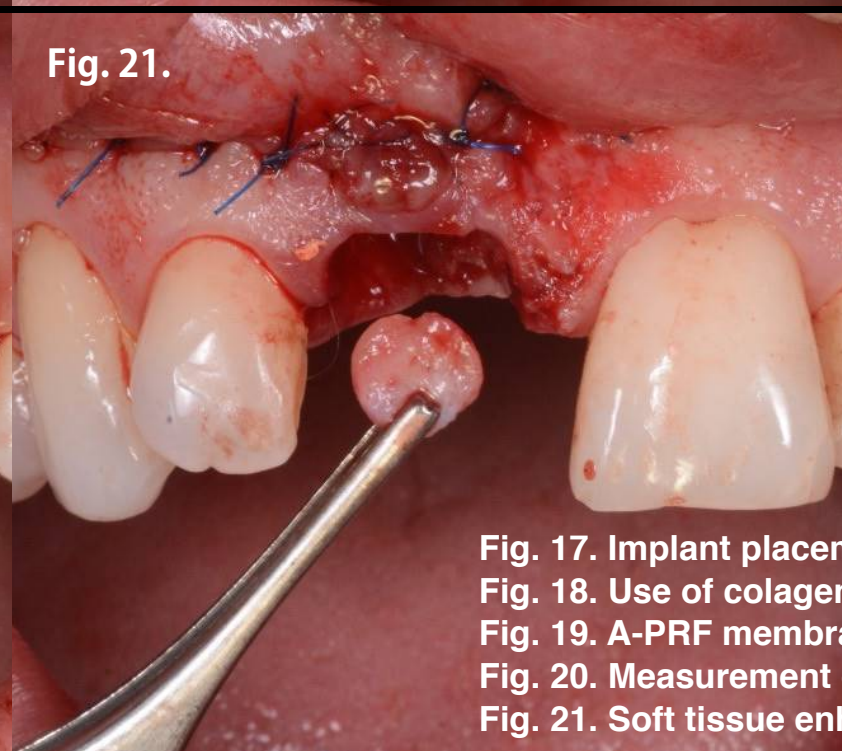


Fig. 21.

Fig. 17. Implant placement
Fig. 18. Use of collagen resorbable membrane (Evolution, Osteobiol) and
Fig. 19. A-PRF membrane
Fig. 20. Measurement of the Implant stability with Mega ISQ (ISQ 75)
Fig. 21. Soft tissue enhancement with gingival graft from the tuberosity

AnyRidge 4 Ø/ L 15 mm T-50 implant placement using AnyRidge surgical Kit (Fig. 17-21)

Fig. 22.



Fig. 23.



Fig. 24.

Case by Dr. Miguel Stanley



Fig. 25.



Fig. 26.

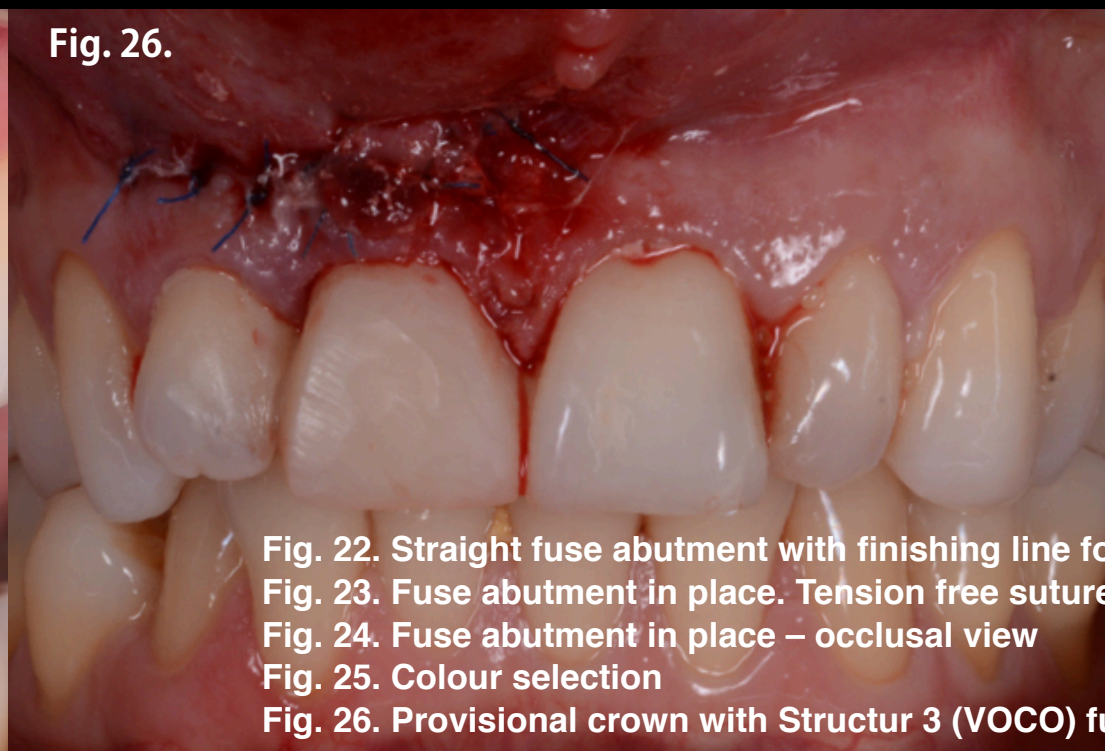


Fig. 22. Straight fuse abutment with finishing line for soft tissue modulation

Fig. 23. Fuse abutment in place. Tension free suture

Fig. 24. Fuse abutment in place – occlusal view

Fig. 25. Colour selection

Fig. 26. Provisional crown with Structur 3 (VOCO) fused to the adjacent teeth

Fuse abutment placement for temporary crown 17/11/2015 (Fig. 22-26)

Fig. 27.

Case by Dr. Miguel Stanley



Fig. 27. Immediate post-op panoramic x-ray 17/11/2015

Fig. 28.



Fig. 28. Follow-up one week after surgery: Intra-oral photo

Fig. 29. Impression coping x-ray: 2 months after surgery (16/1/2016)

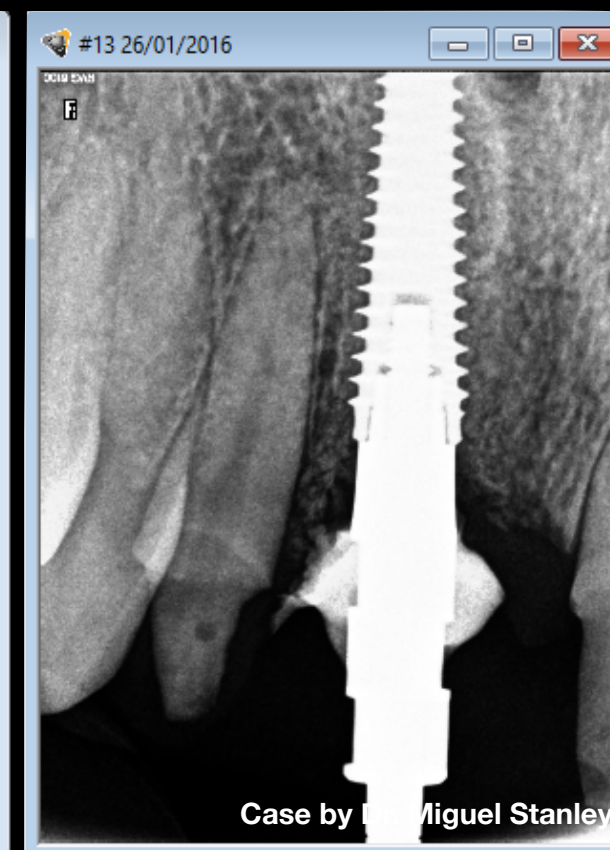
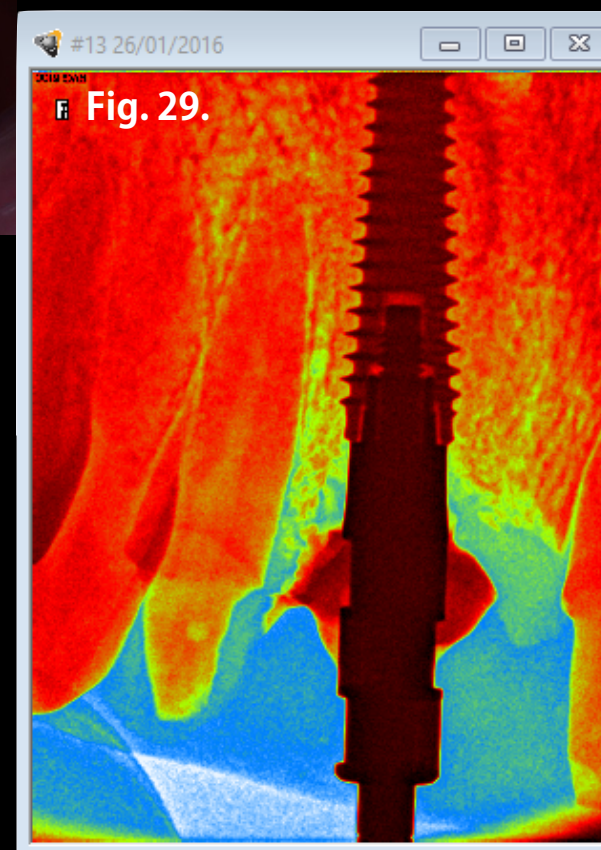


Fig. 30.



Fig. 31.

Case by Dr. Miguel Stanley

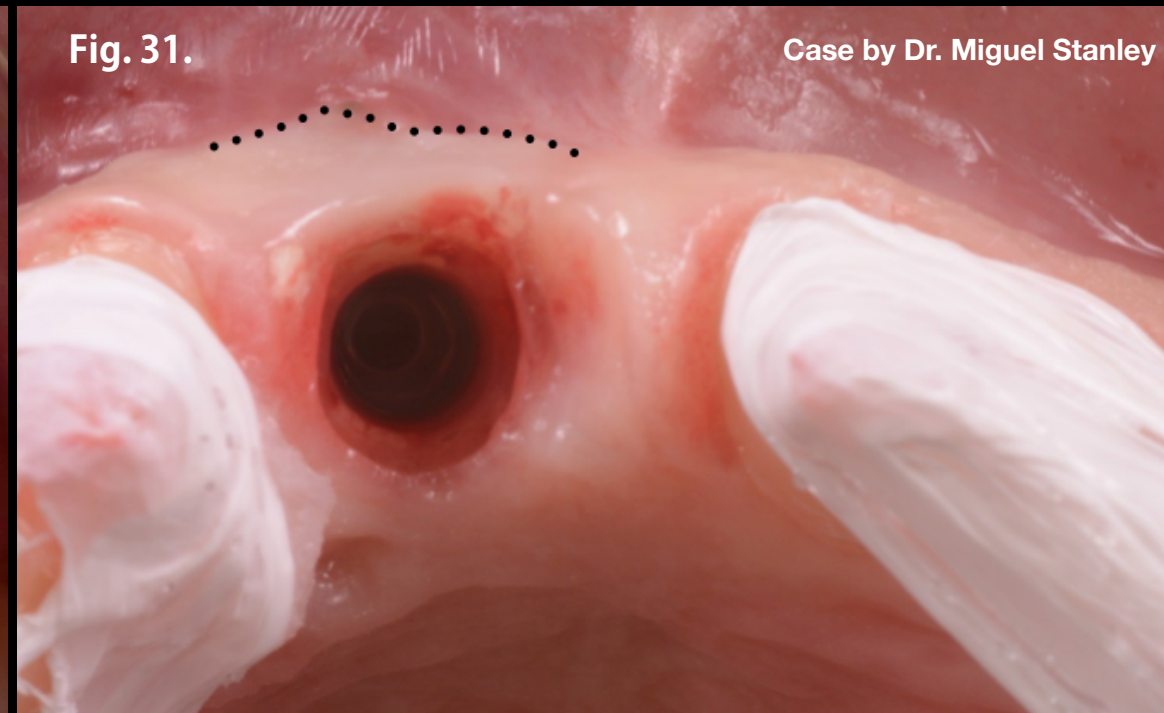


Fig. 32.



Fig. 33.



Fig. 32. First try-in – Frontal view
Fig. 33. First try-in – Lateral view

Occlusal view of the volume of buccal tissues. 3 months after surgery Fig. 29 and 30 (Patient asked for cosmetic enhancement of the anterior region).

Fig. 34.



Fig. 35.

Case by Dr. Miguel Stanley



Fig. 36.



Fig. 34. Customized zirconia abutment

Fig. 35. Cementation dual-cured resin-based adhesive luting system
(Bifix QM, Voco)

Fig. 36. Intra-oral X-ray

Final Cementation

Fig. 37.



Fig. 38.

Case by Dr. Miguel Stanley



Fig. 39.

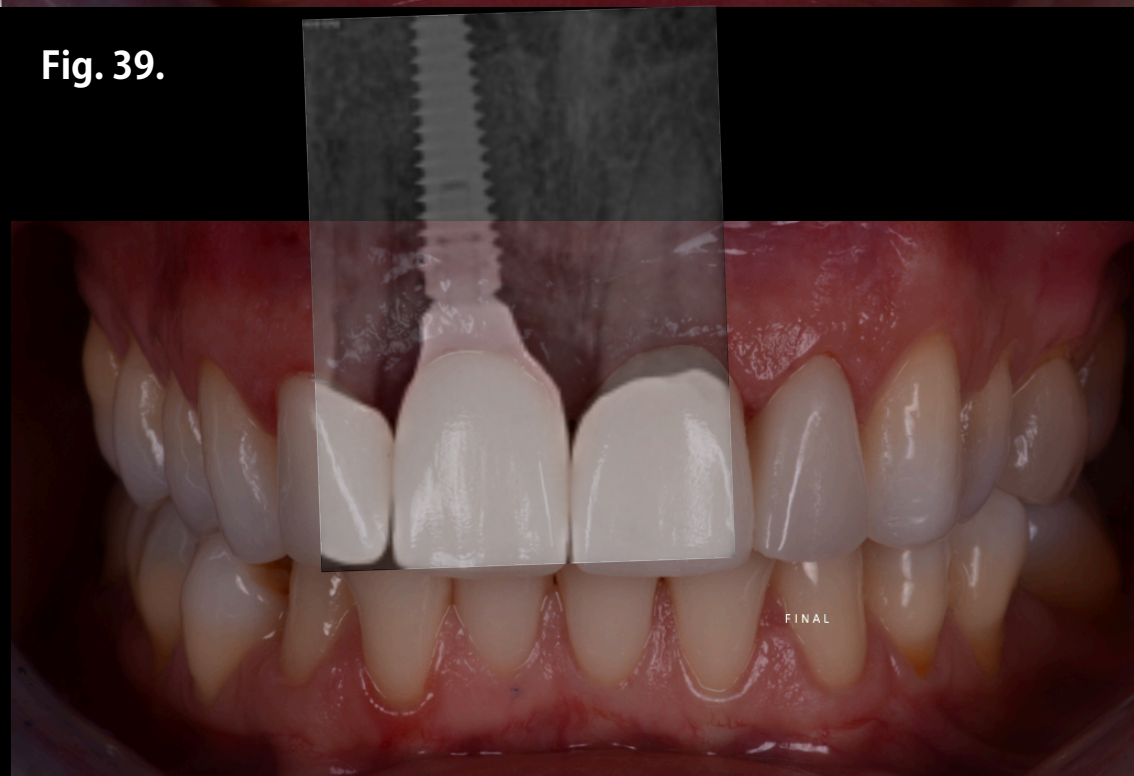


Fig. 40.



Fig. 37 and 38 Final crowns in milled zirconia and stratified lithium disilicate ceramic

Fig. 39. Intra-oral x-ray

Fig. 40. Final Panoramic x-ray (18/7/2017)

Fig. 41.

Case by Dr. Miguel Stanley



Fig.41. Follow-up: 6 months after loading

Fig. 42.

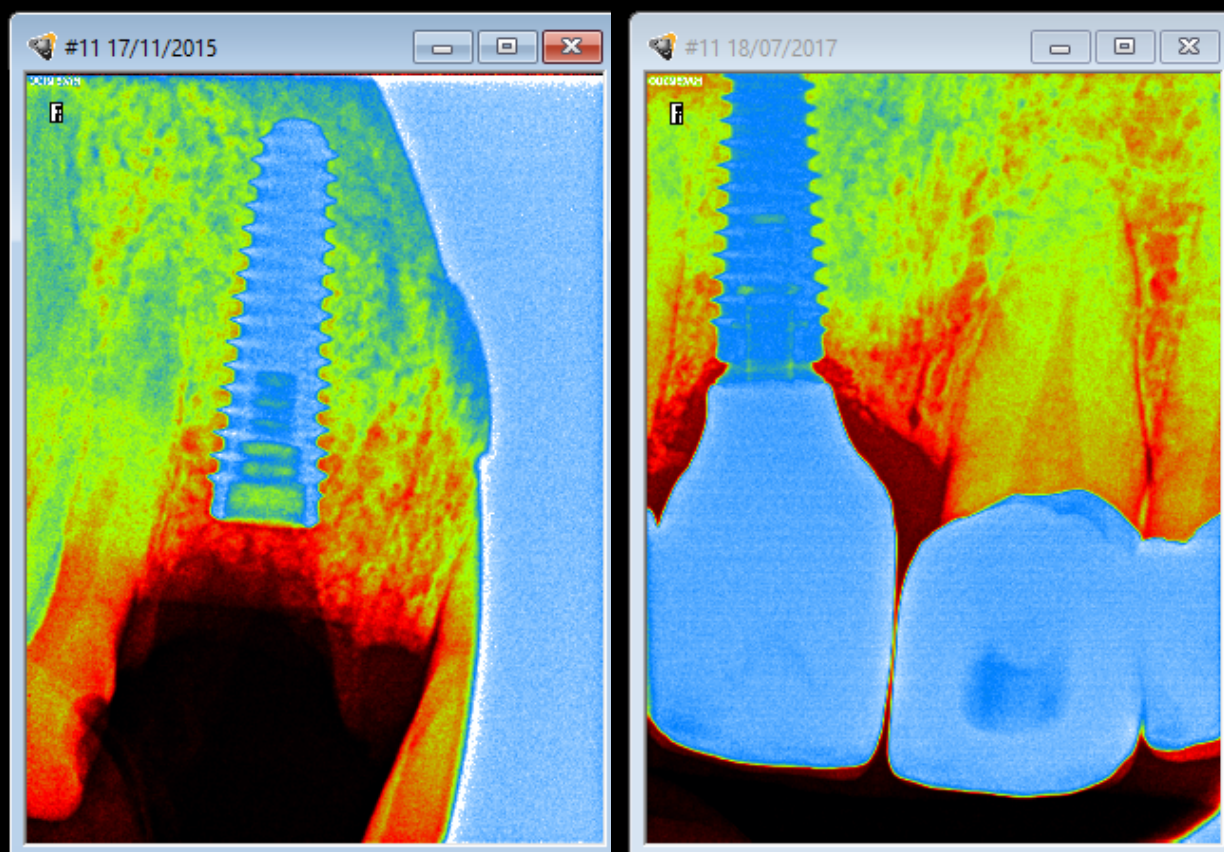


Fig. 43.

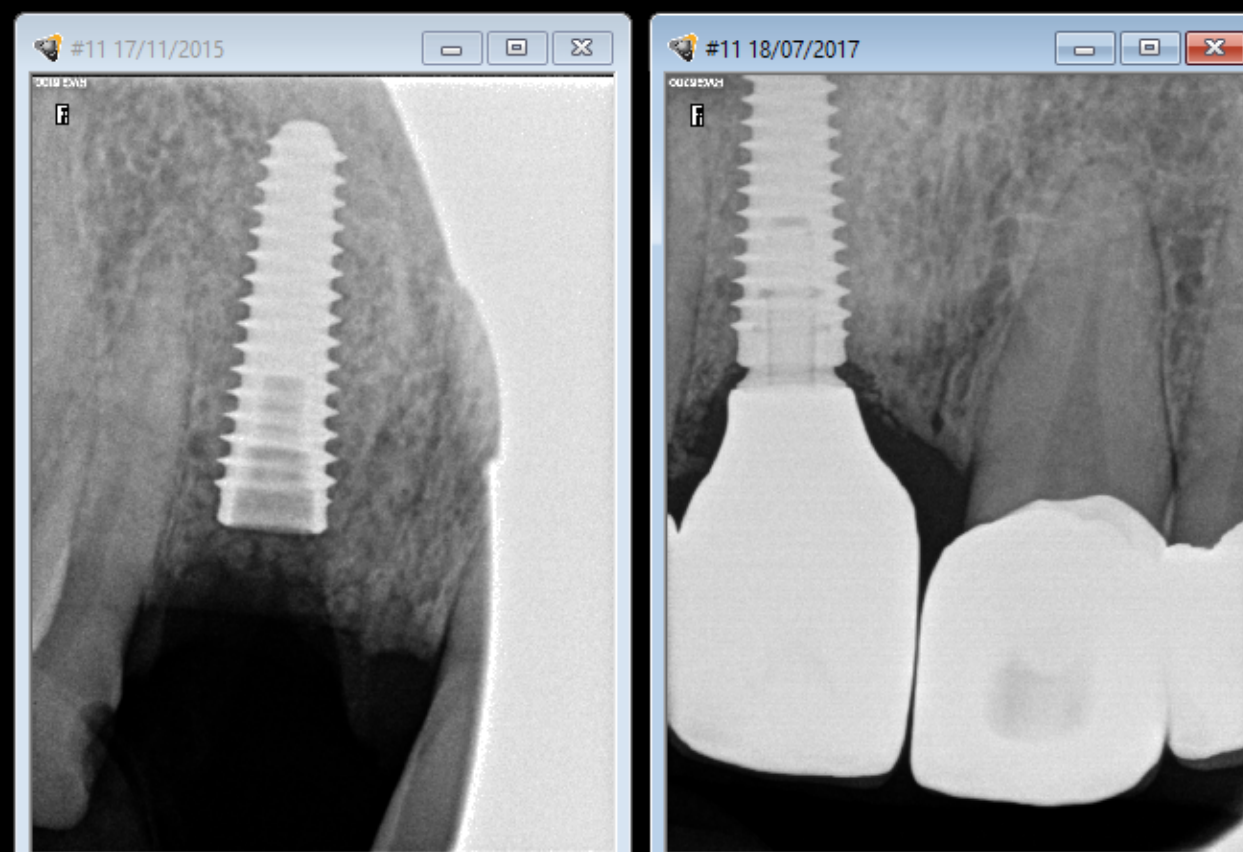


Fig. 42 and 43. Follow-up x-rays: 2 Years after loading (17/11/2015--18/7/2017)

Fig. 44.

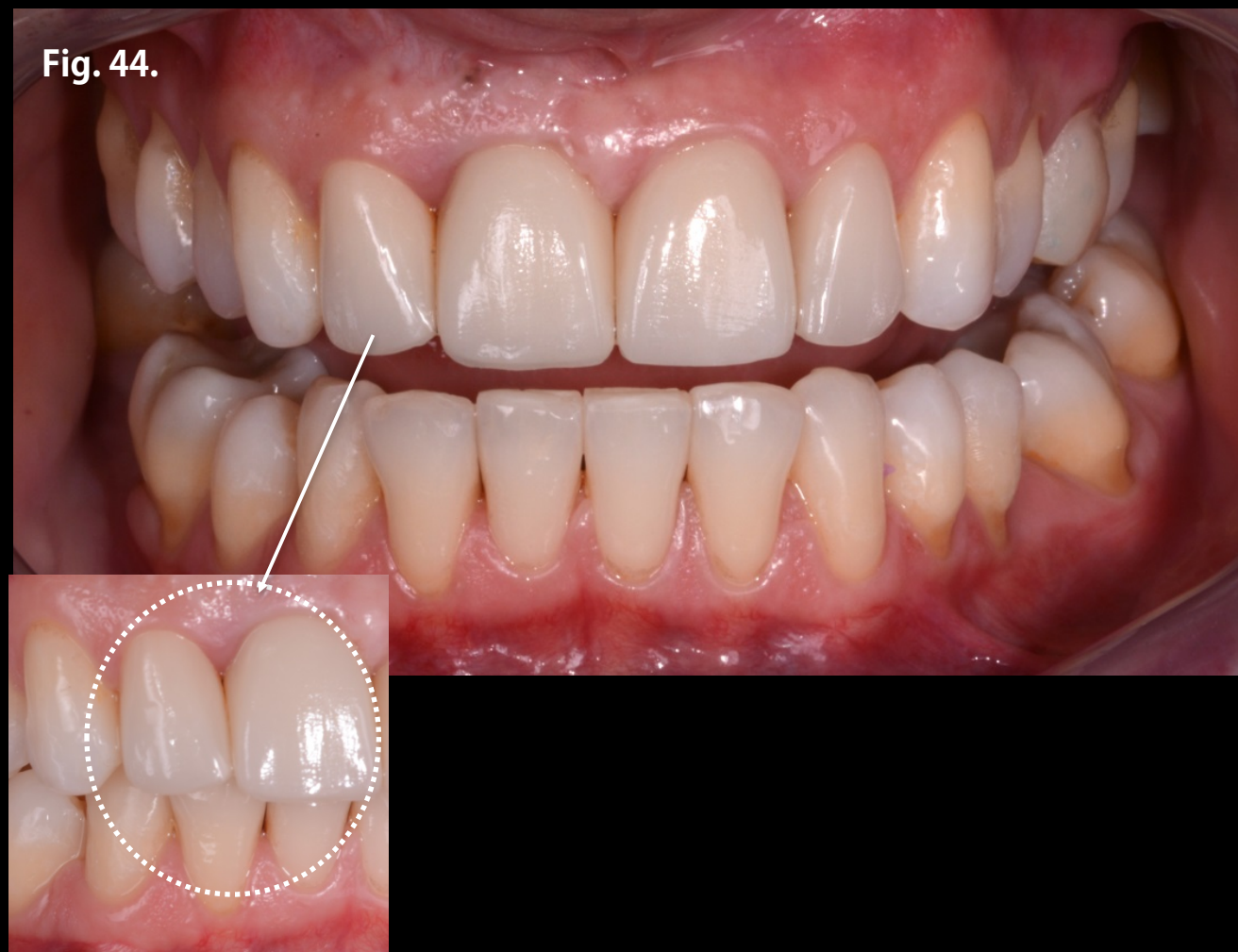


Fig. 45.

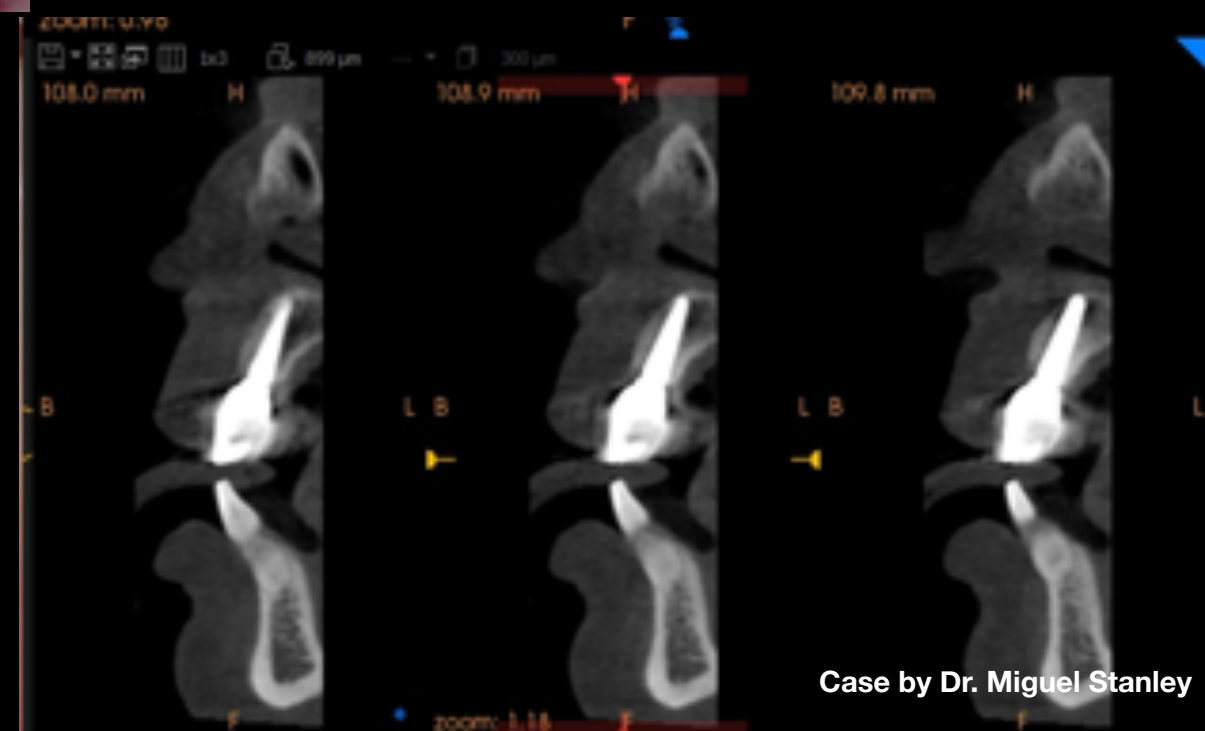


Fig. 44. Follow-up: 2 years after loading (18/7/2017)

Fig. 45. Follow-up: CBCT 2 years after loading (18/7/2017)

