



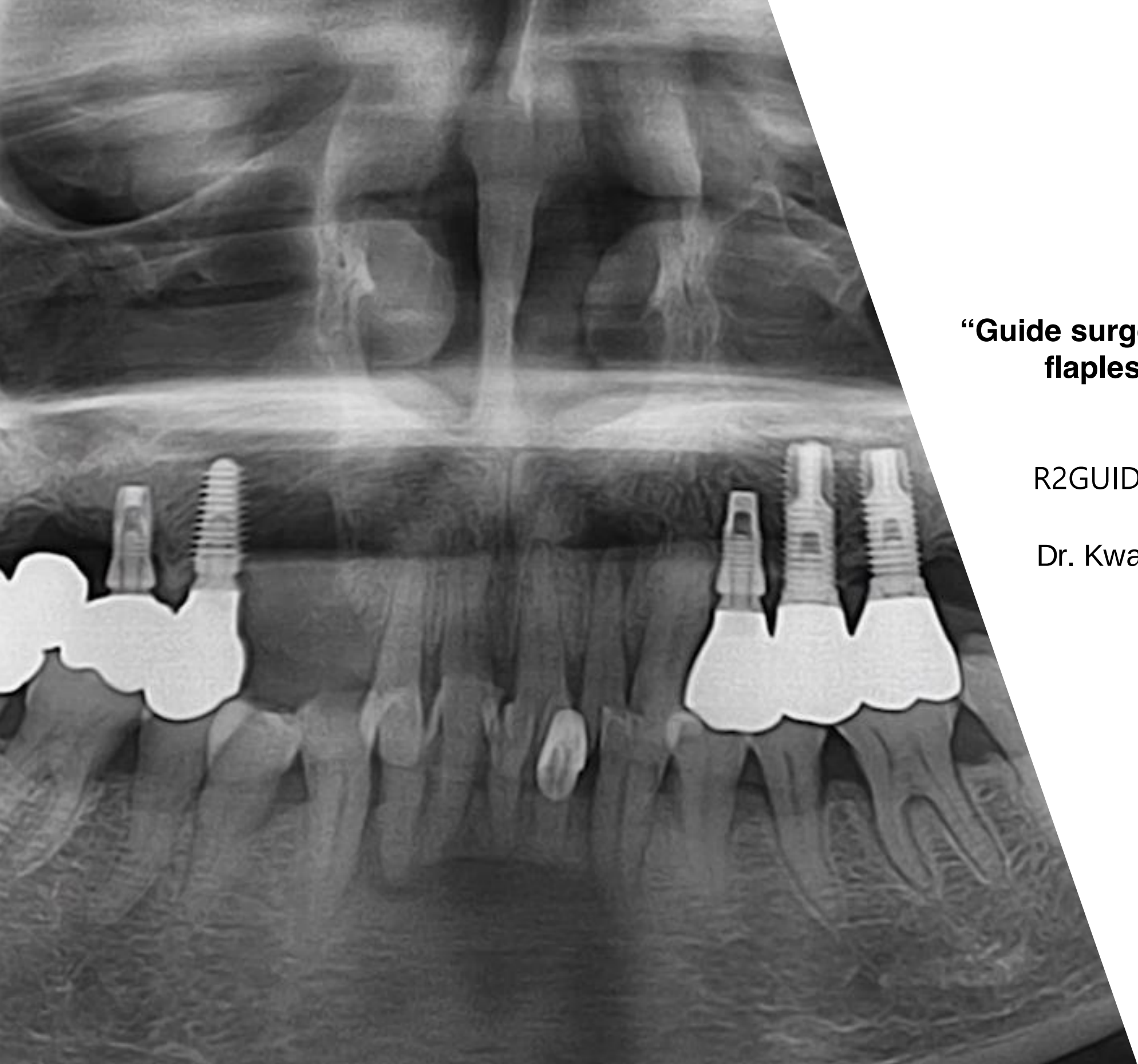

**“Guide surgery always means  
flapless surgery?”**

R2GUIDE surgery with GBR

Dr. Kwang Bum Park

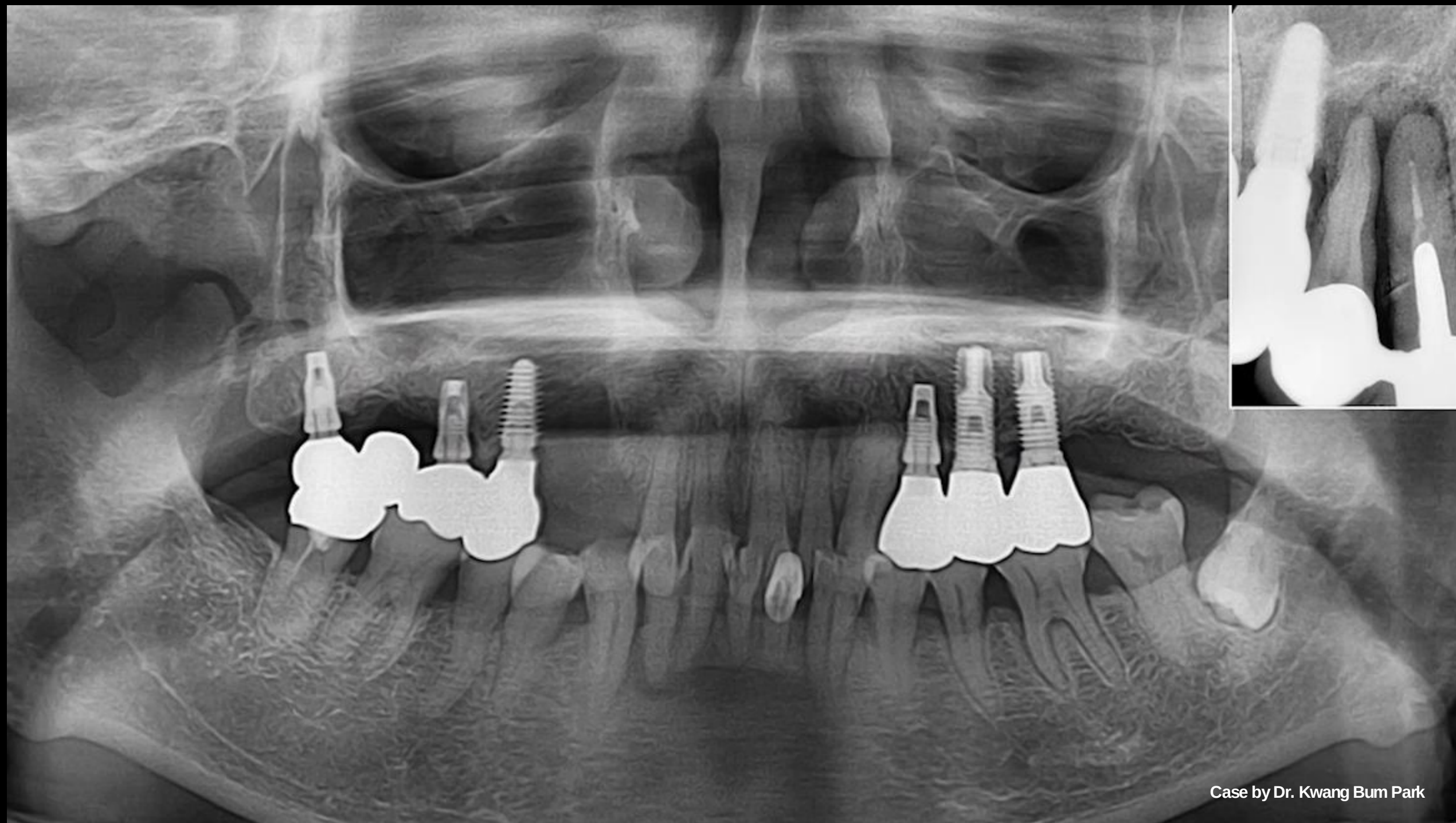


Case by Dr. Kwang Bum Park

**Clinical Case:** # 13,14 Implant placement Case, Surgery using R2GUIDE



Case by Dr. Kwang Bum Park

**Clinical Case: # 13,14 Implant placement Case**

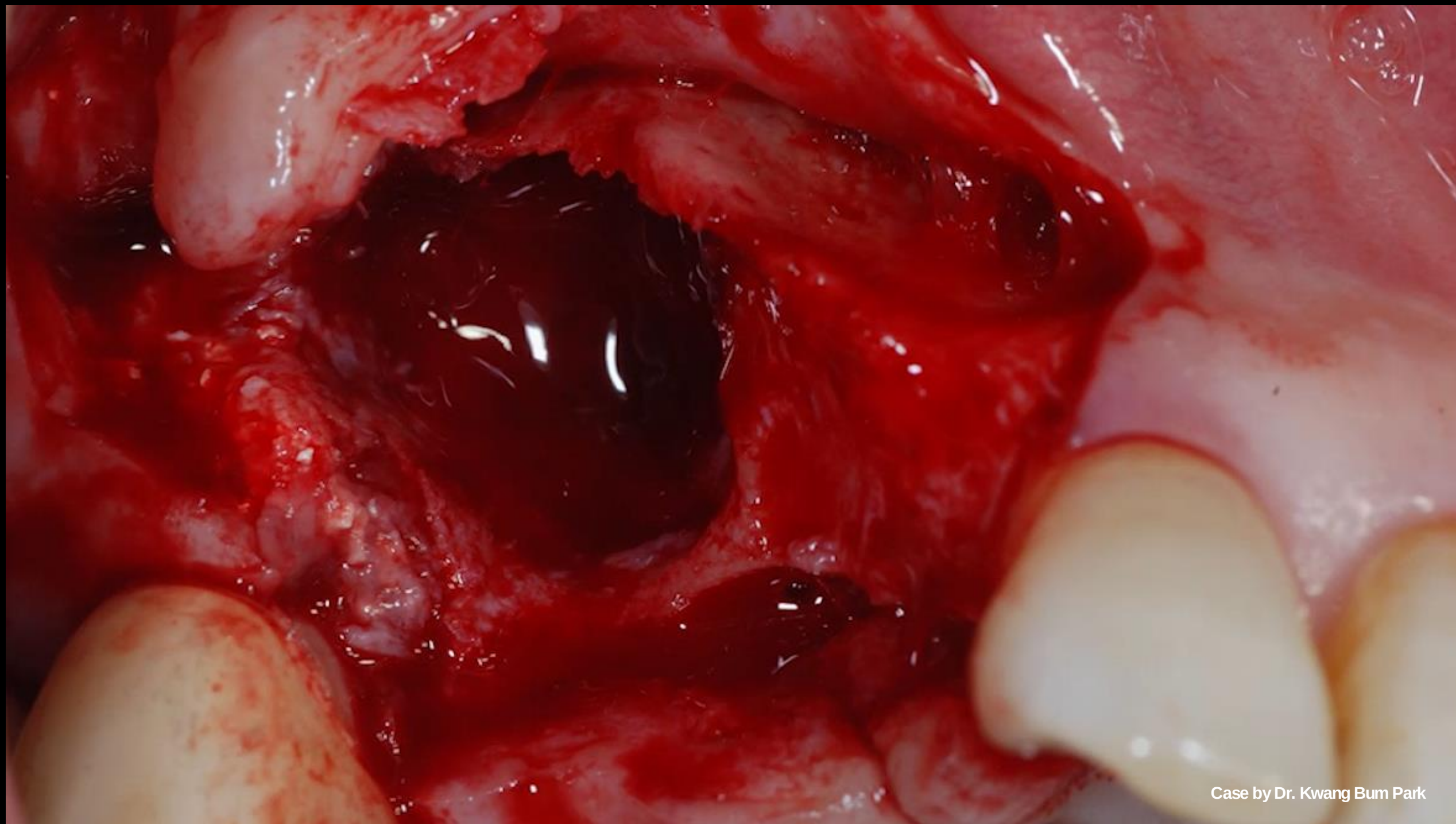


Case by Dr. Kwang Bum Park

**Clinical Case:** Even if you use the guide, open flap as needed.




Case by Dr. Kwang Bum Park

**Clinical Case : # 13,14 Implant placement Case**

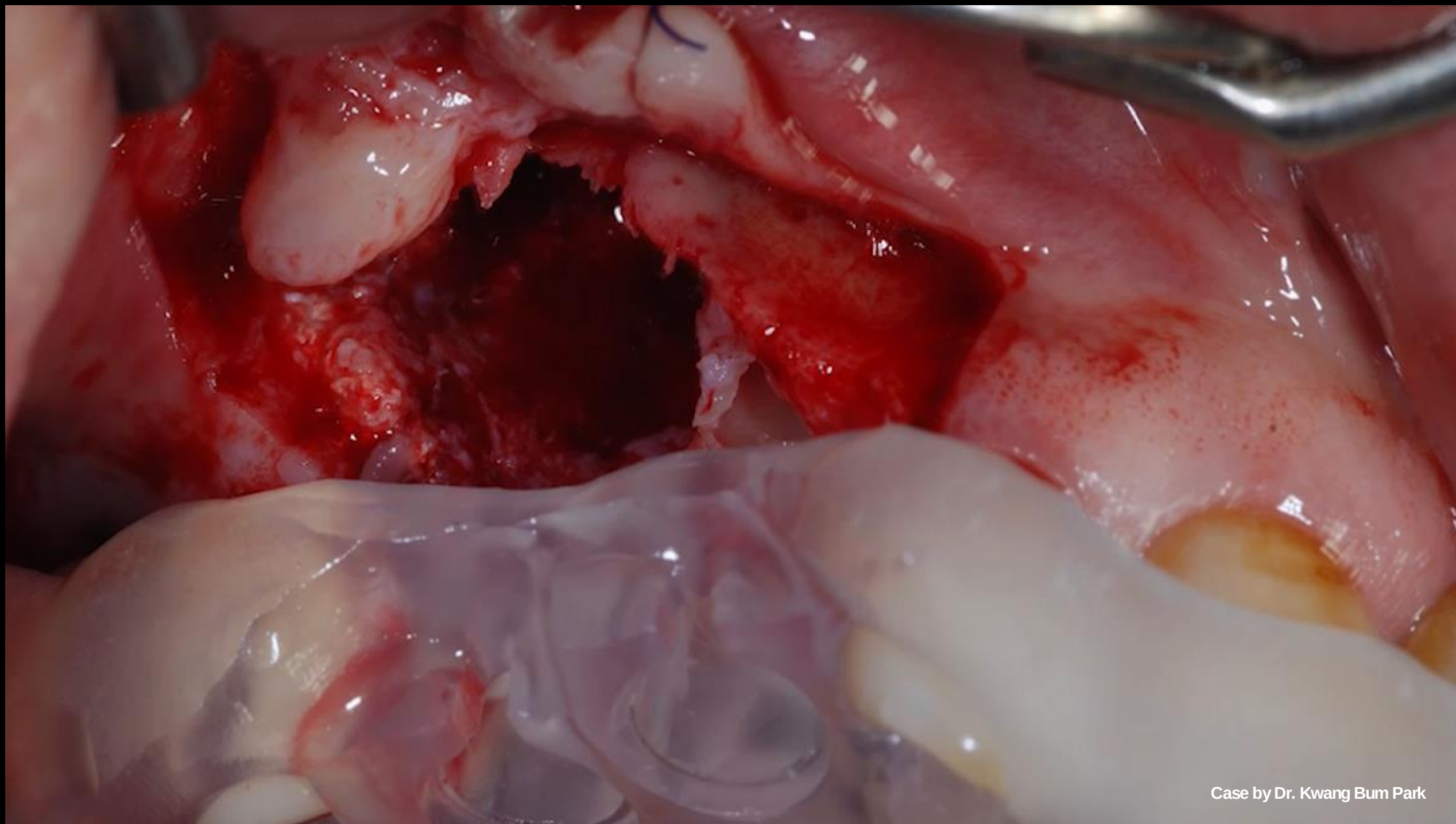


Case by Dr. Kwang Bum Park

**Clinical Case :** Put the R2GATE Guide on the teeth

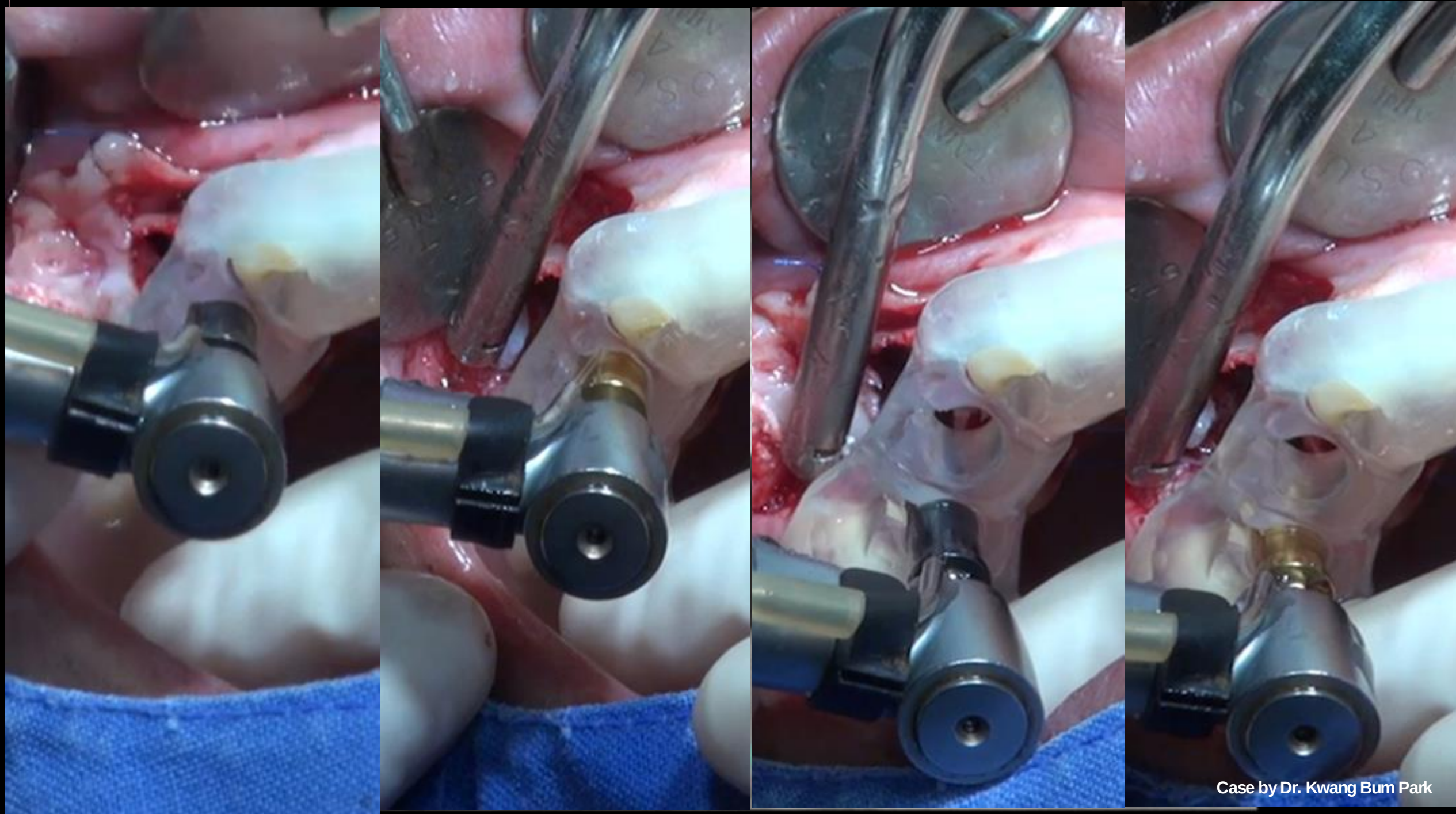

Case by Dr. Kwang Bum Park

**Clinical Case :** Drilling according to drilling sequence using bone quality through digital eye

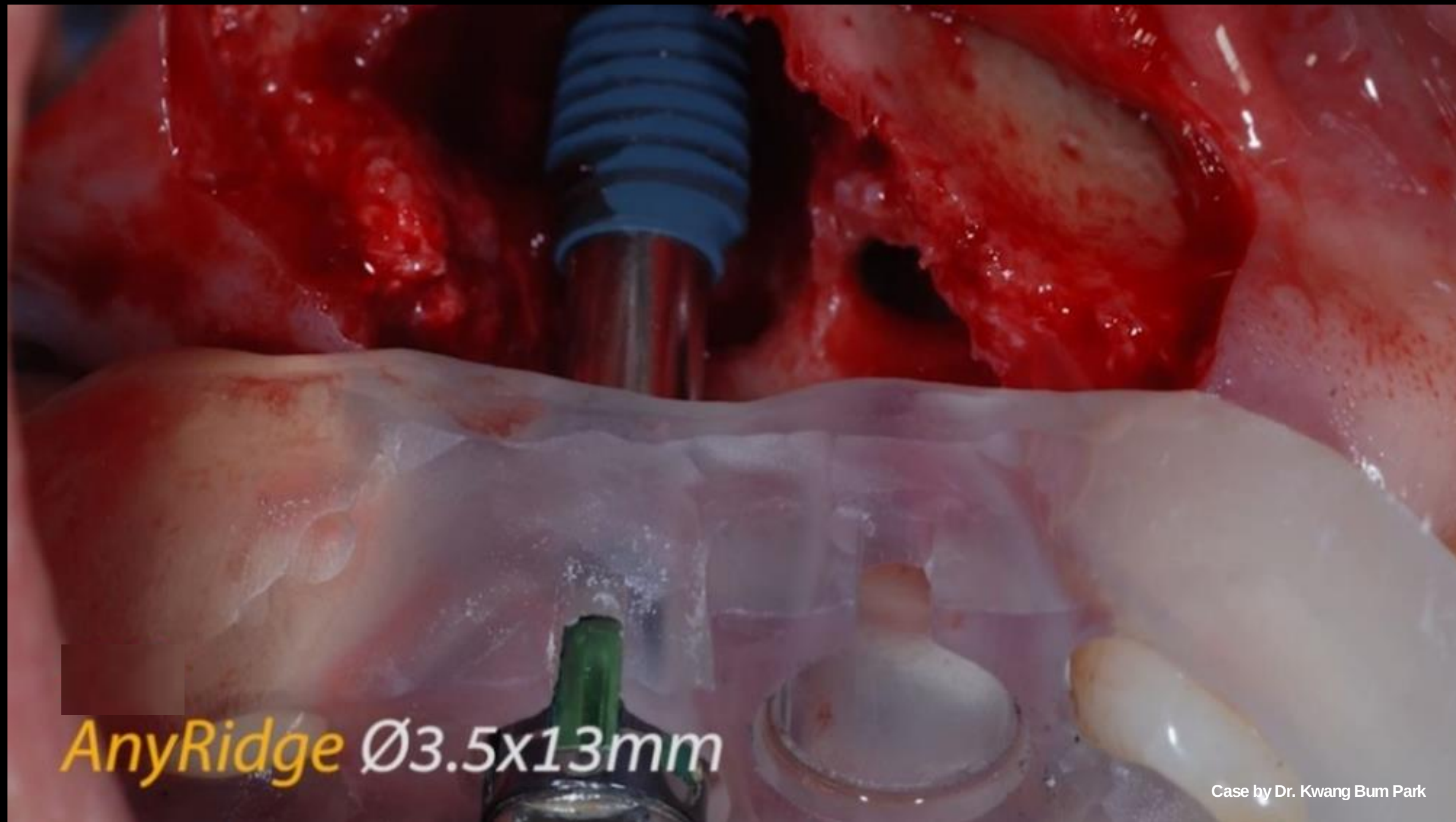


**Clinical Case :** #14, AnyRidge 3.5x13mm Implant placement through the guide



#13  
*AnyRidge* Ø3.5x13mm

Case by Dr. Kwang Bum Park

**Clinical Case :** #13, AnyRidge 3.5x13mm Implant placement through the guide

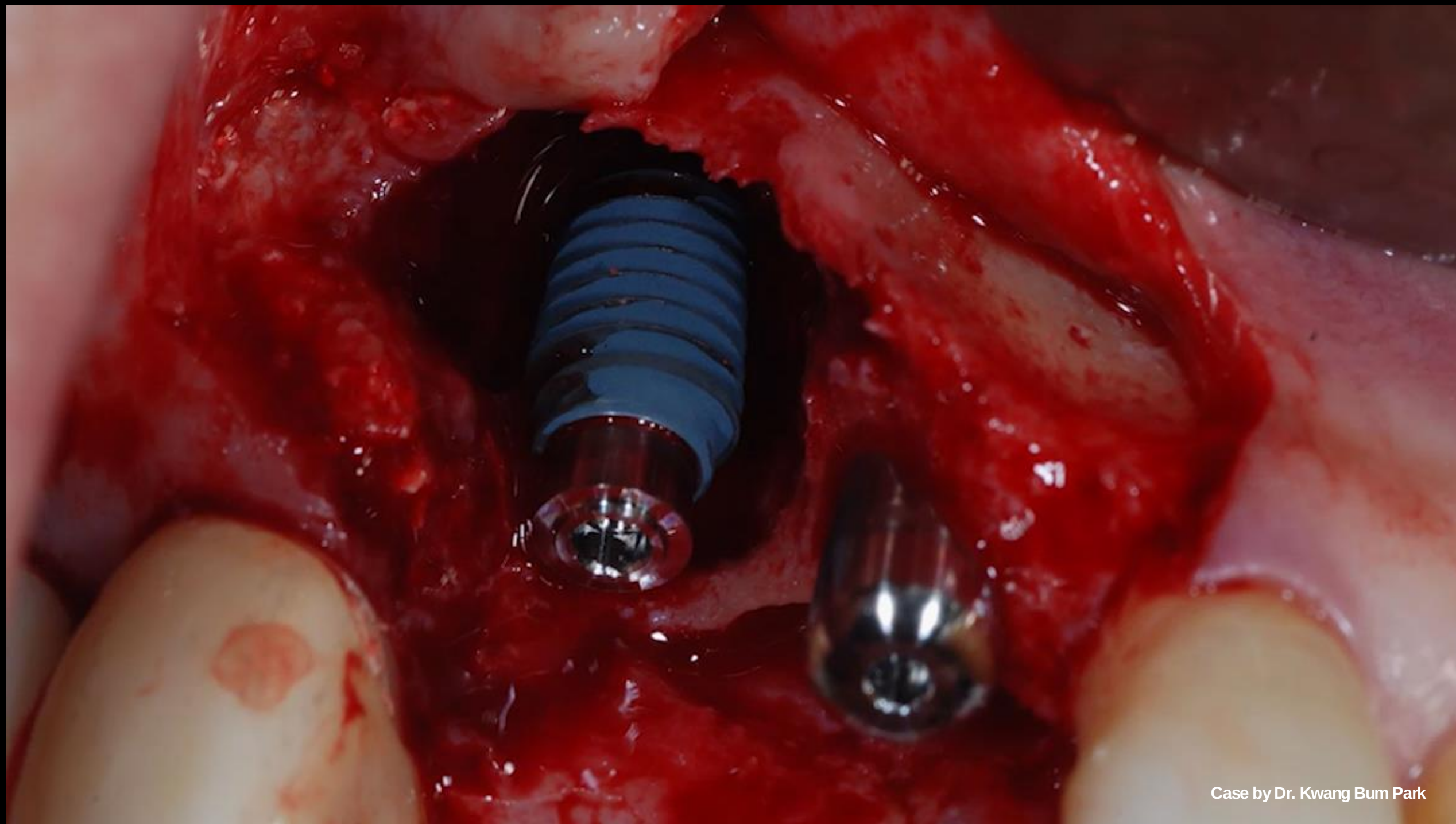

Case by Dr. Kwang Bum Park


**Clinical Case :** Check the implant position



*Use Mega-Oss*

Case by Dr. Kwang Bum Park

**Clinical Case :** GBR using Mega-Oss

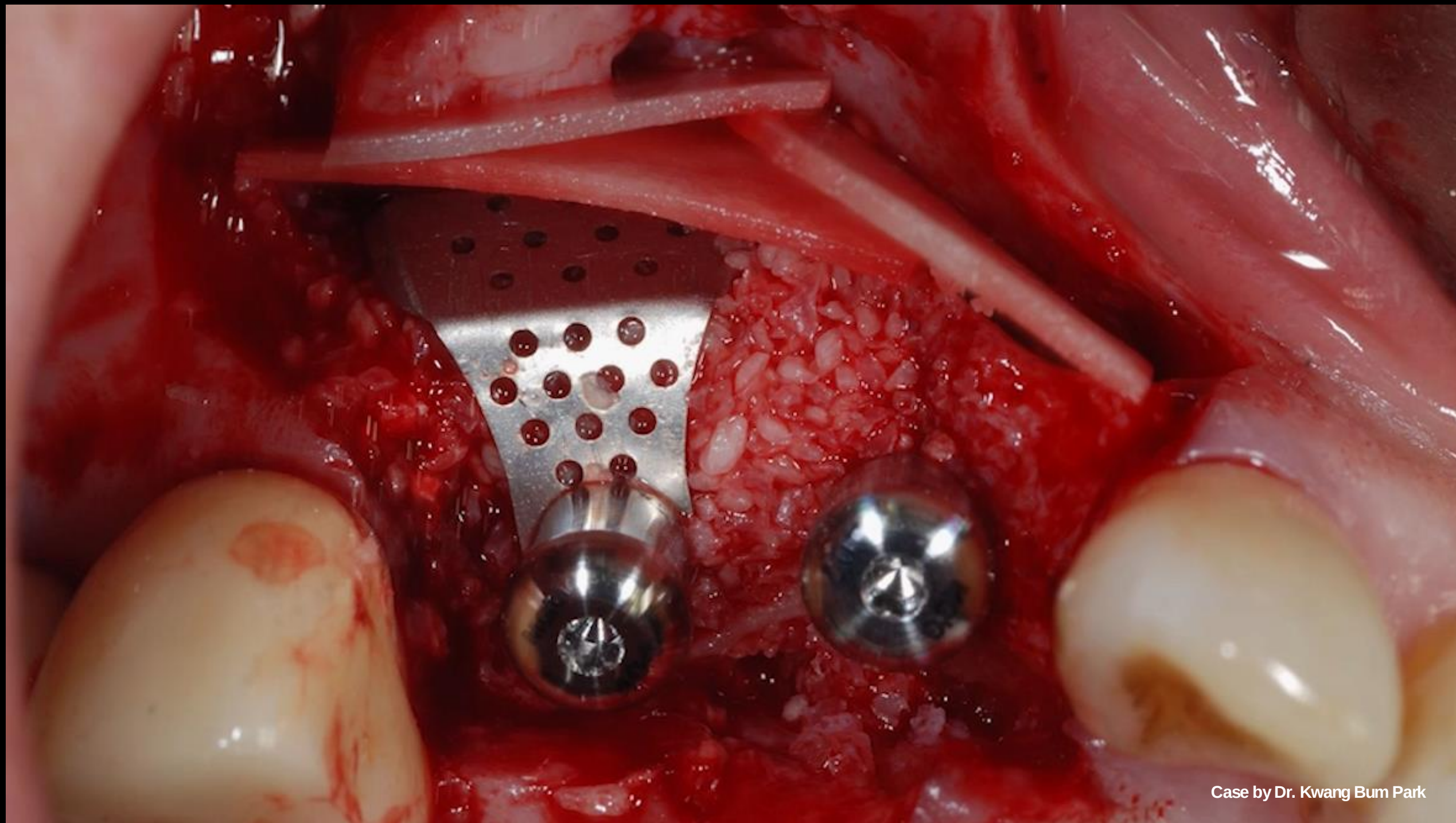

Case by Dr. Kwang Bum Park


**Clinical Case : GBR finished**



*Use i-Gen membrane*

Case by Dr. Kwang Bum Park

**Clinical Case :** Use i-Gen membrane (Titanium mash)

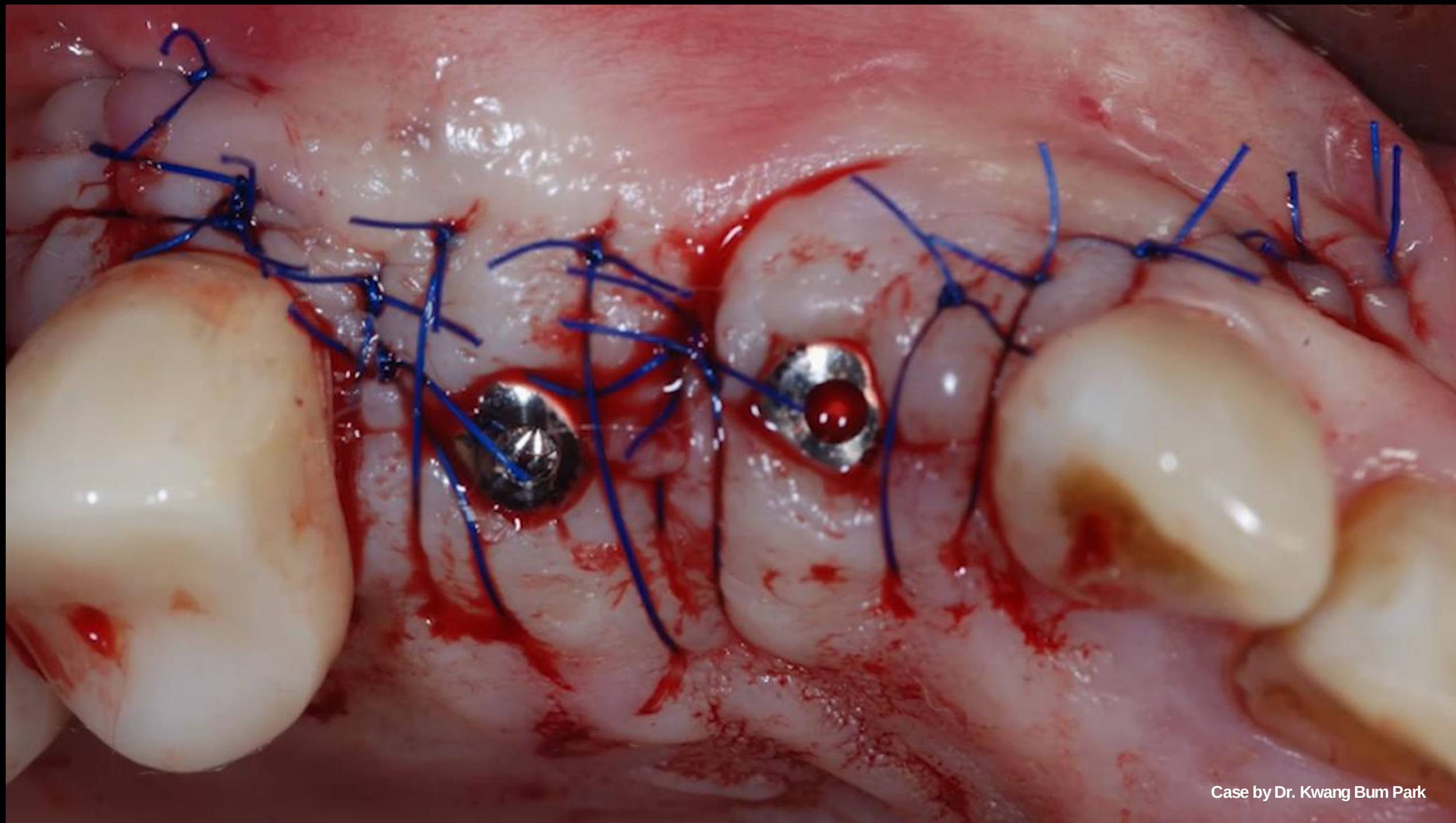

Case by Dr. Kwang Bum Park

**Clinical Case :** Suture finished (occlusal surface)

