


HOW to use *MEGA ISQ™ II*

1. Activate the MEGA ISQ II by picking it up then place a barrier sleeve over the instrument.

The instrument will start-up and be ready for measurement in the BL (Buccal – Lingual) direction.

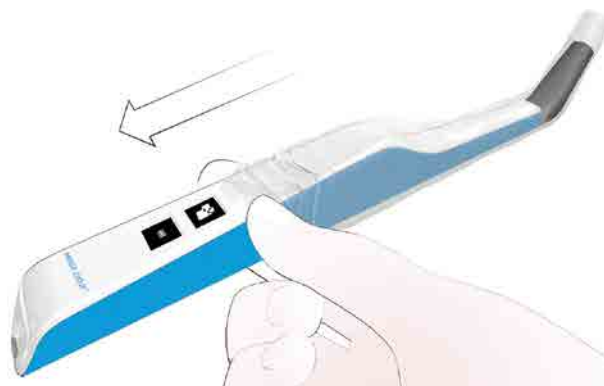

2. Attach the SmartPeg to the implant or abutment by screwing the SmartPeg Mount using finger force of approximately 4-6 Ncm.

3. Hold the instrument tip close (3-5mm) to the top of the SmartPeg without touching it. A sound indicates when measurement has started and the value will be shown in the upper display combined with a colored light indication below the instrument tip.

Bring the instrument out of the mouth to clearly read the ISQ value and the colored indication. The measured ISQ values will be displayed in the upper display for a couple of seconds and then switch to indicate ready for measurement in mesial-distal direction.

