

►► How to use

1 . BonePen



Use suitable BonePen diameter -Drill tip must reach the alveolar bone.
Check the position of the alveolar bone to be cut. Adhere product's cuff side to tooth side and use it while the saline is running.
In order to avoid the BonePen spinning beyond the cutting point, use handpiece with both hands.

2 . BonePin



Select a BonePin that has same numbering as BonePen used.
Insert the bottom of the product as deep as possible into the insertion hole that you have made previously.
If the bone insertion is completed, visibly assess the gap between product's external diameter and adjacent tooth, also cutting height of the alveolar bone and occlusal surfaces during occlusion.

3 . GBR Pen




Decortification

Collected bone inside the cup
(8th drilling)

Bone graft


Suture

4 . BoneShaper

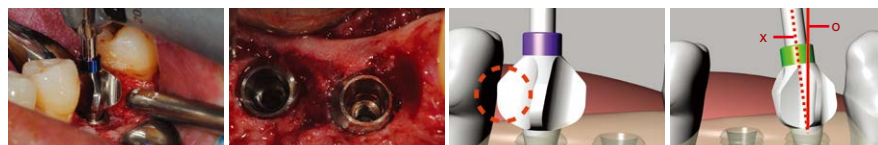

Select a BoneShaper that has same numbering as BonePin used.
Place BoneShaper onto the center of the fixture, and begin countersinking by pushing it downward.

5 . GingivalShaper

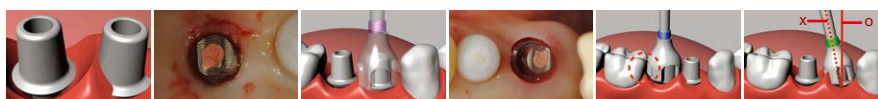

Select a GingivalShaper that is suitable for a diameter of the Abutment used.
Place GingivalShaper onto Abutment Shoulder, and use a begin countersinking using a hand-piece for 3~5 seconds.