




Immediate post-extraction insertion of implant & immediate loading


Prof. Giuseppe Luongo

Case by. Prof. Giuseppe Luongo

Case by. Prof. Giuseppe Luongo

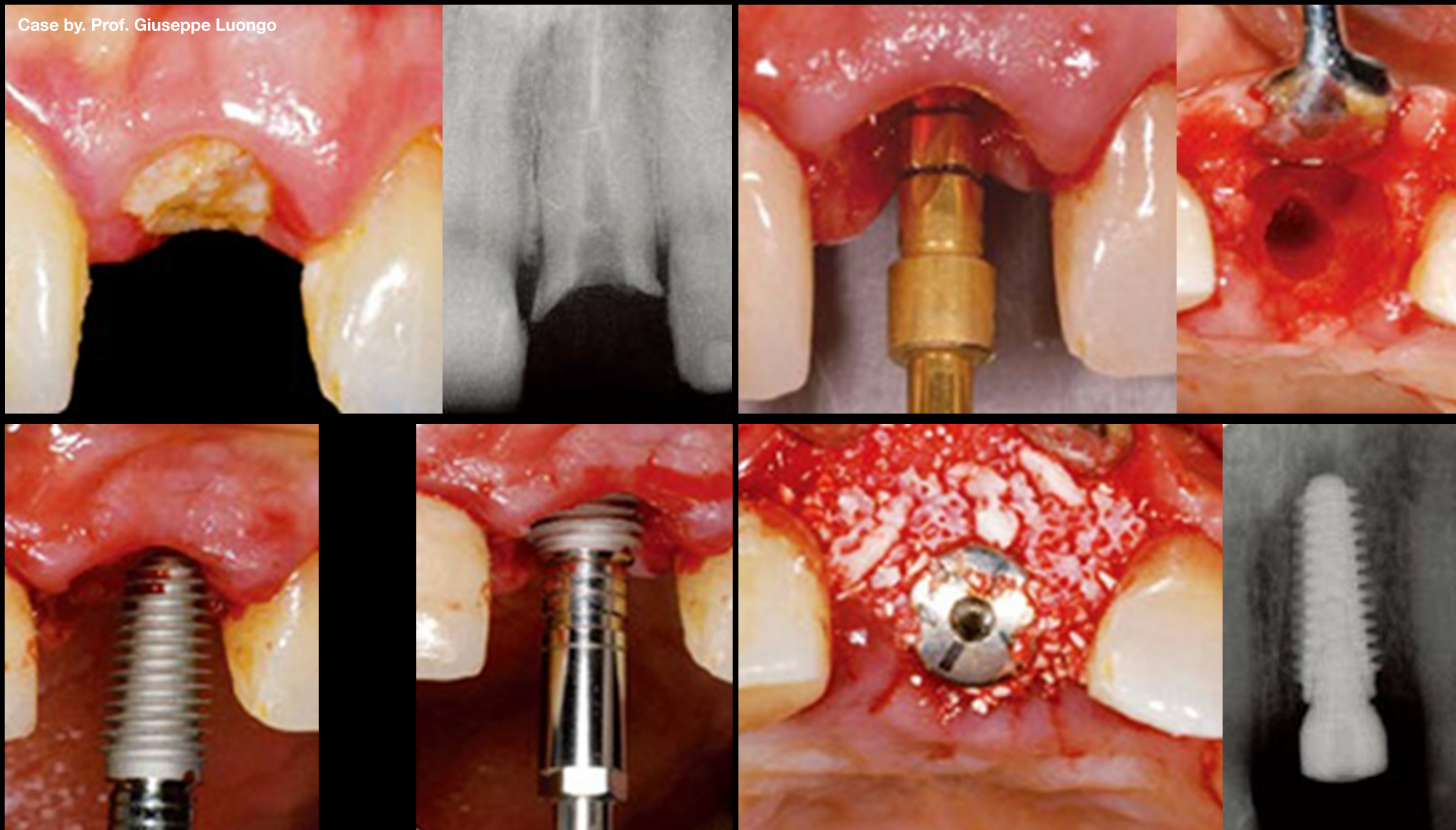

Fracture of #21 tooth / Good stability of hard & soft tissues indicates immediate post-extraction insertion of implant & immediate loading
To protect esthetic outcome, implant site was prepared via slightly palatal alveolar access
4.5×11mm AnyRidge was placed in prepared site
Biomaterial was added to vestibular aspect to improve stability of esthetic

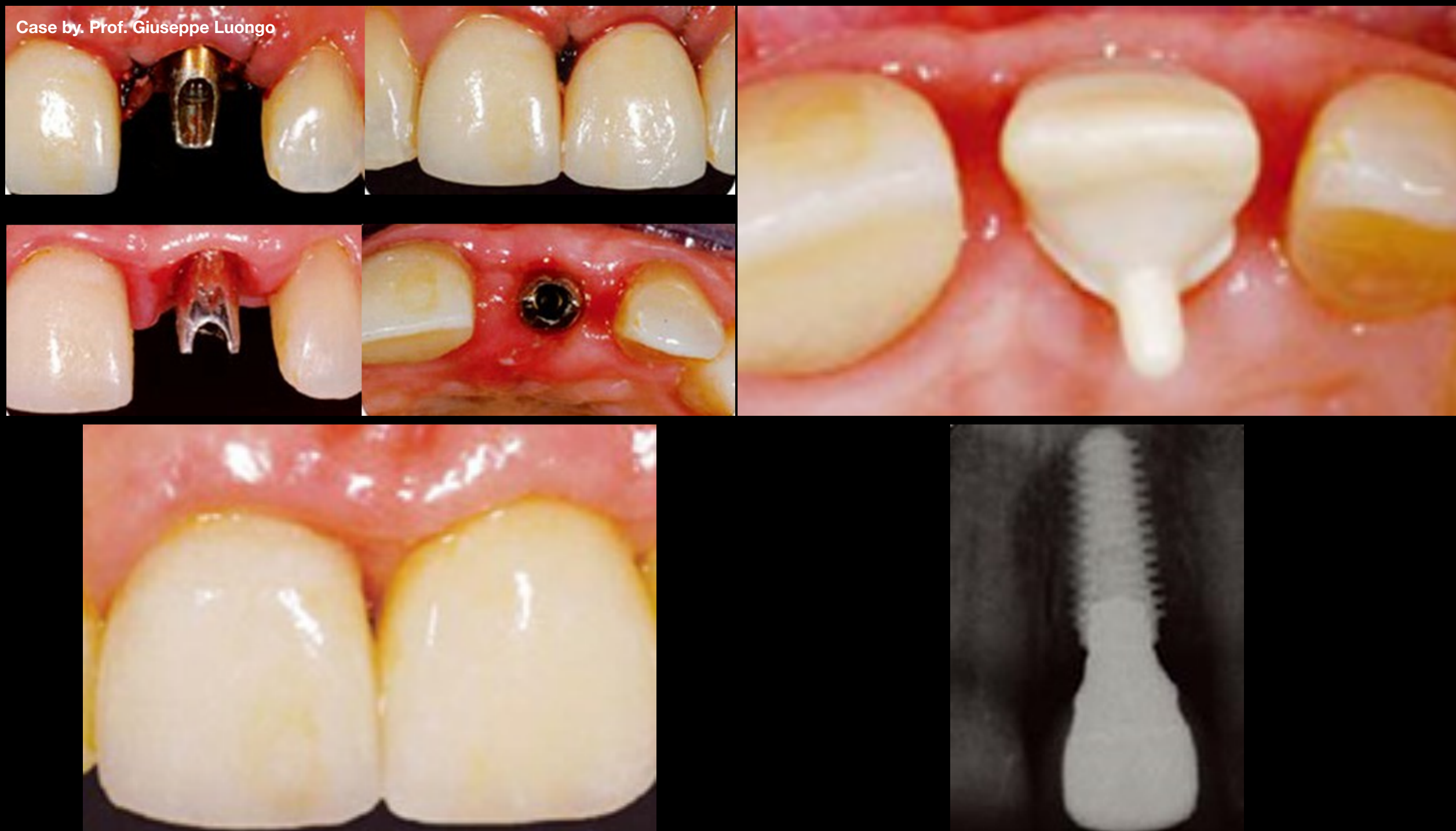

Temporary crown was placed immediately / After 12 wks, tissue was ready to proceed with final abutment & crown
Implant position harmonized with surrounding tissue
Photo of final crown at time of placement
Intraoral scan of final crown at time of placement

Case by. Prof. Giuseppe Luongo



Photo of final crown at 1-yr follow-up

Intraoral scan of final crown at 1-yr follow-up

Photo of final crown at 5-yr follow-up

Intraoral scan of final crown at 5-yr follow-up

