




Installation of dental implants in
complicated anatomic conditions using
crest lifting methods

Dr. Alexander Lysov

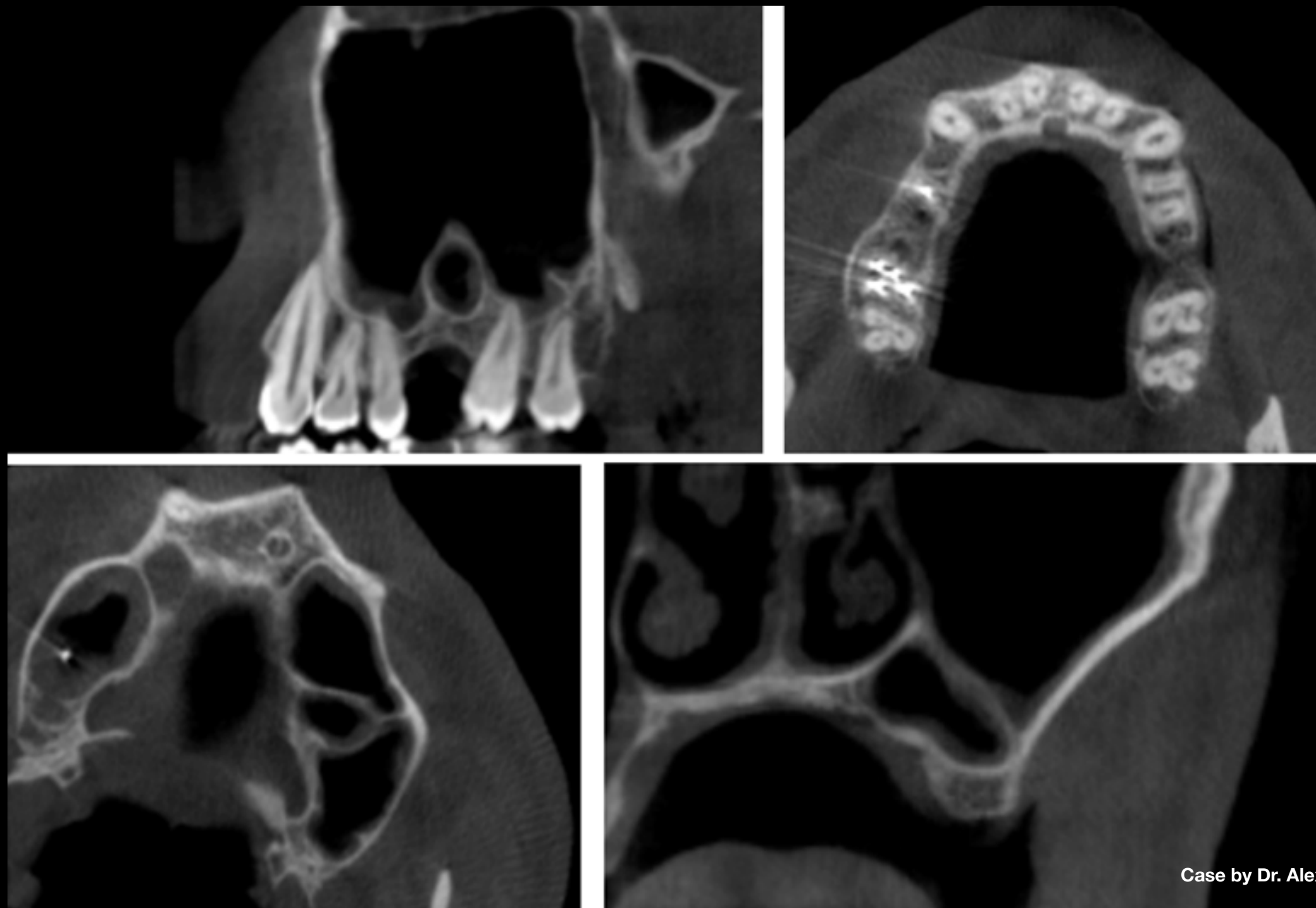


Case by Dr. Alexander Lysovov

Clinical Case: As alveolar ridge thickness of #26 tooth projection was 5mm, crest sinus lifting using MICA Kit was performed for best results





Case by Dr. Alexander Lysovov

Clinical Case: Penetration milling guide marks on trephine

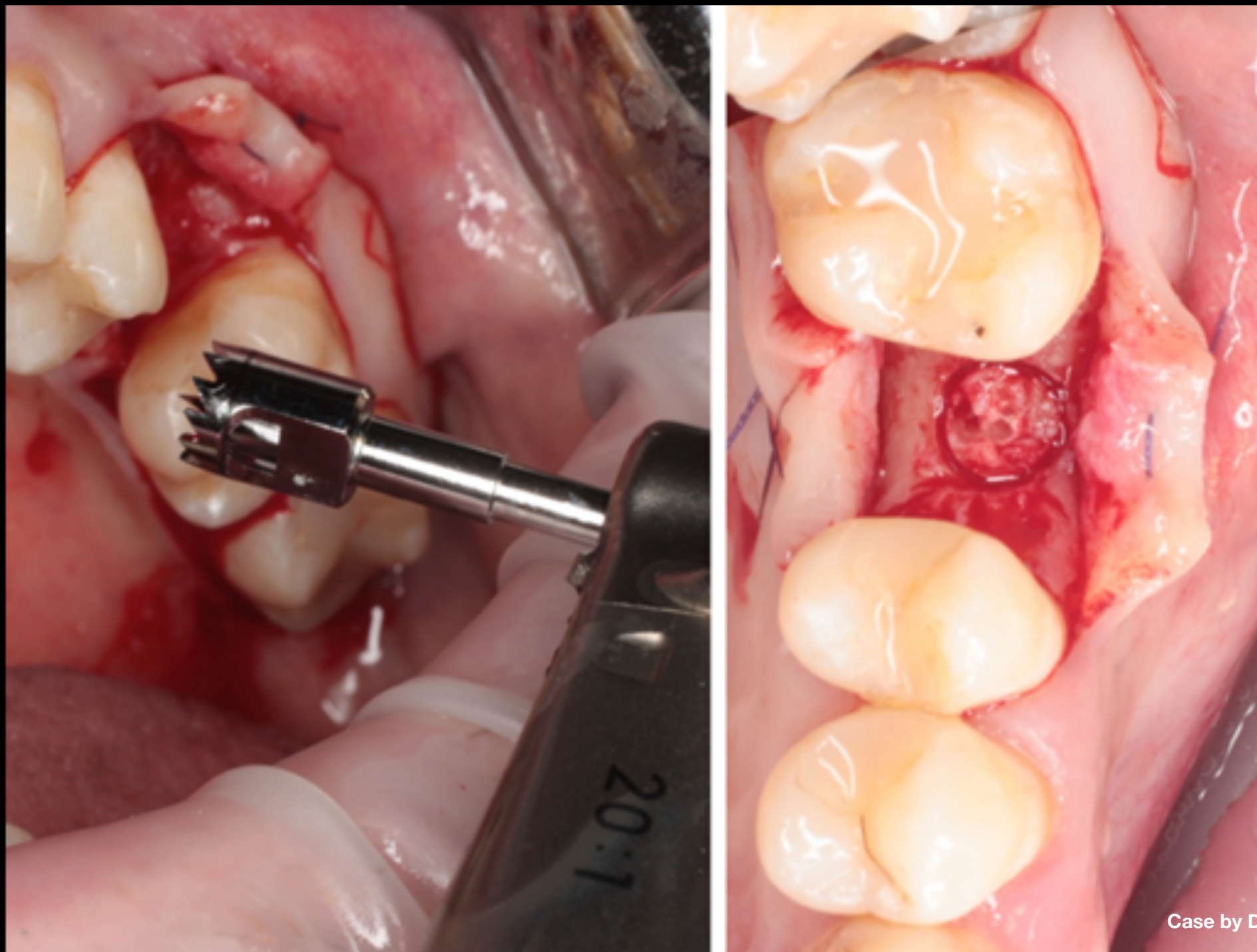

Case by Dr. Alexander Lysovov

Clinical Case: Drilling using point trephine bur (2mm per time until laser marking is reached)

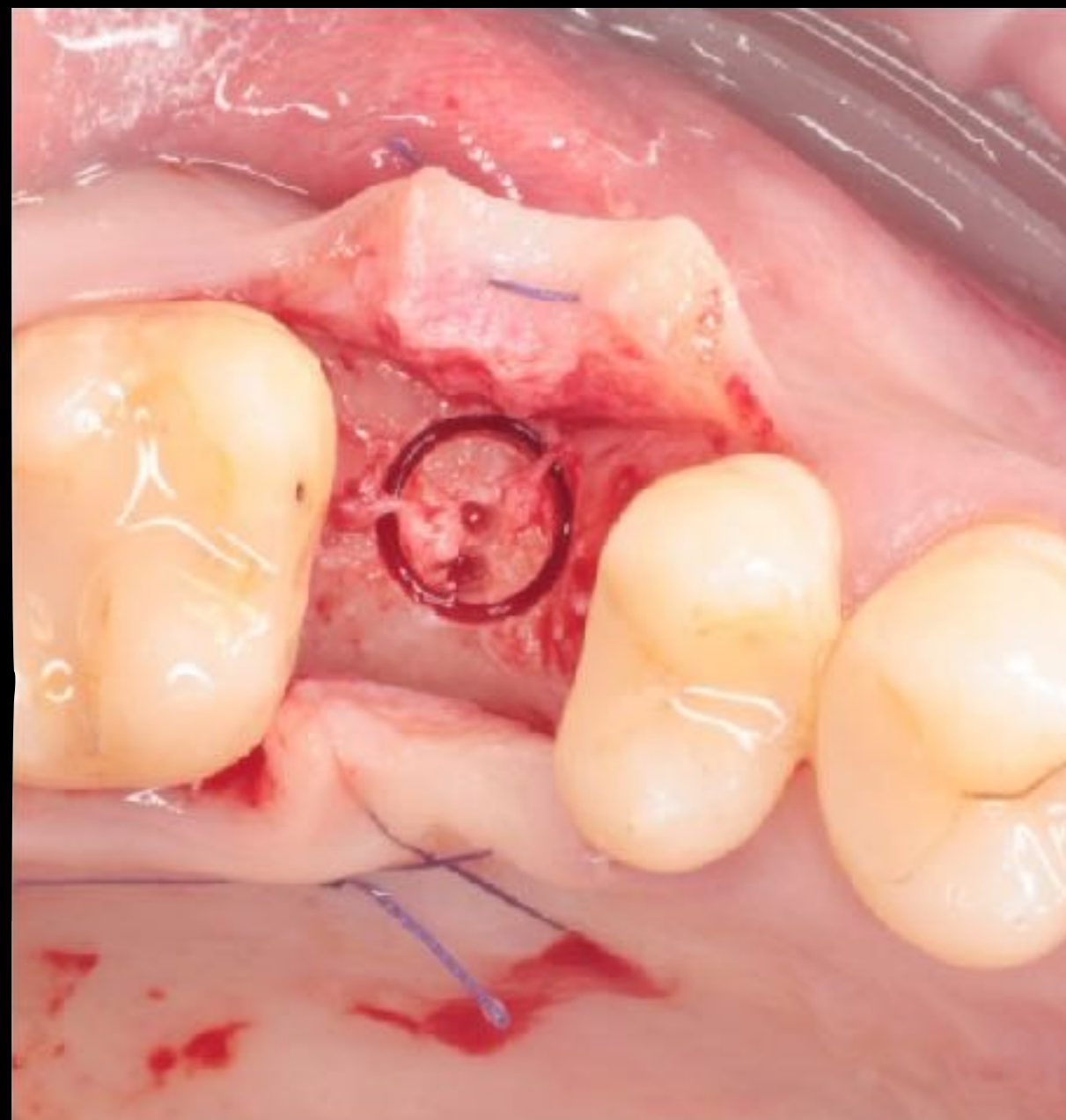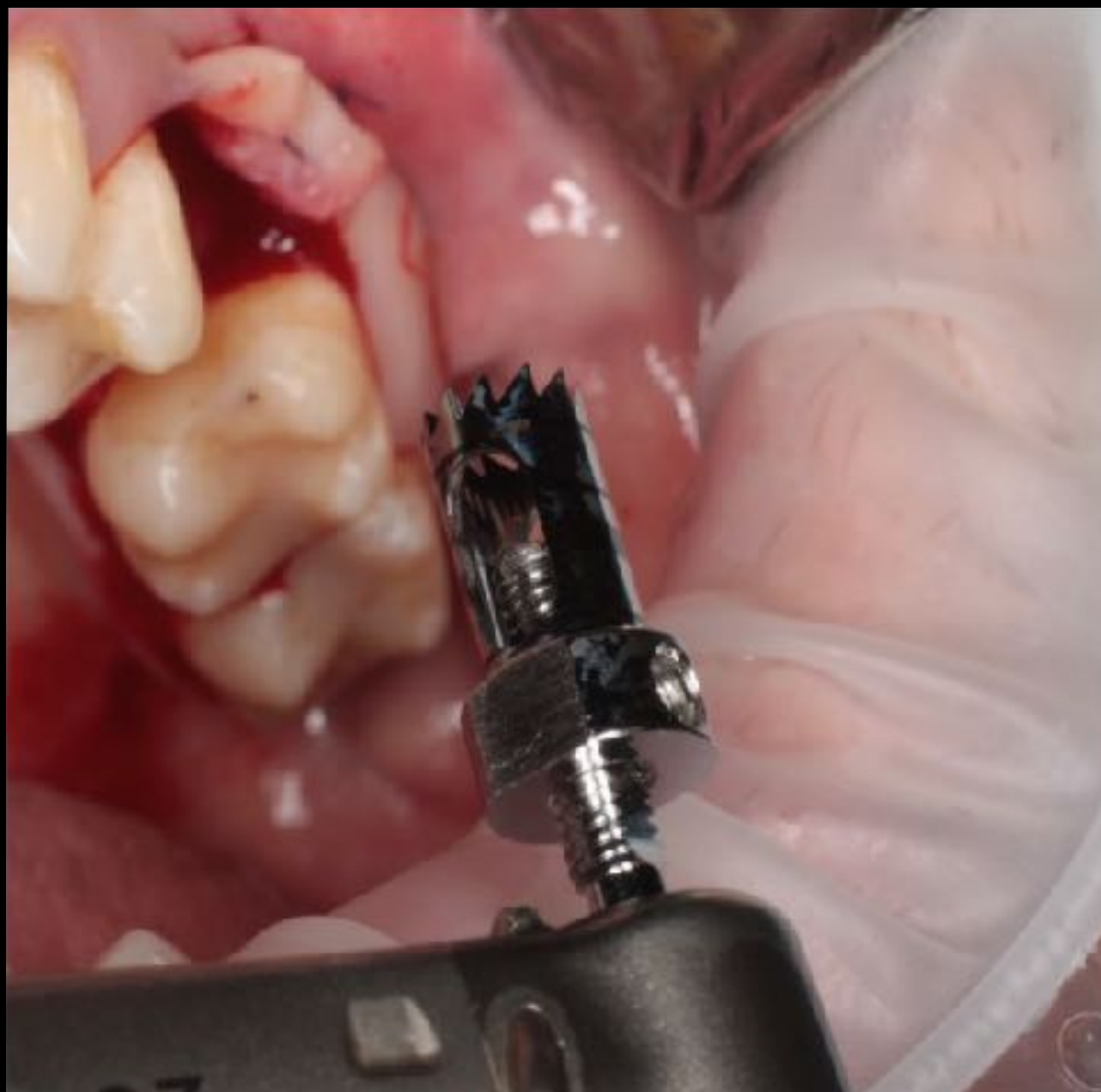


Case by Dr. Alexander Lysovov

Clinical Case: Drilling using ASBE trephine bur until 1-2mm of bone is left




Case by Dr. Alexander Lysovov

Clinical Case: Extraction of bone tissue by slightly tilting bur

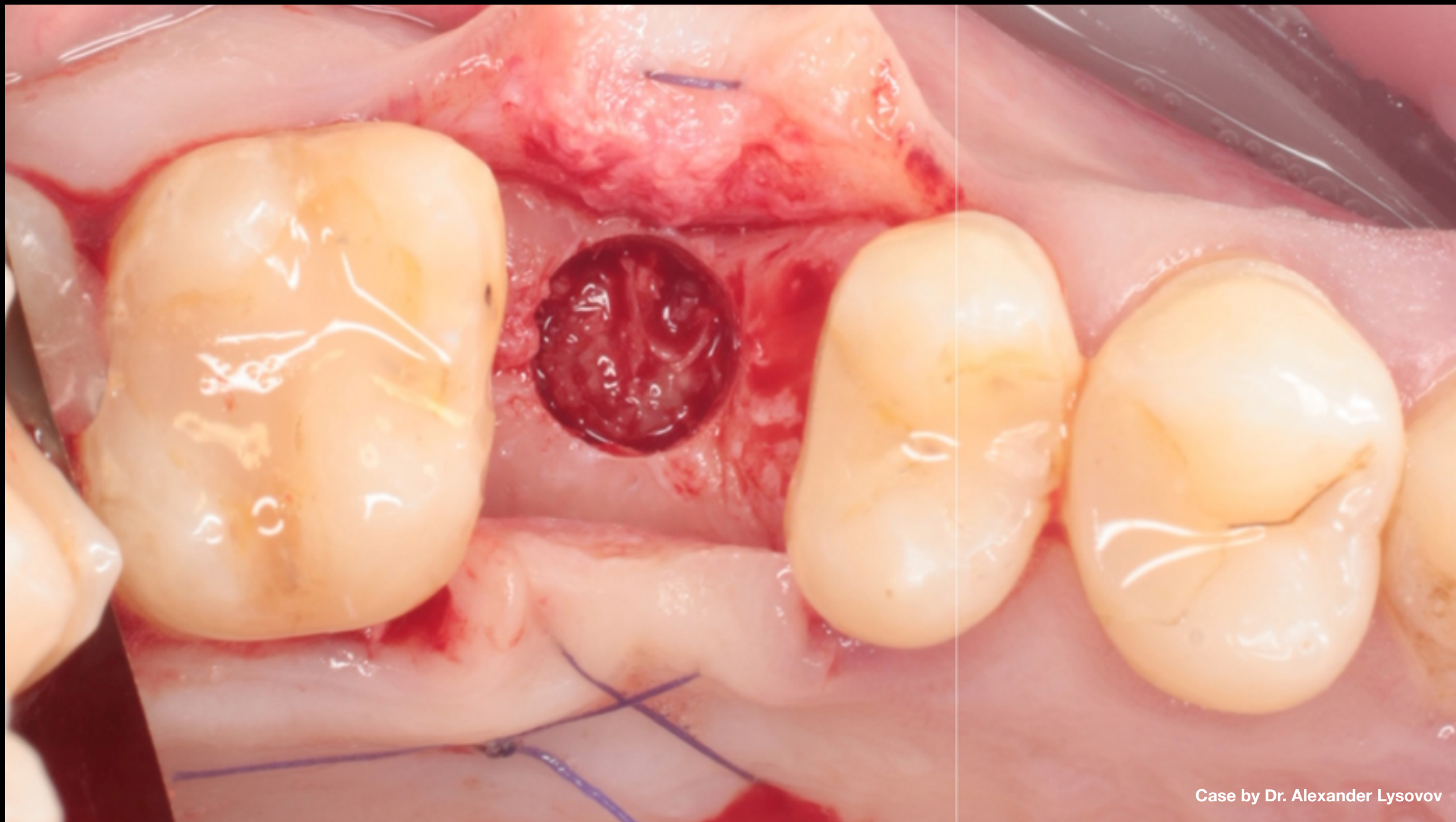

Case by Dr. Alexander Lysovov

Clinical Case: Extraction of bone tissue by slightly tilting bur

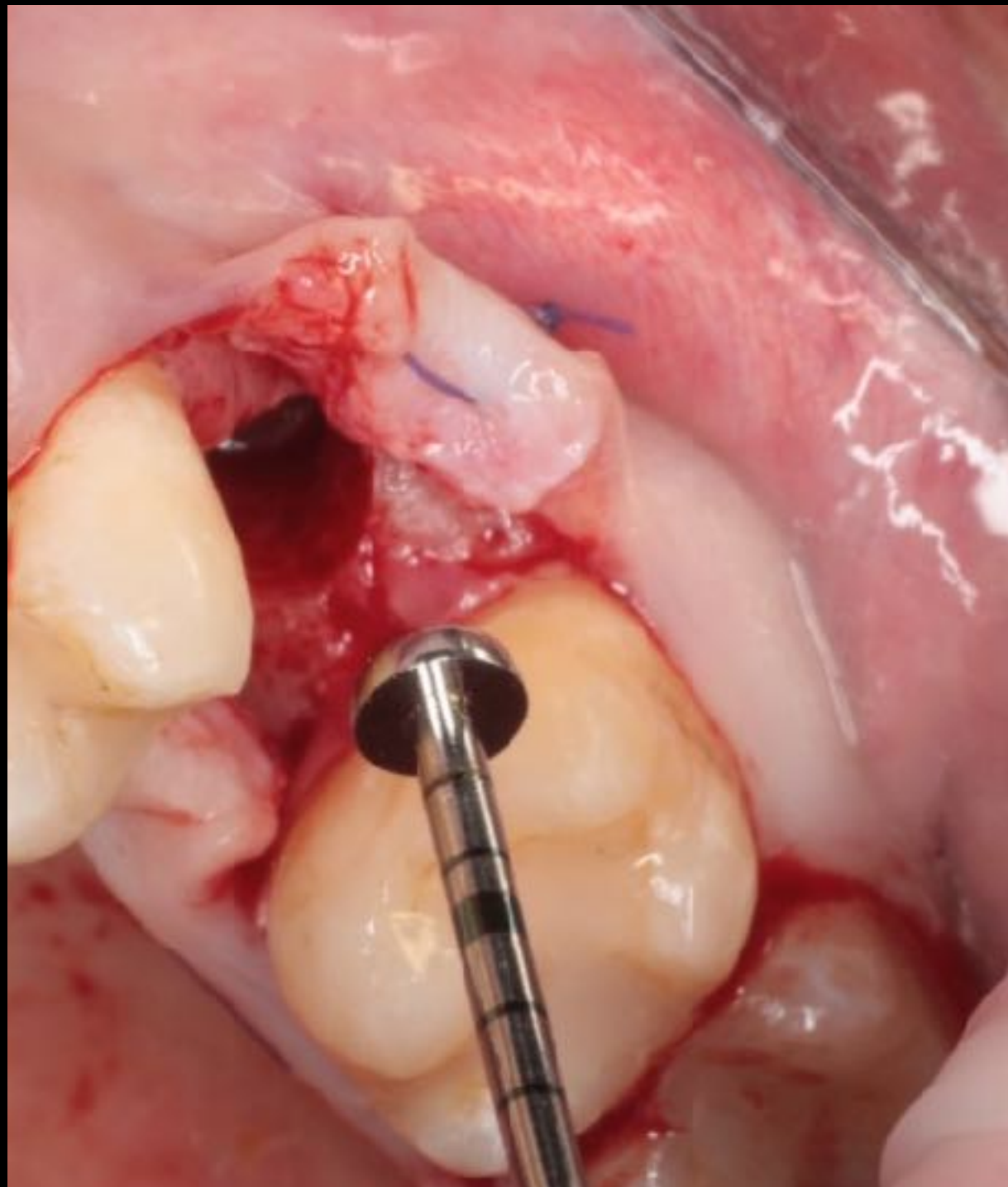



Case by Dr. Alexander Lysovov

Clinical Case: Schneider membrane detachment using Mushroom & Cobra





Case by Dr. Alexander Lysovov

Clinical Case: Schneider membrane detachment using Mushroom & Cobra

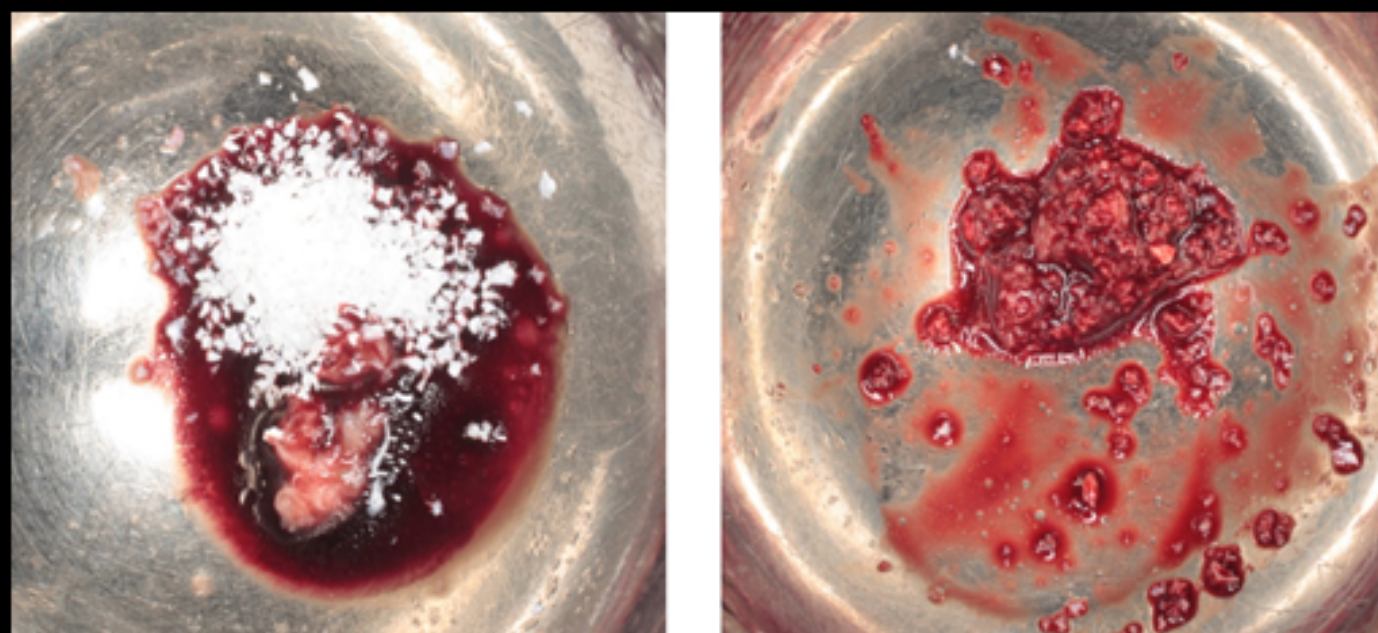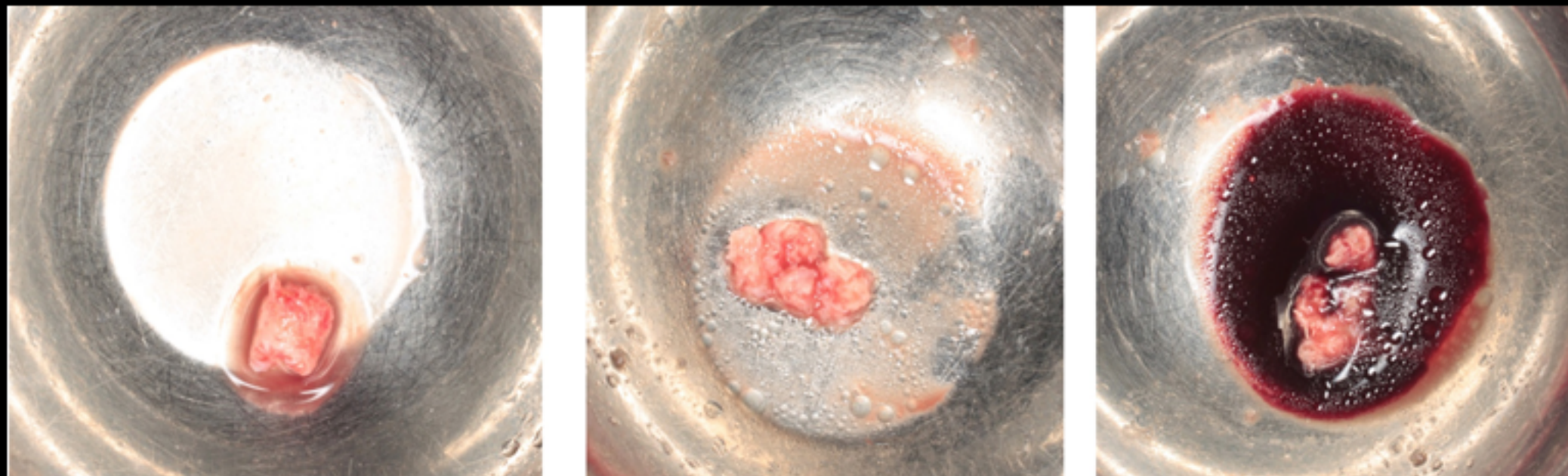

Case by Dr. Alexander Lysovov

Clinical Case: Creation of biologic osteo-plastic material

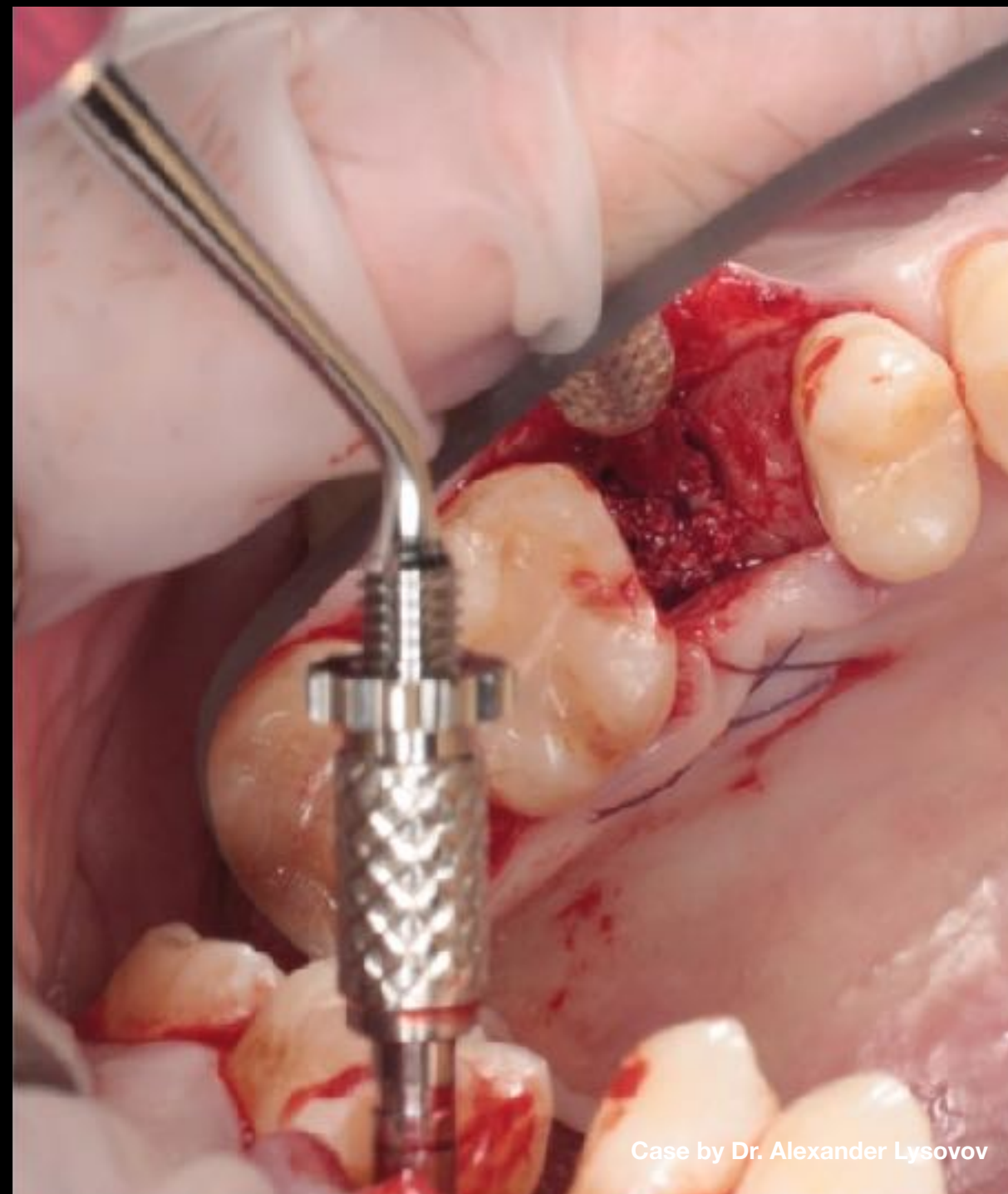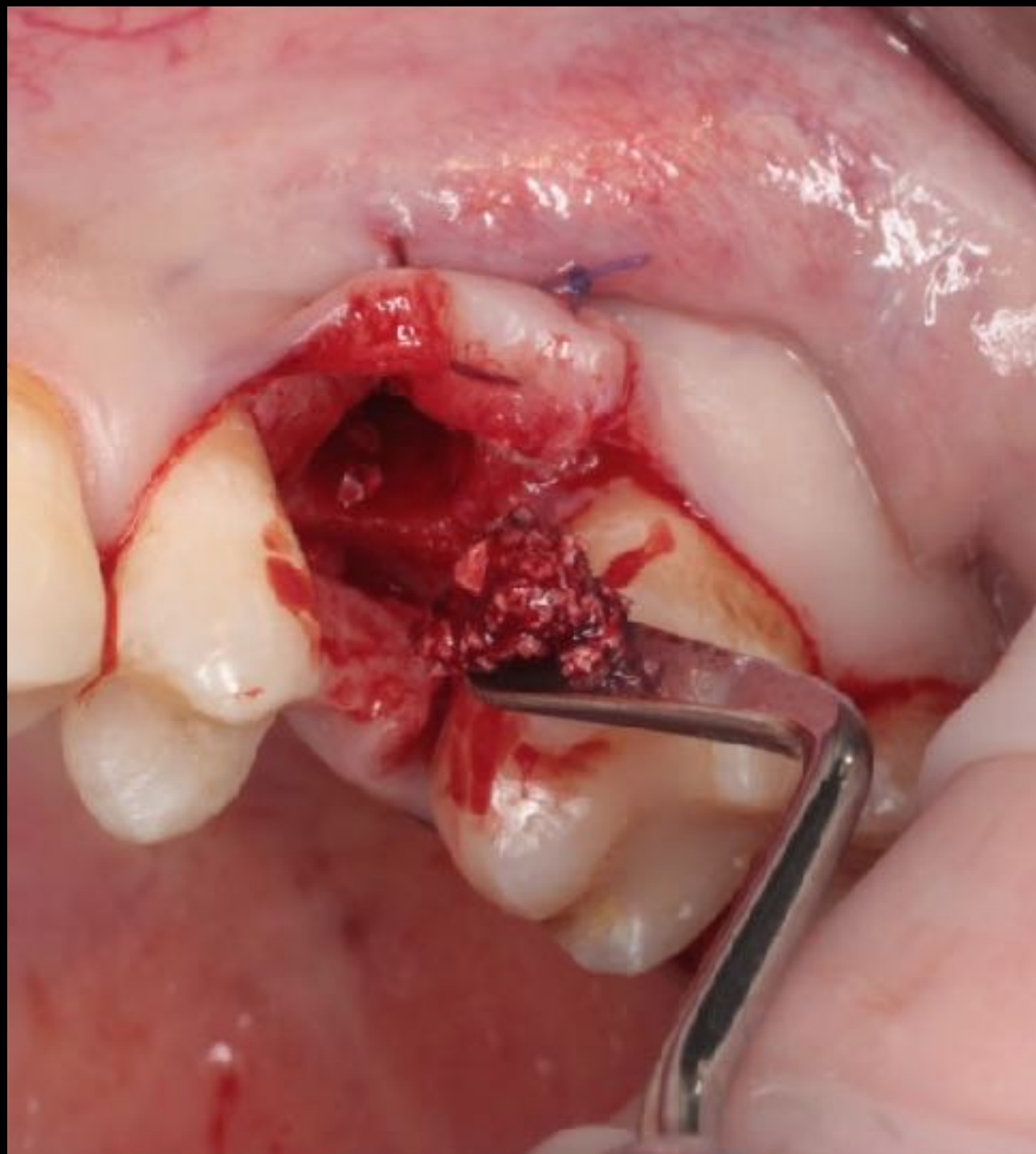

Case by Dr. Alexander Lysovov

Clinical Case:

- Graft harvested bone & alloplastic material using Spreader
- Adjust Condenser stopper & fill bone material to desired depth by pressing




Case by Dr. Alexander Lysovov

Clinical Case: Blue-colored Xpeed surface guarantees cleanliness



Clinical Case: Placement of AnyRidge implant

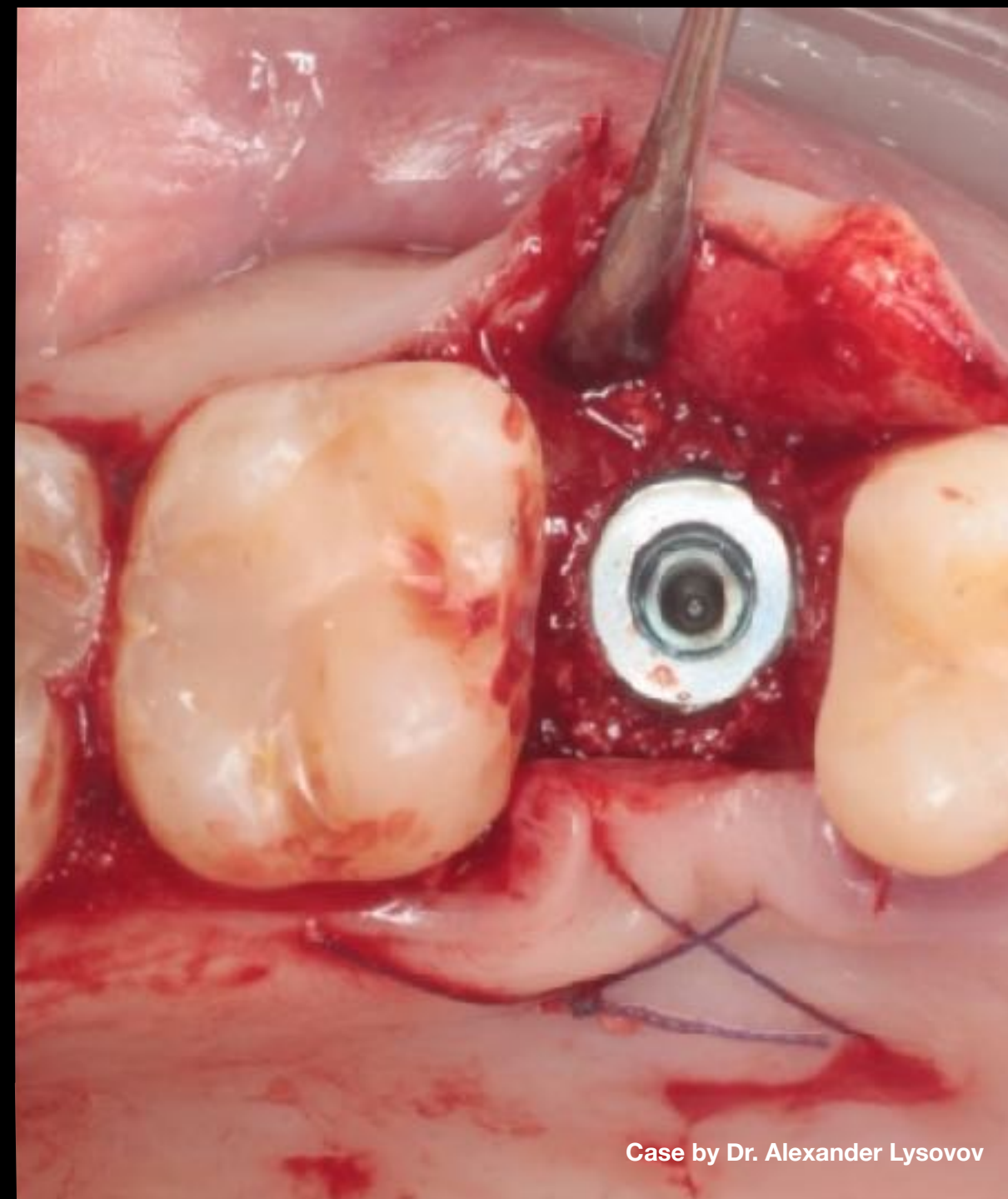


Case by Dr. Alexander Lysovov




Case by Dr. Alexander Lysovov


Clinical Case

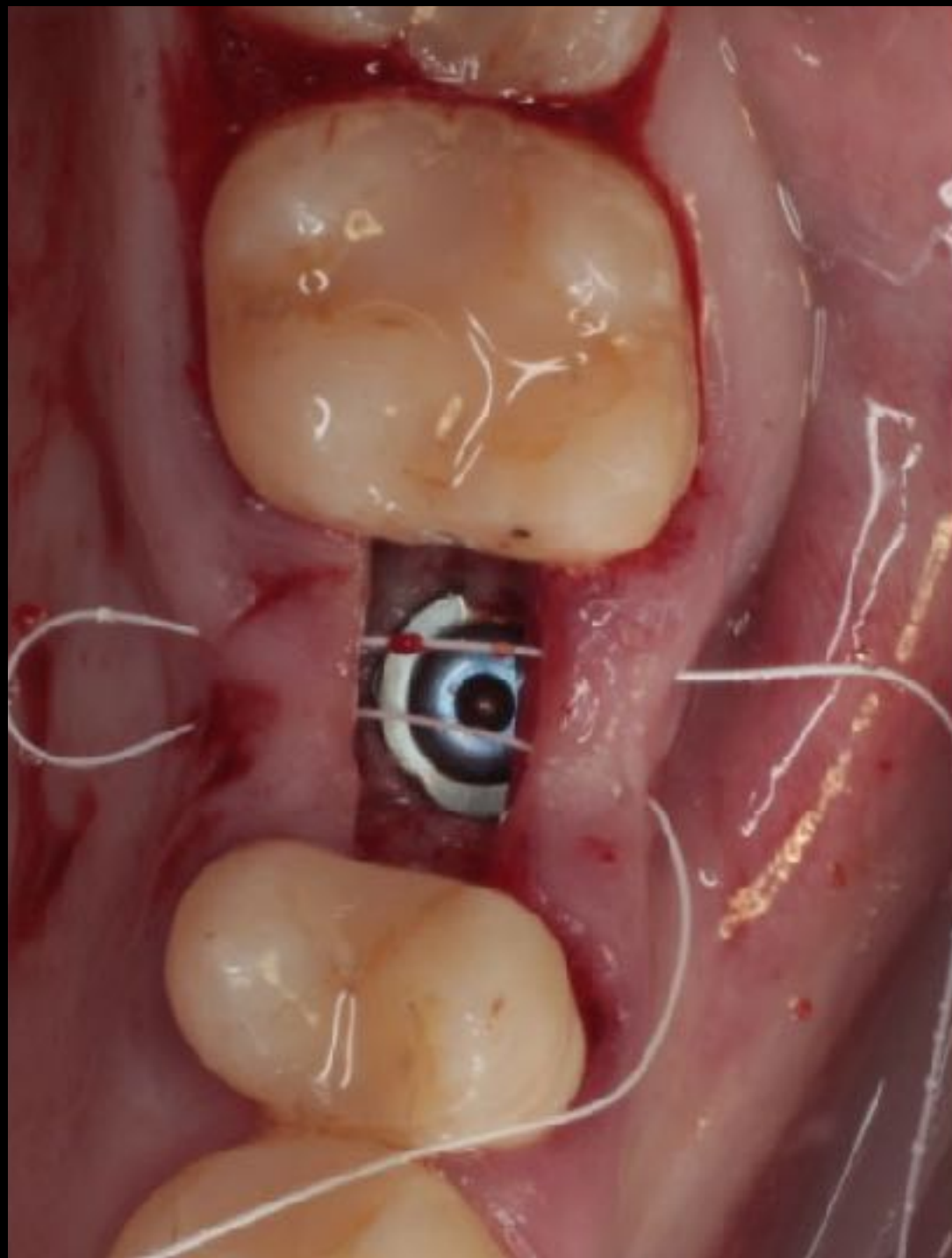

Clinical Case: Wound closure

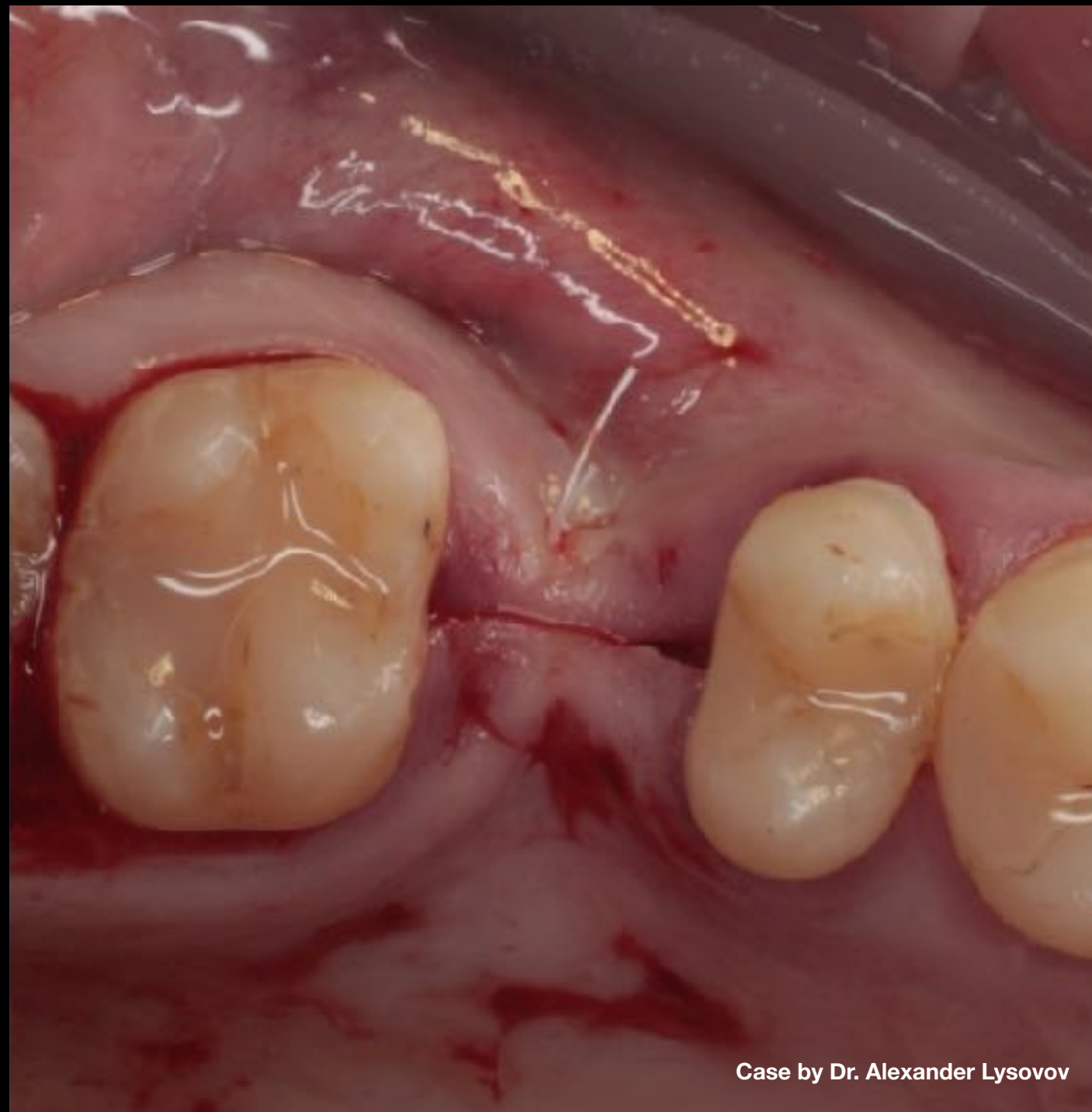


Case by Dr. Alexander Lysovov

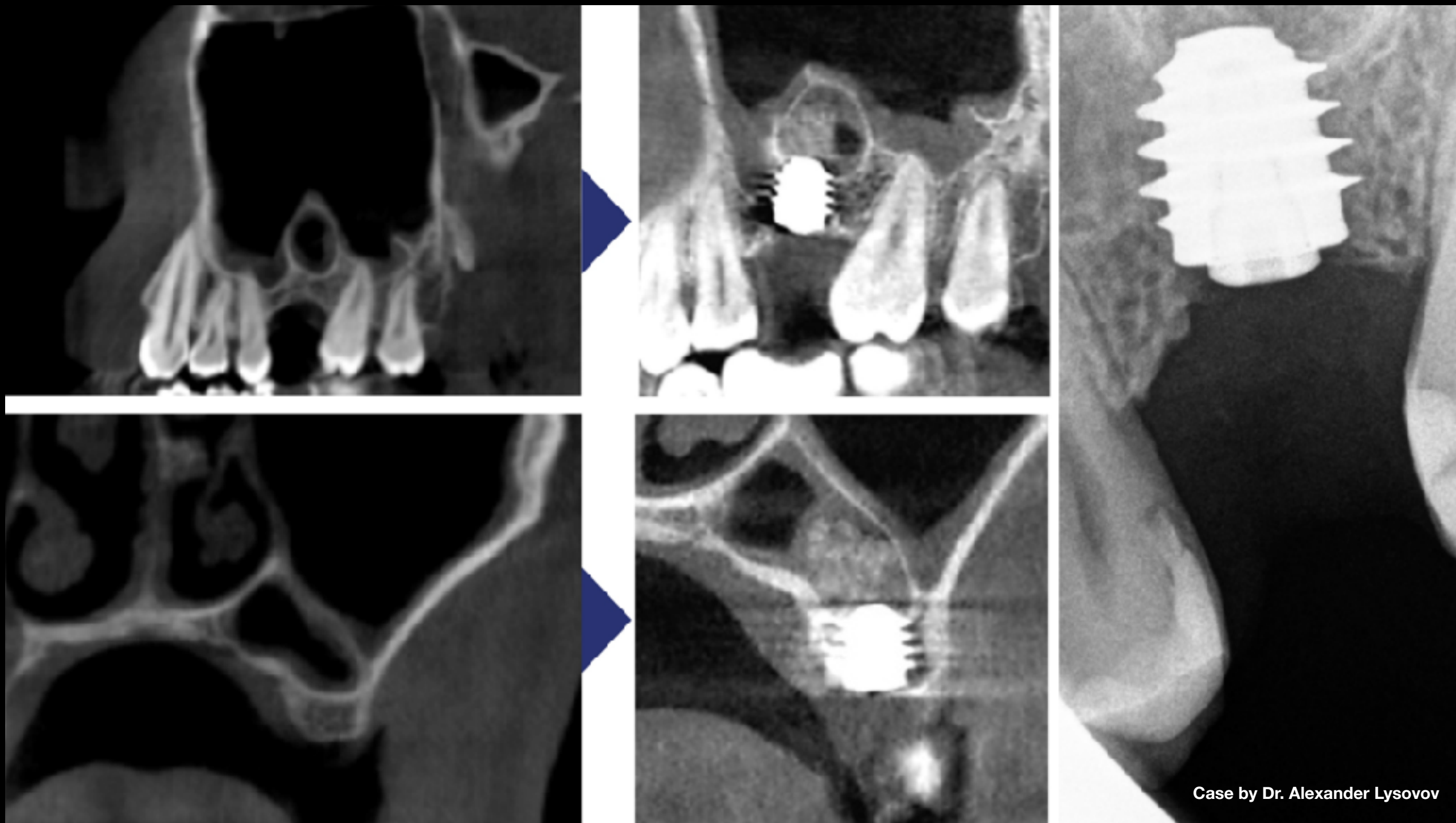

Case by Dr. Alexander Lysovov

Clinical Case: Post-operation tomography

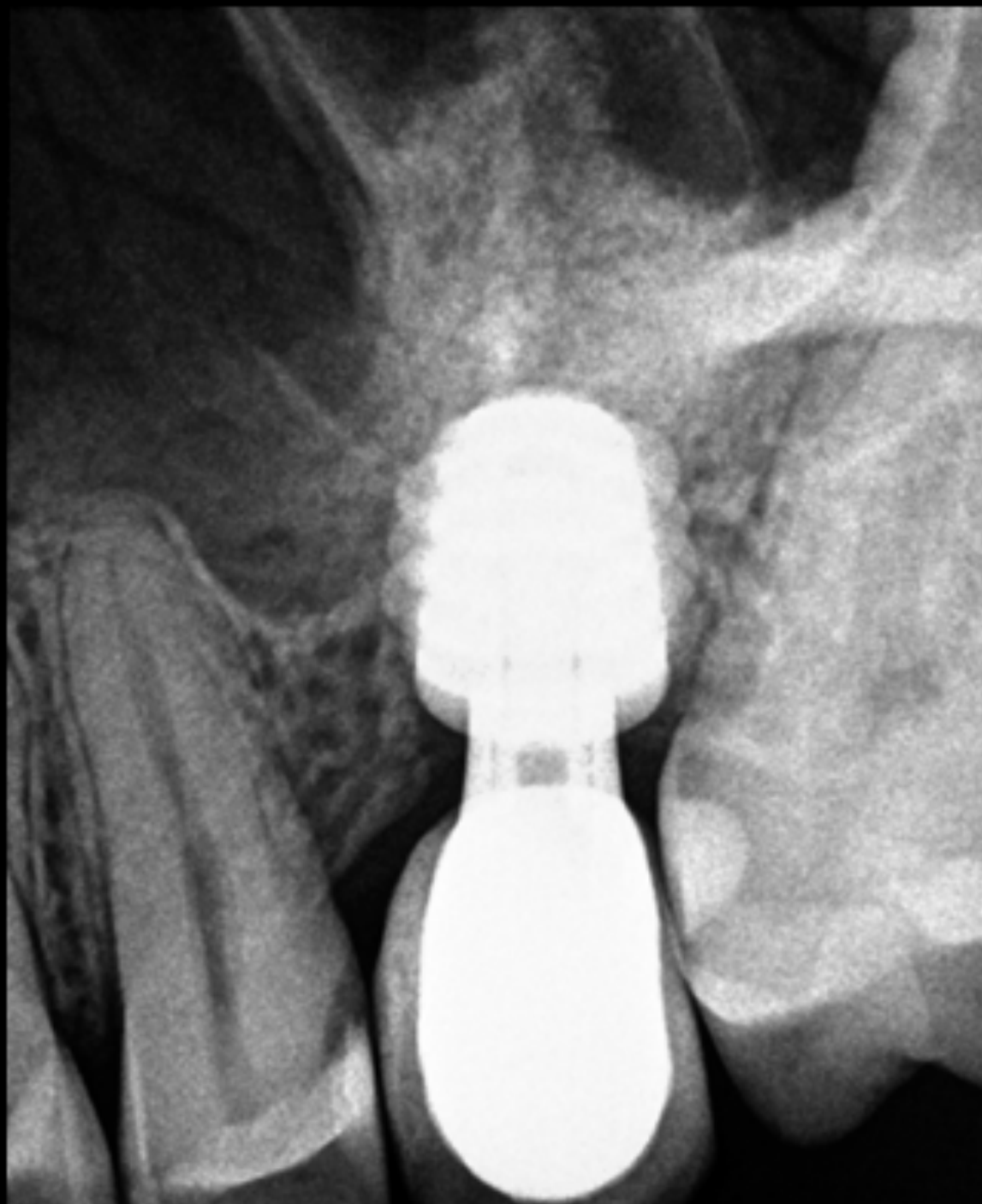

Clinical Case: 1 year after prosthesis

Case by Dr. Alexander Lysovov

