





Surgery using # 35 implant placement case, R2GATE GUIDE and ONE-DAY Prosthesis (Zirconia customized abutment and PMMA Crown).

Dr. Kwang Bum Park

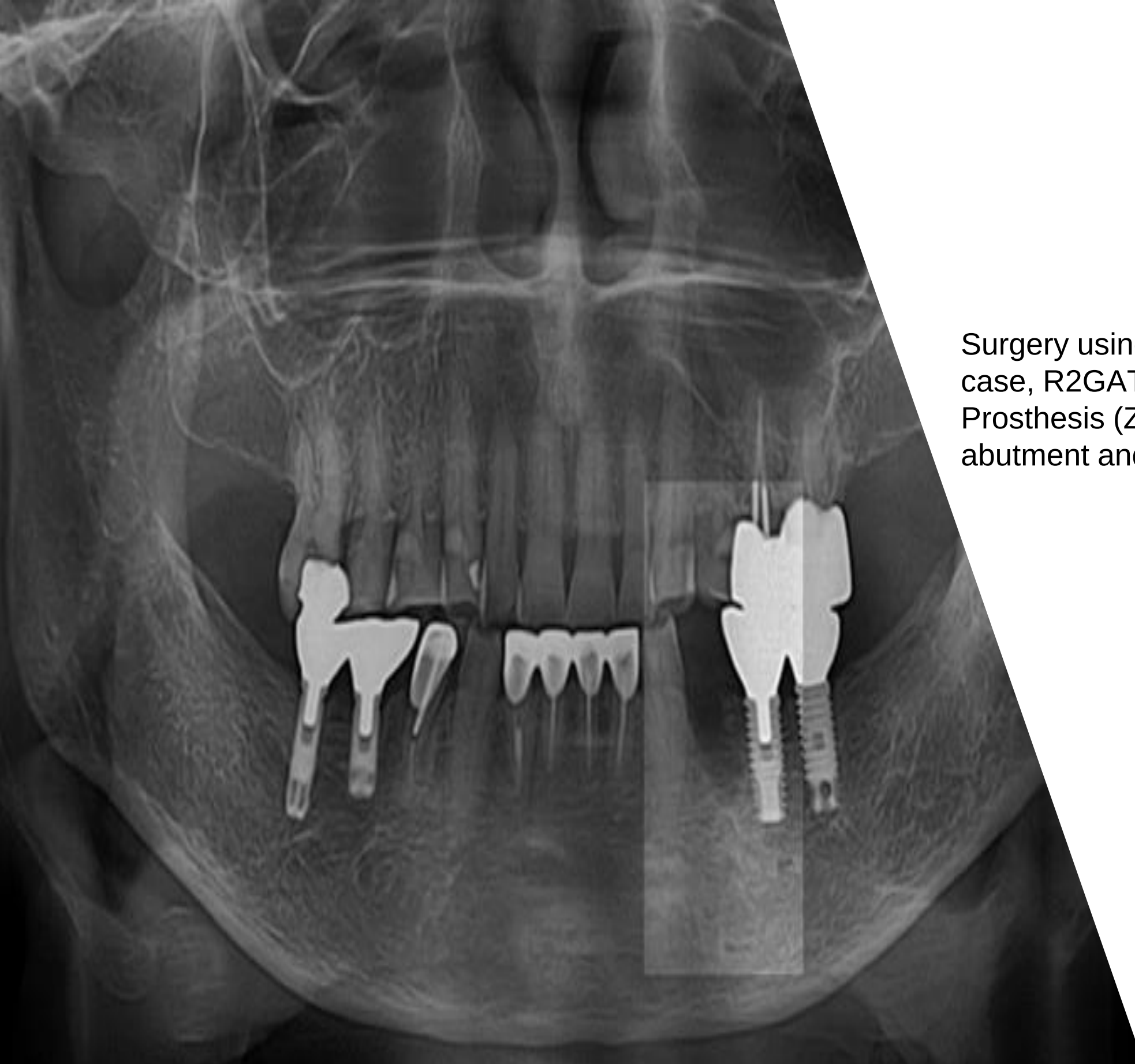

Case by Dr. Kwang Bum Park

Clinical Case: # 35 with tooth extraction socket

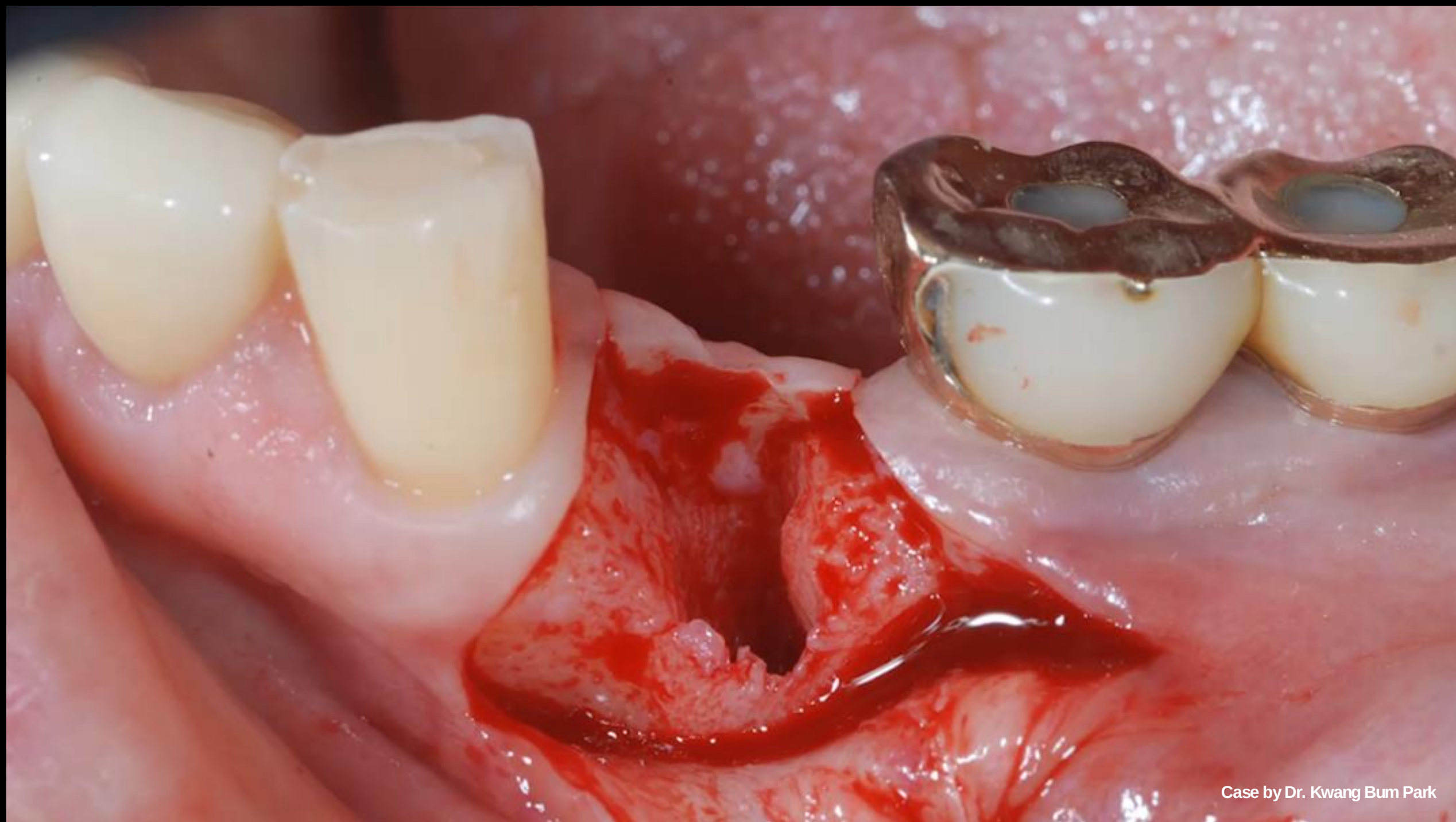



Case by Dr. Kwang Bum Park

Clinical Case: Flap open




Case by Dr. Kwang Bum Park

Clinical Case: ONE-DAY Prosthesis (PMMA SCRP TYPE Crown)

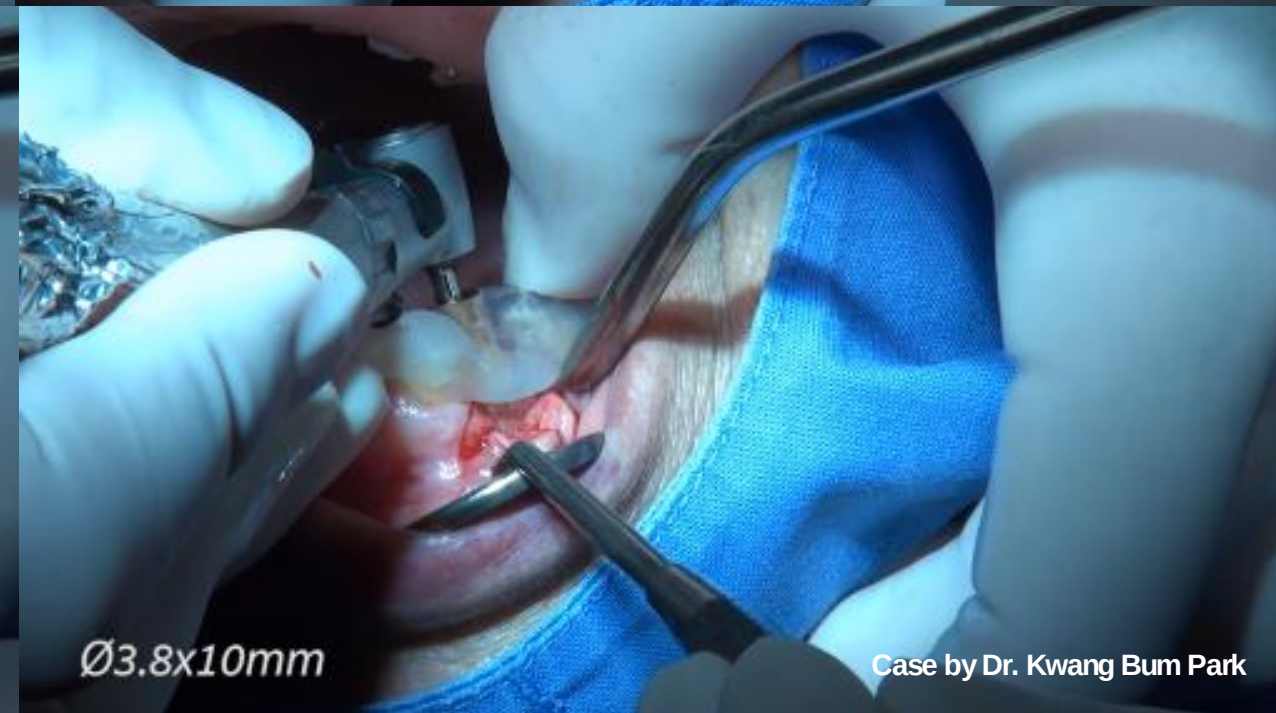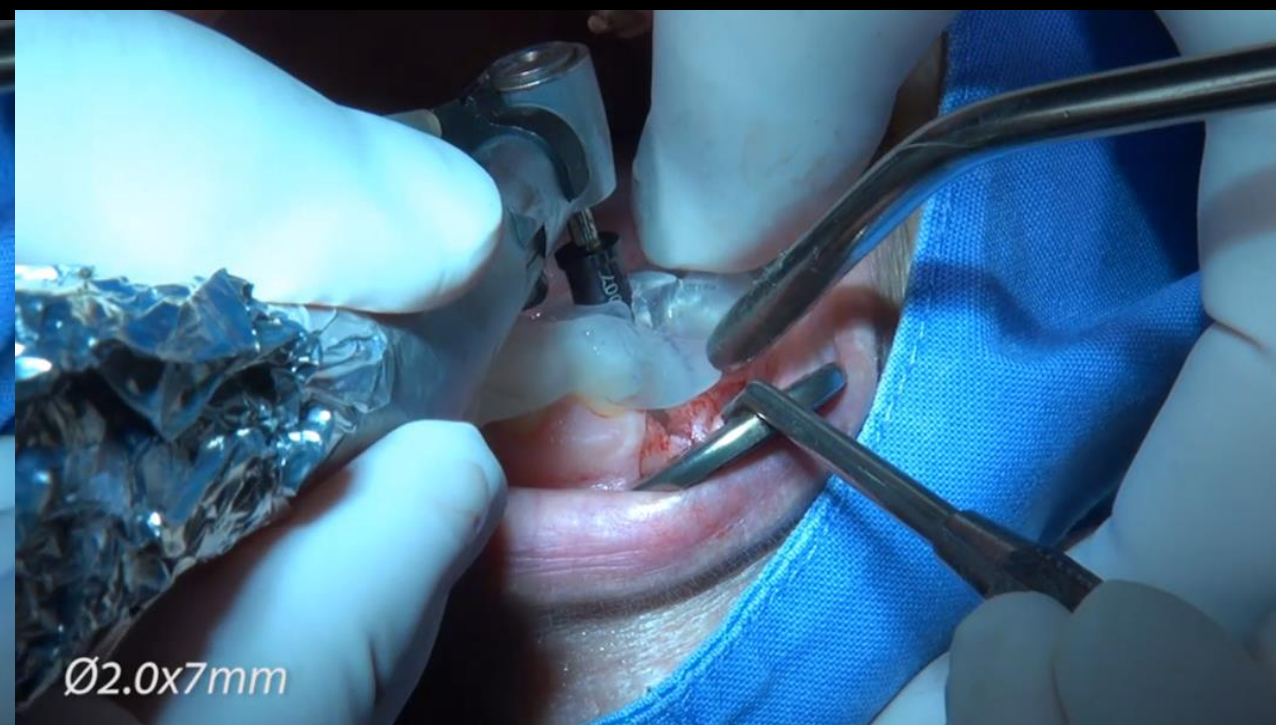


Clinical Case: Drilling according to drilling sequence using bone quality through digital eye: Drilling sequence(Initial drill, Second drill, 2.0x7mm, 2.0x10mm, 2.5x10mm, 2.8x10mm, 3.3x10mm, 3.8x10mm)

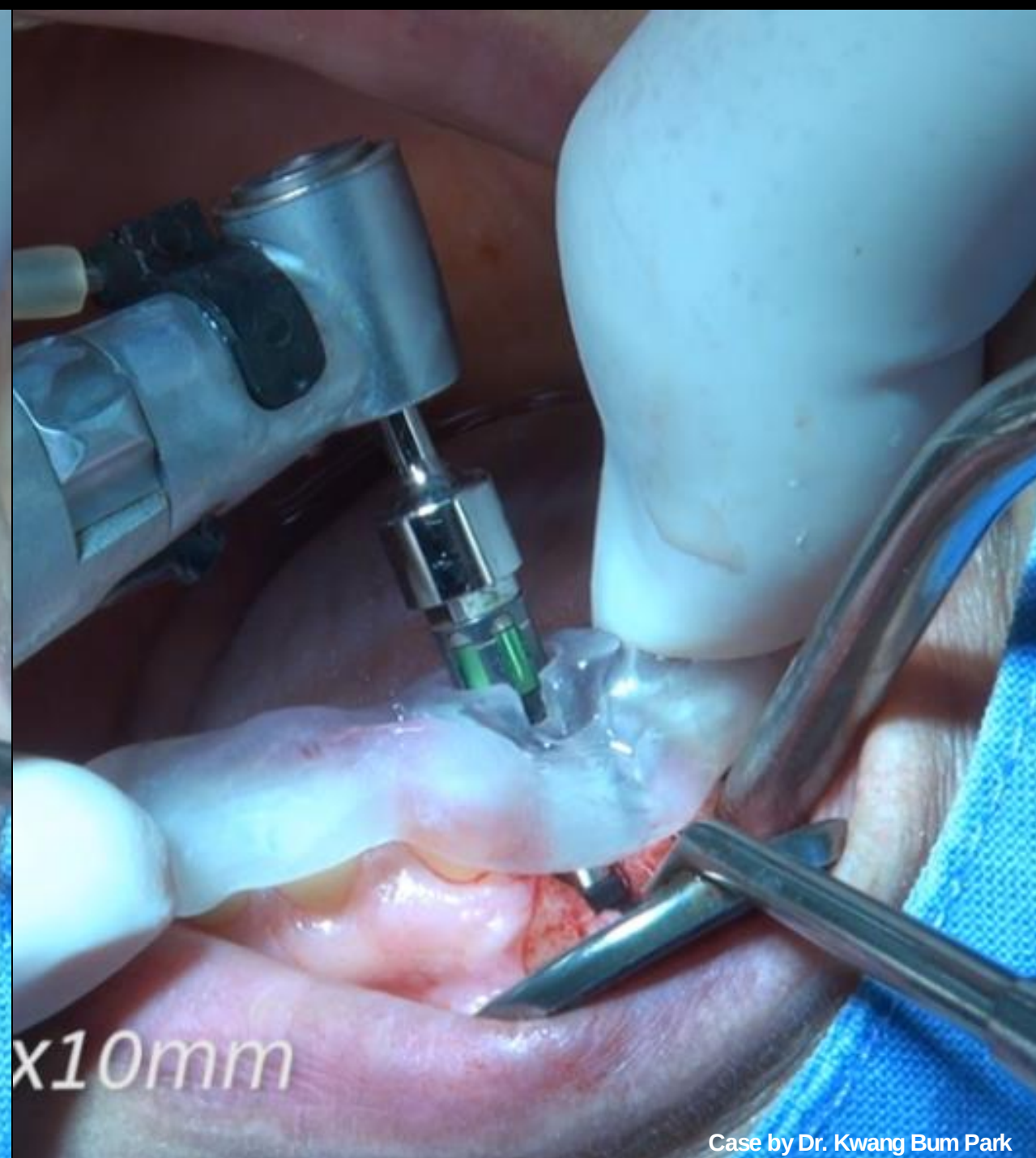


Case by Dr. Kwang Bum Park

Clinical Case: AnyRidge 4.5x10mm Implant Placement

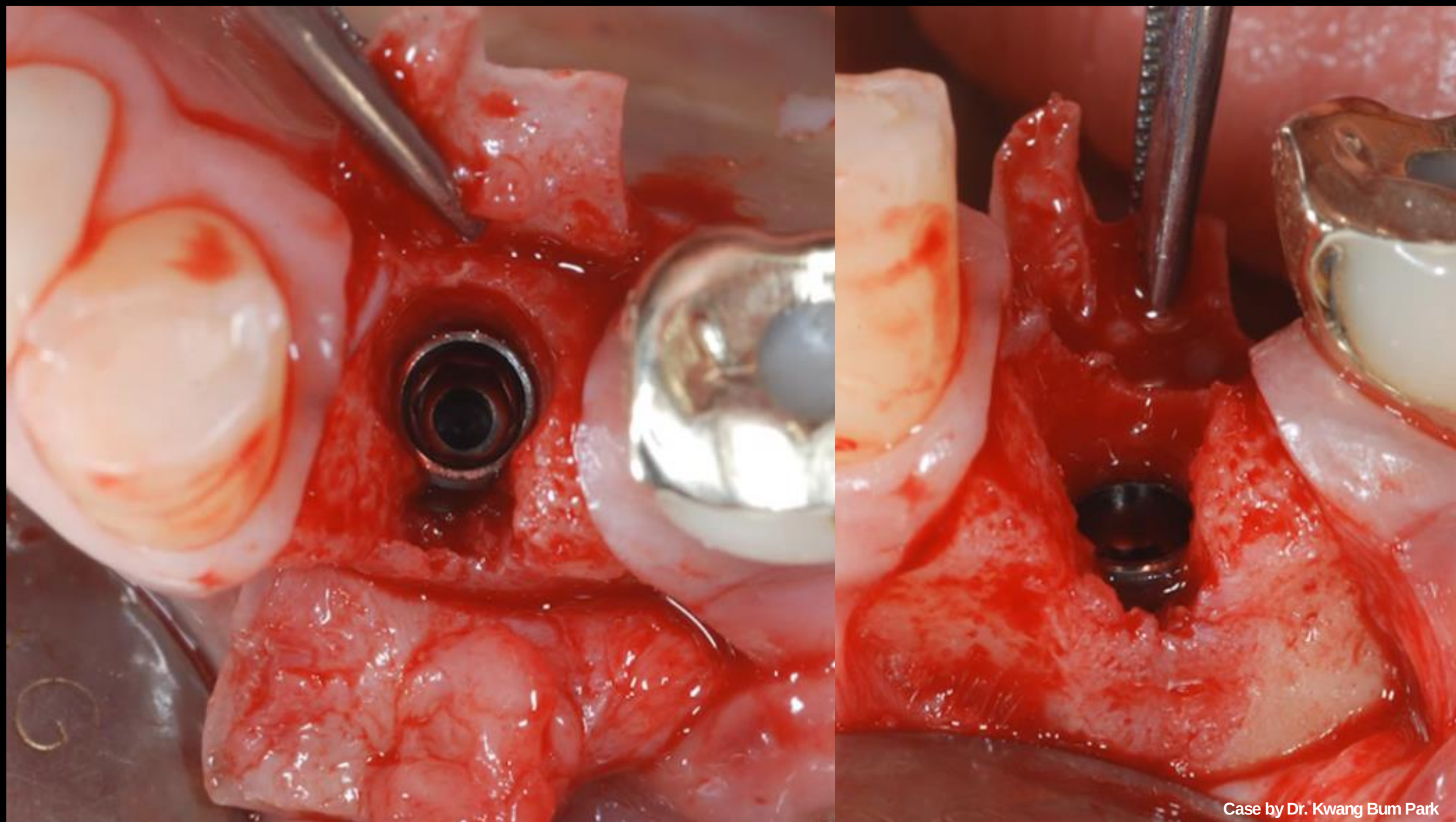

Case by Dr. Kwang Bum Park

Clinical Case: Check the implant placement.

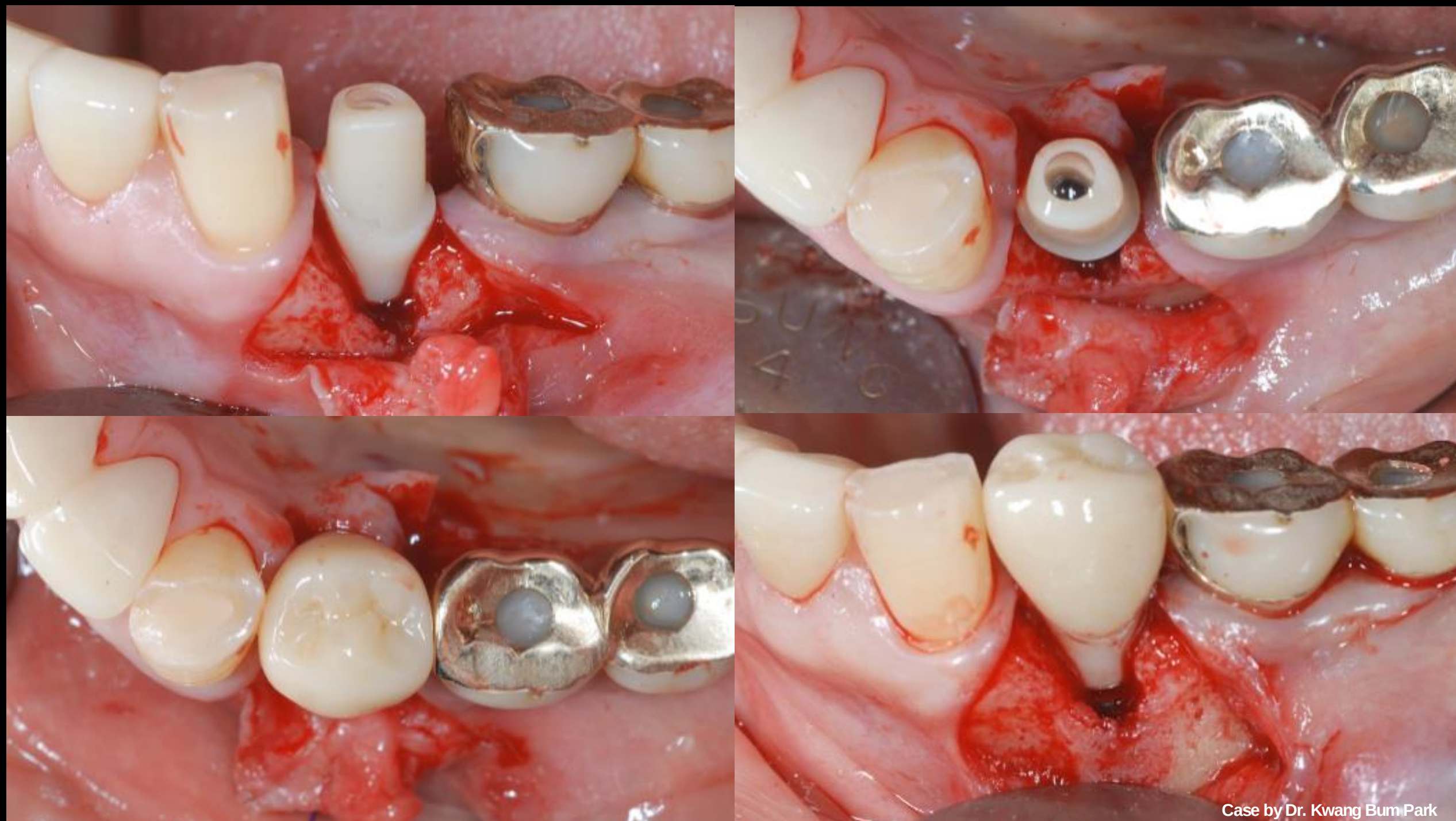


Clinical Case: Set the Zirconia customized abutment and PMMA Crown




Case by Dr. Kwang Bum Park

Clinical Case: GBR using Mega-Oss

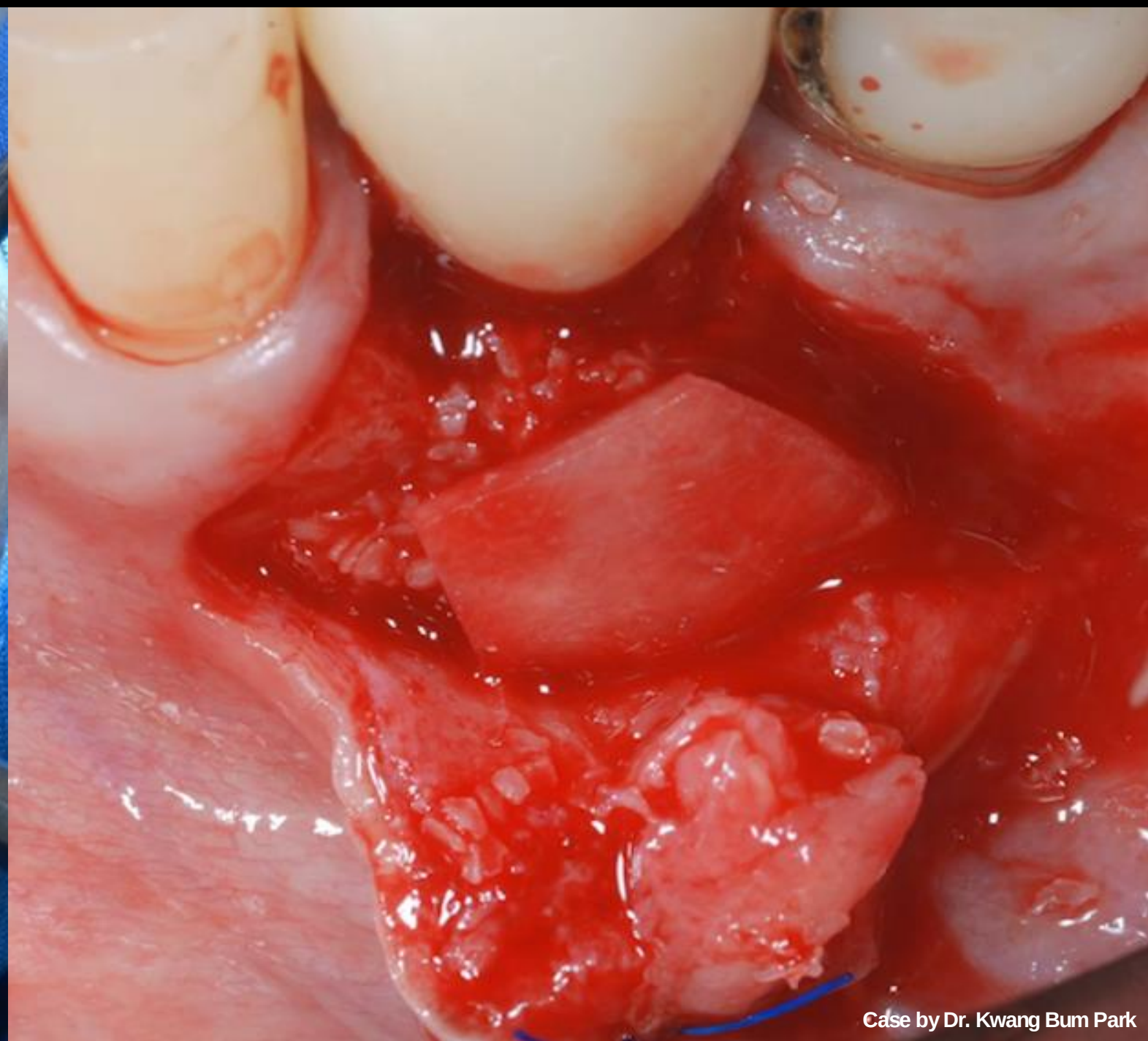

Case by Dr. Kwang Bum Park

Clinical Case: Protect bones implanted with lyoplants

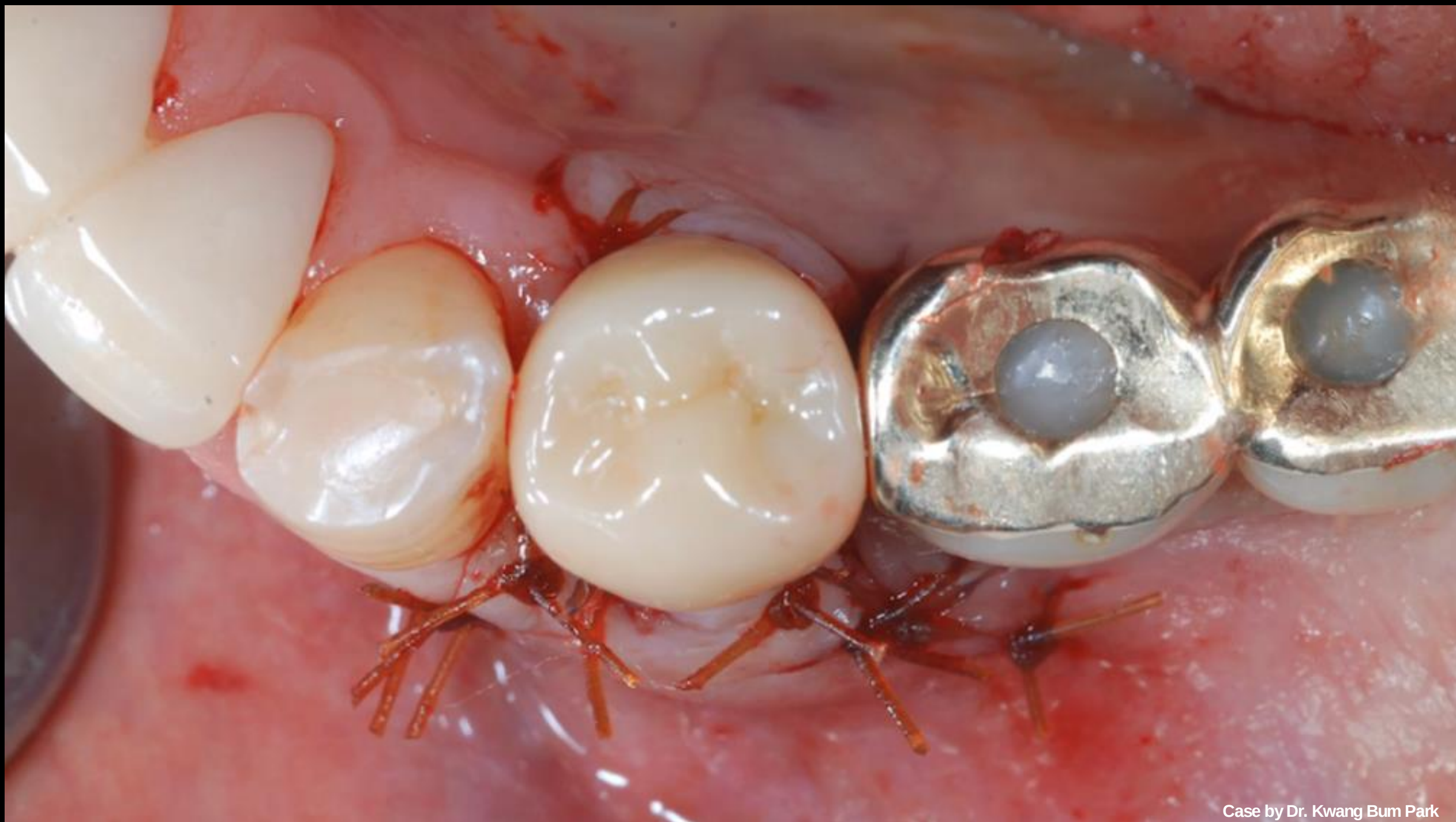


Case by Dr. Kwang Bum Park

Clinical Case: Suture finished (occlusal surface)



Case by Dr. Kwang Bum Park

Clinical Case: Suture finished (narrowing)

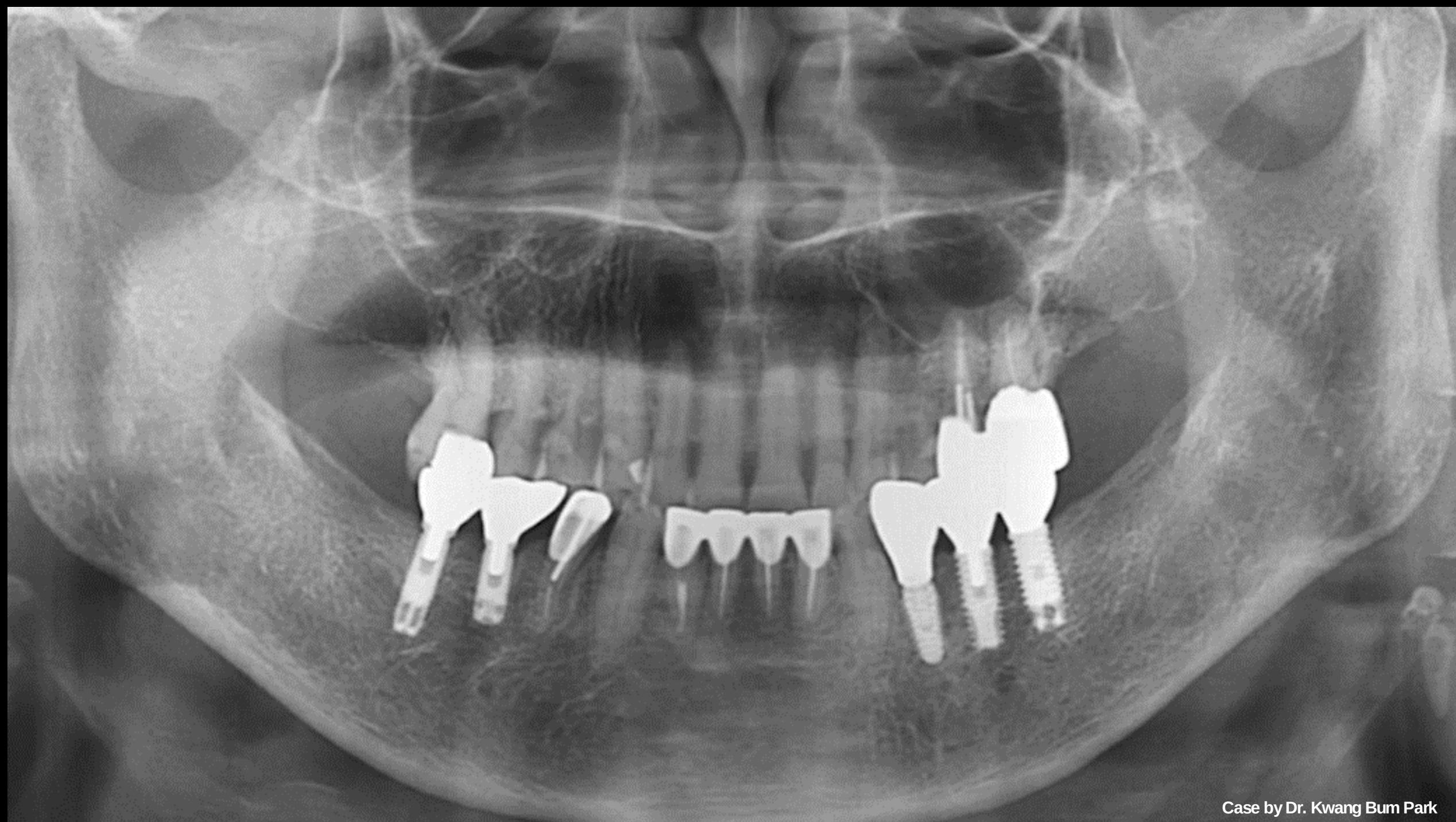

Case by Dr. Kwang Bum Park

Clinical Case: Panoramic picture after surgery

